Ziegelei 1 72336 Balingen Germany

Tel. +49 7433 9933-0 Fax +49 7433 9933-146

info@kern-sohn.com www.kern-sohn.com







# **MICROSCOPES & REFRACTOMETERS**

for laboratory, industry and food





















2007














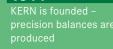















manufacturer (2009/23/EC)



Significant anniversary year! Construction of Ziegelei 3.0, 25 years of accredited DKD extension of administration







**Taras Mikitisin** Tel. +49 7433 9933-143 Fax +49 7433 9933-29143 Mobil +49 171 5590115

mikitisin@kern-sohn.com

#### KERN Models A-Z

OBB-C	0.0	ORT-1
OBD-0		OSE-4
OBE-10 · OBE-11	24	OSE-41
OBE-12 · OBE-13	15 16	OSE-42
OBE-S	73	OSF-4G
OBF-1		OZB-H
OBF-S		OZB-IR
OBL-12 · OBL-13	20	OZB-M
OBL-14 · OBL-15	22	OZB-UE
OBL-S	74	OZB-UP
OBN-13	26	OZC-5
OBN-14	28	OZG-4
OBN-S		OZL-44
OBS-1		OZL-45
OBT-1		OZL-45R
OCM-1		OZL-46
ODC-2		OZL-9
ODC-8		OZL-S
ODC-9		OZM-5
OKM-1	34	OZM-9
OKN-1		OZM-S
OKO-1		OZP-5
OLM-1		OZP-S
OPM-1	41	OZR-5
OPN-1		OZS-5
OPO-1		VIS (OXM-9)
ORA		Checklist for your microscop
ORF		Checklist for your refractome
ORF-E	106	•

## Do you have questions about our products? Our customer consultants will be pleased



Maren Möwert Tel. +49 7433 9933-132 Fax +49 7433 9933-29132

Mobil +49 151 46143240 maren.moewert@kern-sohn.com



Fax +49 7433 9933-29118 Mobil +49 [0]151 44611168 florian.baechle@kern-sohn.com

GR, CY, BG, HU, RO, SK, CZ, AL, Ex-Yugoslavia, GUS Aleksandar Delić

Tel. +49 7433 9933-120 Fax +49 7433 9933-29120 Mobil +49 151 11196197 aleksandar.delic@kern-sohn.com

> Nordamerika, Near- & Middle East, Africa, Asia, Oceania, TR Corinna Klaass

Tel. +49 7433 9933-215 Fax +49 7433 9933-29215 Mobil +49 [0]151 44568364 klaass@kern-sohn.com

> **Technical Service** Ralf Kurschat

technical service and maintenance:

For special questions please contact our

Manager KERN Optics

Tel. +49 7433 9933-155

Fax +49 7433 9933-29195

daniel.junger@kern-sohn.com

Tel. +49 7433 9933-173

samuel.haasis@kern-sohn.com

Fax +49 7433 9933-29173

Technical sales KERN Optics &

**SAUTER Measuring instrumente** 

Tel. +49 7433 9933-306

Fax +49 7433 9933-29306

ralf.gutbrod@kern-sohn.com

Produkt specialist SAUTER

Tel. +49 7433 9933-205

Mobil +49 171 3059946

Quick help for all questions regarding our

Fax +49 7433 9933-29205

nadine.wenzler@kern-sohn.com

Nadine Wenzler

Measuring Instruments & KERN

Product management KERN Optics

Daniel Junger

Samuel Haasis

Ralf Gutbrod

**OPTICS-Team:** 

Tel. +49 7433 9933-166 Fax +49 7433 9933-195 ralf.kurschat@kern-sohn.com



60

\_78/79

### KERN Pictograms




360° rotatable microscope head



Monocular Microscope For the inspection with one eve



Binocular Microscope For the inspection with both eyes

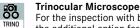

For the inspection with both eyes and the additional option for the connection of a camera



With high numerical aperture for the concentration and the focusing of light

For pictures bright and rich in contrast



LED illumination Cold, energy saving and especially

Halogen illumination



Incident illumination For non-transparent objects



Transmitting illumination For transparent objects



Fluorescence illumination For stereomicroscopes



Fluorescence illumination for compound microscopes With 100 W mercury lamp and filter



Fluorescence illumination for compound microscopes With 3 W LED illumination and filter



Darkfield condenser/unit

For a higher contrast due to indirect

Phase contrast unit



Polarising unit To polarise the light



Infinity corrected optical system



Zoom magnification For stereomicroscopes



**Abbreviations** 


Adapter for the connection of a C-Mount camera to a trinocular microscope

high eye point for wearers of glasses)

Long Working Distance

Numerical Aperture

Parallel optical system

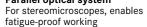



Integrated scale In the eyepiece



SD card For data storage



USB 2.0 digital camera For direct transmitting of the picture to a PC



USB 3.0 digital camera For direct transmitting of the picture to a PC



WLAN data interface: For transmitting of the picture to a mobile display device



HDMI digital camera For direct transmitting of the picture to a display device



PC software

To transfer the measurements from the device to a PC.



Automatic temperature compesation For measurements between

10 °C and 30 °C



Protection against dust and water splashes IPxx



The type of protection is shown by the pictogram.



**Battery operation** Ready for battery operation. The battery type is specified for each device.



Battery operation rechargable Prepared for a rechargable battery operation



Mains adapter 230V/50Hz in standard version for EU. On request GB, AUS or USA version.



Power supply

Integrated in microscope. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.



Package shipment

The time required to manufacture the product internally is shown in days in the pictogram.

**SLR camera** Single-Lens Reflex camera Super Wide Field (Field number at least Ø 23 mm for 10× evepiece) Frames per second High (Super) Wide Field (Eyepiece with W.D. Working Distance

Wide Field (Field number up to Ø 22 mm for 10× eyepiece)

### Why you should choose a KERN microscope now!

weighing and measuring technology. This claim is the driving force for the development of our microscope and refractometer ranges.

Thanks to consistent customer focus paired with smart ideas and the latest available technology we are proud to be suppliers of high-quality, durable top microscopes and refractometers, which help you to be as efficient as possible in your daily work.

When developing our microscopes we have concentrated on the very best optical quality and have used only high-quality glass and the latest technologies to achieve this. The high-quality Philips halogen and modern LED illumination produce razor-sharp images with high contrast and which will impress you with their brilliant true-colour display - you must have noticed this yourself.

#### Your advantages:

- all mechanical parts have been designed for a long service life
- special attention has been given to the ergonomy of our microscopes, as this allows the user to work for several hours in a comfortable position which does not cause fatigue
- our microscopes are fully-equipped and can be used immediately
- Highlight for 2020: die KERN camera software you will be amazed at how user-friendly and intuitive it is, a high-quality tablet camera as well as a comprehensive range of calibration services for refractometers
- · and much more...

For 175 years, KERN & SOHN has been synonymous with high precision Use our practical "Check list for microscopes and refractometers", which may help you to quickly determine specifications for the future instrument. Together with our KERN product specialist you can choose the right product for you.

> If there is no suitable product in the standard range, for example, then we will of course configure an individual microscope for you.

Our aim is to develop a market-driven product solution, so with our microscope and refractometer range, the saying holds true: good quality at a competitive price! This is what we stand for and work towards, every day!

With our current 2020 product range you can benefit from improved quality and a clear reduction in price, which we have been able to achieve through more efficient production methods and increased global sales of our microscopes and refractometers and of course we pass this straight on to you

Do you have any questions about our range of microscopes and

Your KERN customer consultants are available at any time to help you further.

I hope that you enjoy working efficiently with our KERN Optics

## Your advantages

#### fast

- · 24 hour dispatch service for products in stock - order today, on its way tomorrow
- · Sales & service hotline from 8:00 am to 5:00 pm

+49 7433 9933-0

Service hotline

Online shop

+49 7433 9933-199

www.kern-sohn.com

#### reliable

- Up to 3 years warranty
- Certified QM system **DIN EN ISO 9001**

#### versatile

- · One-stop shopping: from microscope through to refractometer - everything from one supplier
- Quick as a flash, find the product you want with the "Quick-Finder" at www.kern-sohn.com



#### E-mail order info@kern-sohn.com



Fax Order +49 7433 9933-146



Our team of consultants will assist you from Monday to Friday

from 8:00 am to 5:00 pm

www.kern-sohn.com

Information on current product availability, product data sheets, user instructions, useful knowledge, technical glossary and much for you to download, practical topic areas, which will guide you to the right product in your industry as well as a microscopes and refractometers search engine.

#### Retailer information

#### Sales conditions



All prices are valid as of January 1st 2020 until a new version of the KERN catalogue is released. In Europe, all prices do not include the applicable V.A.T.

At KERN there is no minimum order value. For orders less than € 15.00 there is no re-sale discount available.

**Delivery Conditions**: we supply ex works Balingen, i.e. the transport costs are invoiced. With our contract logistics partners we always offer you shipping at a good rate. When shipping to islands or other more distant areas, or areas which are difficult to access, the transport costs may possibly be higher, please ask for details. Any goods supplied remain KERN's property until complete payment for the goods sold has been



Delivery is usually via courier service.



When you see this symbol by truck, please ask for prices.

#### Extract from general terms and conditions: :

Court of jurisdiction/Legal domicile: 72336 Balingen, Germany; Commercial register N°: HRB 400865, AG Stuttgart; Managing director: Albert Sauter. For the full Terms and Conditions, please refer to the website. www.kern-sohn.com/en/kern/agbs.html

Price changes and product changes are likely in individual cases due to product modifications as well as error.

Sale or return: within 14 days of purchase. Not valid for custom processes, such as, for example, special orders.

#### Warrantv:

- 3 years for products with list price ≥ € 400.00
- 2 years for products with list price < € 400.00</li>
- · Does not apply to consumables such as batteries, rechargeable battery packs, illumination units, etc.




KERN DirectCash: The quick, secure COD procedure for protection against non-payment. With the KERN DirectCash COD system, you can safely deliver orders to end customers with unknown credit rating, with no risk of non-payment. Please request further details on this procedure.

Services

Financing is available using KERN hire purchase - easy and convenient. Hire Purchase gives you the option of purchasing any product from the range against a simple monthly installment. The product value is financed over the period of the agreement. On payment of the last installment, the ownership of the contract item automatically transfers from the contractor to the contractee. The Hire Purchase Agreement can - if you so choose – be set for a period of between one and five years. This package includes the transfer of items as well as the guarantee for the entire transfer period. Compared with buying the product, KERN hire

purchase offers the advantage that the initial financial investment is largely unnecessary. This is particularly relevant when purchasing a number of products, for example when refitting a laboratory, a company department or a hospital ward. In addition the monthly installments constitute a direct cost and the item does not have to be capitalised by the purchaser.

#### After-Sales-Service

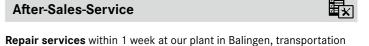

costs are additional. In urgent cases, if requested by customer, a replacement device will be provided until all repairs are completed;

Price reduction on a new device: if repair costs exceed the current value of the defective device, a new device will be offered at a discount price. This offer is valid only up to 2 years after warranty expiration.

**Spare parts service** within 48 hours, transportation costs are additional.

#### Marketing support



#### KERN catalogues, brochures, branch prospectuses your own personalised marketing tools

This catalogue and branch prospectuses are available free of charge. A neutral version of the catalogue, without the KERN address imprint, is also available for your marketing activities free of charge.

On demand, your company logo and address can be printed on the back of the KERN catalogue free of charge. In this way you have your own individual marketing tool. The KERN catalogue branch prospectuses are available in the following languages: DE, GB, FR, IT, ES.

#### Important notice



Our models are not suitable for rooms with a high level of air humidity (condensing). Please observe the applicable electrical regulations.

#### Miscellaneous

Product pictures printed in catalogue

All product pictures contained in our catalogue show devices similar to our products. Please note that possible technical innovations might be the cause of such deviations.

#### Accessory for optical instruments

Further extensive accessories for our optical instruments you can also find in the KERN Online shop on www.kern-sohn.com

# **Product group index 2020**



## Microscopes

Compound microscopes	08	01
Metallurgical microscopes	33	02
Polarising microscopes	40	03
Stereomicroscopes	45	04
Digital microscope sets	72	05
Stereo microscope sets	77	06
Stereo microscope modular system	80	07
External light sources for stereomicroscopes	89	80
Microscope cameras & Software	91	09



#### Refractometers

Retractometers		
Analogue refractometers – type: hand-held	100	10
Digital refractometers – type: hand-held	106	11
Abbe refractometers – type: desktop	112	12

# **NEWS**

# Innovative products in familiar KERN quality:

#### KERN microscope cameras and software

The new generation of KERN C-mount cameras. You now have an even larger selection of professional C-mount cameras for your documentation









Further details ▶ page 92/93 from € 1060,-

#### **School microscopes KERN OBT-1**

The new, high-quality compound microscope for use in school education

Further details ➤ page 11/12

from € 195,-

#### Basic laboratory microscope KERN OBE-12 OBE-13

Microscopes with modern design, which ensure even greater ergonomics, have now been added to the topseller Kern OBE range



Further details ➤ page 16/17 from € 425,-

#### Professional darkfield condensor

The darkfield method, thanks to the cardioid system, can now also be used for 1000× magnification



## Spring loaded arm for stereomicroscopes

Through additional coarse focussing, KERN spring loaded arms give even greater flexibility of use



Further details ➤ page 87 only € 420,-

Further details ► page 27, 29 only € 580,-

# HIGHLIGHTS 2020

#### **KERN Calibration service**

Your partner for calibration services, management of test equipment and support



Further details ► page 114 only € 99,-

#### **KERN Microscope VIS software**

The perfect software for measuring, counting and documenting your samples



#### Further details ▶ page 97

Included with every KERN camera delivery

#### KERN Checklist for "My ideal microscope"

The clever helper when you are looking for the right microscope for your application



Further details ▶ page 116

#### **Digital ECO refractometers series KERN ORF-E**

Digital refractometry through additional Eco variants is now more affordable



Further details ➤ page 106 only € 290,-

#### **KERN Tablet Camera ODC-2**

Digital microscopy brought up to date – tablet with integrated camera for optimal observation and digital documentation of the sample



Further details ► page 96 only € 890,-



# Microscopes



1	Compound microscopes	08
	Compound, Phase contrast, Digital,	
	Fluorescence and Inverted microscopes	
2	Metallurgical microscopes	33
3	Polarising microscopes	40
Ū	r oldriging interescopes	
4	Stereomicroscopes	45
	Stereo, Stereo-Zoom, Coaxial and	
	Gem microscopes	
5	Digital misusages acts	72
3	Digital microscope sets	12
6	Stereo microscope sets	77
	· ·	
7	Stereo microscope modular system	80
8	External light sources for stereomicroscopes	89
	Ring illumination and cold light sources	
9	Microscope cameras & Software	91
-		, ,




**Daniel Junger** Manager KERN Optics

Tel. +49 7433 9933-155 Fax +49 7433 9933-29195 daniel.junger@kern-sohn.com

# Compound microscopes Compound, Phase contrast, Digital, Fluorescence and Inverted microscopes



Objectives OBS

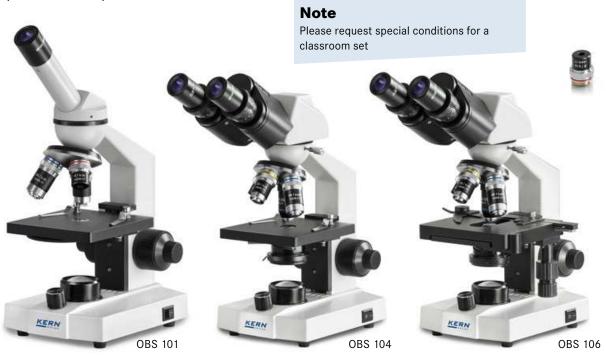

#### **EDUCATIONAL LINE**

### The school microscope – For the first steps in microscopy and for use in biology lessons

#### **Features**

- The KERN OBS range is a solid and simple school microscope range, which is easy to use due to its intuitive control elements
- The continuously dimmable 0.5W LED guarantees optimum illumination of the samples and also ensures long service life. Mobile use is also no problem through the use of rechargeable batteries
- The simple 0.65 condenser on the OBS 101 (condenser disc) and the OBS 102 (fixed condenser) ensures the very best concentration of light and illumination of the sample. The OBS 103, 104, 105 and 106 models have a 1.25 Abbe condenser which

is height-adjustable and can therefore be focussed and has an aperture diaphragm, which ensures the very best concentration of light

- · To focus the object, all models have a coarse and fine focusing knob on both sides. The mechanical stage enables you to work with the samples and move them rapidly (only for OBS 105, 106)
- · A large selection of different eyepieces and objectives is also available
- · Please find detailed information in the following model outfit list

#### Scope of application

· Primary school, secondary school, training, hobby use

#### Applications/Samples

• Translucent, thin, high-contrast, less complex samples (e.g. plant tissue, coloured cells/ parasites)

#### **Technical data**

- Finite optical system (DIN)
- Triple (OBS 101, 102) or quadplex (OBS 103, 104, 105, 106) nosepiece
- Tube 45° (OBS 101, 102, 103, 105) or 30° (OBS 104, 106) inclined/360° rotatable
- · Diopter adjustment: Both-sided (for binocular models)
- · Overall dimensions W×D×H 130×300×310 mm
- · Net weight approx. 3 kg

#### STANDARD



















Model	Standard configuration								
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	Stage	ex works €		
OBS 101	Monocular	WF 10×/Ø 18 mm	Achromatic		0,5W LED (transmitted) (battery incl., rechargeable)	fix	160,- 🕛		
OBS 102 W	Monocular	WF 10×/Ø 18 mm	Achromatic		0,5W LED (transmitted) (battery incl., rechargeable)	fix	170,-		
OBS 103 W	Monocular	WF 10×/Ø 18 mm	Achromatic	4×/10×/40×	0,5W LED (transmitted) (battery incl., rechargeable)	fix	175,-		
OBS 104	Binocular	WF 10×/Ø 18 mm	Achromatic	4^/10^/40^	0,5W LED (transmitted) (battery incl., rechargeable)	fix	245,- 🕛		
OBS 105 🔤	Monocular	WF 10×/Ø 18 mm	Achromatic		0,5W LED (transmitted) (battery incl., rechargeable)	mechanical	210,-		
OBS 106	Binocular	WF 10×/ø 18 mm	Achromatic		0,5W LED (transmitted) (battery incl., rechargeable) mechanical		280,- 🕛		

New model

Price reduction

Model outfit			Мо	del KE	RN			Order number	Price/piece excl. of VAT
		OBS 101	OBS 102	OBS 103	OBS 104	OBS 105	OBS 106		ex works €
	WF 10×/ø 18 mm	1	✓	✓	44	✓	11	OBB-A1473	30,-
Eyepieces	WF 16×/Ø 13 mm	0	0	0	00	0	00	OBB-A1474	30,-
(23,2 mm)	WF 20×/ø 11 mm	0	0	0	00	0	00	OBB-A1475	30,-
	WF 10×/Ø 18 mm (with Pointer)	0	0	0	0	0	0	OBB-A1561	30,-
	4×/0,10 W.D. 18,0 mm	1	✓	✓	✓	✓	✓	OBB-A1476	35,-
	10×/0,25 W.D. 7,0 mm	✓	✓	1	✓	✓	✓	OBB-A1477	40,-
Achromatic objectives	40×/0,65 (spring-loaded) W.D. 0,53 mm	1	✓	✓	✓	✓	1	OBB-A1478	45,-
02,000	60×/0,85 (spring-loaded) W.D. 0,1 mm	0	0	0	0	0	0	OBB-A1479	60,-
	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	0	0	0	0	0	0	OBB-A1480	65,-
	4×/0,10 W.D. 14,5 mm	0	0	0	0	0	0	OBB-A1562	40,-
	10×/0,25 W.D. 5,65 mm		0	0	0	0	0	OBB-A1563	55,-
E-Plan objectives	40×/0,65 (spring-loaded) W.D. 0,85 mm	0	0	0	0	0	0	OBB-A1564	80,-
	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	0	0	0	0	0	0	OBB-A1565	85,-
	100×/0,80 (dry) (spring-loaded) W.D. 0,15 mm	0	0	0	0	0	0	OBB-A1442	170,-
	Plan 100×/1,0 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	0	0	OBB-A1441	175,-
Monocular tube	45° inclined/360° rotatable	<b>✓</b>	✓	✓		✓		OBB-A1471	
Binocular tube	30° inclined/360° rotatable     Interpupillary distance 55-75 mm     Diopter adjustment: Both-sided				<b>✓</b>		<b>4</b>	OBB-A1472	
Fixed stage	Stage size W×D 110×120 mm     Coaxial coarse and fine focusing knobs, scale: 2,5 µm	~	~	~	~				
Mechanical stage	Stage size W×D 115×125 mm Travel 75×18 mm Coaxial coarse and fine focusing knobs, scale: 2,5 µm					<b>✓</b>	<b>~</b>		
C d	Simple condenser N.A. 0,65	✓	✓						
Condenser	Abbe N.A. 1,25 (aperture diaphragm)			✓	✓	✓	✓		
Illumination	0,5 W LED illumination system (transmitted) (rechargeable)	<b>✓</b>	<b>✓</b>	✓	✓	<b>✓</b>	1		
	Blue	✓	✓	✓	✓	✓	✓	OBB-A1466	19,-
Colour filters	Green	0	0	0	0	0	0	OBB-A1467	19,-
for transmitted illumination	Yellow	0	0	0	0	0	0	OBB-A1468	19,-
	Grey	0	0	0	0	0	0	OBB-A1184	19,-

<sup>✓ =</sup> Included with delivery

O = Option





#### Note

Please request special conditions for a classroom set



Monocular version







Objectives OBT

#### **EDUCATIONAL LINE**

### The modern compound microscope for teaching in your class room

#### **Features**

- The KERN OBT range is a high-quality school microscope, which will impress you with its intuitive control elements, sturdy construction and modern design
- The infinitely dimmable 1W LED guarantees optimum illumination of the samples and also ensures long service life. Mobile use is also no problem through optional battery operation
- The simple 0.65 condenser lens with adjustable aperture diaphragm on the OBT 101 ensures the very best concentration of light and illumination of the sample. The OBT 102, 103, 104, 105, 106 models have a 1.25 Abbe condenser which is height-adjustable and can therefore be focussed and has an aperture diaphragm, which ensures the very best concentration of light
- To focus the object accurately, all models have a coarse and fine focusing knob on both sides. The mechanical angle table enables you to work with the samples and move them rapidly (for OBT 103, 104, 105, 106 models)
- A large selection of different eyepieces and objectives is also available
- · A dust cover as well as user instructions are included with the delivery
- · Please find detailed information in the following model outfit list

#### Scope of application

• Primary school, secondary school, training, hobby use

#### Applications/Samples

• Translucent, thin, high-contrast, less complex samples (e.g. plant tissue, coloured cells/ parasites)

#### **Technical data**

- · Finite optical system (DIN)
- Triple (OBT 101) or quadplex (OBT 102, 103, 104, 105, 106) nosepiece
- Tube 45° inclined/360° rotatable
- · Diopter adjustment: Both-sided (for binocular models)
- · Overall dimensions W×D×H 195×147×325 mm
- · Net weight approx. 2,5 kg





STANDARI	)						OPTION	
Ø	0	00	7	Ð				uuu
360°	MONO	BINO	ABBE	LED	230 V	1 DAY	BATT	SCALE
			not					

	OBI 101								
Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	Stage	ex works €		
OBT 101	Monocular	HWF 10×/Ø 18 mm	Achromatic		1W LED (transmitted)	fix	195,-		
OBT 102	Monocular	HWF 10×/Ø 18 mm	Achromatic		1W LED (transmitted)	fix	215,-		
OBT 103	Monocular	HWF 10×/Ø 18 mm	Achromatic	4×/10×/40×	1W LED (transmitted)	mechanical	245,-		
OBT 104	Binocular	HWF 10×/Ø 18 mm	Achromatic	4^/10^/40^	1W LED (transmitted)	mechanical	325,-		
OBT 105	Monocular	HWF 10×/Ø 18 mm	Achromatic		1W LED (transmitted)	mechanical	270,-		
OBT 106	Binocular	HWF 10×/Ø 18 mm	Achromatic		1W LED (transmitted)	mechanical	345,-		

Model outfit				Model	KERN		Order number	Price/piece excl. of VAT	
		OBT 101	OBT 102	OBT 103	OBT 104	OBT 105	OBT 106		ex works €
	WF 10×/ø 18 mm	✓	✓	1	44	✓	11	OBB-A3200	29,-
Eyepieces (23,2 mm)	WF 10×/ø 18 mm (with Pointer)	0	0	0	0	0	0	OBB-A3201	29,-
(20,2)	WF 10×/ø 18 mm (reticule 0,1 mm)	0	0	0	0	0	0	OBB-A3202	29,-
	4×/0,10 W.D. 27 mm	✓	✓	✓	✓	✓	✓	OBB-A3203	25,-
	10×/0,25 W.D. 7 mm	<b>✓</b>	✓	✓	✓	✓	✓	OBB-A3204	29,-
Achromatic objectives	40×/0,65 (spring-loaded) W.D. 0,6 mm		1	1	✓	✓	✓	OBB-A3205	39,-
	100×/1,25 (oil) (spring-loaded) W.D. 0,2 mm		0	0	0	✓	✓	OBB-A3206	49,-
	60×/0,85 (spring-loaded) W.D. 0,4 mm		0	0	0	0	0	OBB-A3207	59,-
Monocular tube	45° inclined/360° rotatable	<b>✓</b>	✓	1	0	✓	0	OBB-A3221	
Binocular tube	Siedentopf 45° inclined/360° rotatable     Interpupillary distance 48-75 mm     Diopter adjustment: One-sided	0	0	0	<b>✓</b>	0	<b>✓</b>	OBB-A3222	
Fixed stage	Stage size W×D 115×110 mm     Coaxial coarse and fine focusing knobs, scale: 2 μm	<b>✓</b>	~						
Mechanical stage	Stage size W×D 115×110 mm     Travel 52×20 mm     Coaxial coarse and fine focusing knobs, scale: 2 μm     One slide holder			<b>✓</b>	<b>√</b>	<b>✓</b>	~		
	Simple condenser N.A. 0,65	<b>✓</b>							
Condenser	Abbe N.A. 1,25 (aperture diaphragm)		✓	<b>✓</b>	✓	✓	<b>✓</b>		
Illumination	1 W LED spare bulb (transmitted)	✓	~	~	<b>~</b>	✓	<b>✓</b>	OBB-A3208	19,-
	Blue	0	0	0	0	0	0	OBB-A3212	19,-
Colour filters	Green	0	0	0	0	0	0	OBB-A3210	19,-
for transmitted illumination	Yellow	0	0	0	0	0	0	OBB-A3211	19,-
	Grey	0	0	0	0	0	0	OBB-A3209	19,-





Trinocular version

#### Note

Please request special conditions for a classroom set



Monocular version









Objectives OBE





Simple polarising unit

#### **EDUCATIONAL LINE**

### The fully equipped all-round compound microscope for school, training and laboratories

#### **Features**

- The KERN OBE series is a range of high-quality, fully-equipped compound microscopes, which can't be beaten in terms of ease of use and ergonomic design
- The strong and continuously dimmable 3 W LED guarantees optimum illumination of the samples and also ensures long service life. Mobile use of several models is also no problem through the use of rechargeable batteries
- · The height-adjustable and thereby focusable 1,25 Abbe condenser with aperture diaphragm is a further quality feature of the OBE series and ensures the very best concentration of light
- · Height adjustment of the fully-equipped mechnical stage is carried out using a coarse and fine focusing knob on both

- sides. The ergonomically designed coaxial drive enables you to work with the samples and move them rapidly
- A large selection of different eyepieces and objectives, a simple polarising unit and a darkfield kit are available to you as
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

• Training, haematology, sediment investigation, doctor's practise

#### Applications/Samples

• Translucent, thin, high-contrast, less complex samples (e.g. plant tissue, coloured cells/ parasites)

#### **Technical data**

- · Finite optical system
- · Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: One-sided (for binocular and trinocular models)
- · Overall dimensions W×D×H 320×180×365 mm
- · Net weight approx. 5,5 kg

#### STANDARD

























9 MONO	OO BINO	TRINO	ABBE	LED	RECHARGE OBE 103, OBE 113	230 V	1 DAY	DF	POLAR	SCA
/lodel								Stand	ard co	nfigu
	T .			Б			01: : 1			

		OBE 113				
Model			Standar	d configuration		Price excl. of VAT
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €
OBE 101	Monocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	370,- 🕛
OBE 102	Binocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	480,-
OBE 103	Binocular	HWF 10×/Ø 18 mm	Achromatic	4×/10×/40×	3 W LED (transmitted) (battery incl., rechargeable)	550,-
OBE 104	Trinocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	590,-
OBE 111	Monocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	450,-
OBE 112	Binocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	540,-
OBE 113	Binocular	HWF 10×/Ø 18 mm	Achromatic	4×/10×/40×/100×	3 W LED (transmitted) (battery incl., rechargeable)	610,-
OBE 114	Trinocular	HWF 10×/Ø 18 mm	Achromatic		3 W LED (transmitted)	650,-

Model outfit		Model	KERN		Order number	Price/piece excl. of VAT	
		OBE 101	OBE 102	OBE 103	OBE 104		ex works €
	HWF 10×/ø 18 mm	✓	11	11	11	OBB-A1403	45,-
Eyepieces	WF 16×/ø 13 mm	0	00	00	00	OBB-A1354	45,-
(23,2 mm)	HWF 10×/Ø 18 mm (with Pointer)	0	0	0	0	OBB-A1348	40,-
	HWF 10×/ø 18 mm (reticule 0,1 mm) (non-adjustable)	0	0	0	0	OBB-A1349	60,-
	4×/0,10 W.D. 18,6 mm	✓	✓	✓	✓	OBB-A1111	35,-
	10×/0,25 W.D. 6,5 mm	✓	✓	✓	✓	OBB-A1108	45,-
	40×/0,65 (spring-loaded) W.D. 0,47 mm	✓	✓	✓	✓	OBB-A1112	70,-
Achromatic	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	0	0	0	0	OBB-A1109	105,-
objectives	20×/0,40 (spring-loaded) W.D. 1,75 mm	0	0	0	0	OBB-A1110	85,-
	60×/0,85 (spring-loaded) W.D. 0,1 mm	0	0	0	0	OBB-A1113	105,-
	E-Plan 100×/0,80 (dry) (spring-loaded) W.D. 0,15 mm	0	0	0	0	OBB-A1442	170,-
	Plan 100×/1,0 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	OBB-A1441	175,-
Monocular tube	30° inclined/360° rotatable	✓				OBB-A1227	
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Diopter adjustment: One-sided		~	~		OBB-A1123	
Trinocular tube	see binocular tube     Light distribution 20:80				✓	OBB-A1341	
Mechanical stage	<ul> <li>Stage size W×D 125×115 mm</li> <li>Travel 50×70 mm</li> <li>Coaxial coarse and fine focusing knobs, scale: 2 μm</li> </ul>	<b>✓</b>	<b>✓</b>	<b>✓</b>	<b>√</b>		
Condenser	Abbe N.A. 1,25 (aperture diaphragm)	✓	✓	✓	✓	OBB-A1101	85,-
Darkfield unit	Usable for 4× – 40× objectives	0	0	0	0	OBB-A1148	85,-
Polarising unit	Analyser/Polariser	0	0	0	0	OBB-A1276	265,-
Illumination	3 W LED illumination system (transmitted) (non-rechargeable)	✓	✓		✓		
Indinination	3 W LED illumination system (transmitted) (rechargeable)			✓			
	Blue	0	0	0	0	OBB-A1466	19,-
Colour filters for transmitted	Green	0	0	0	0	OBB-A1467	19,-
illumination	Yellow	0	0	0	0	OBB-A1468	19,-
	Grey	0	0	0	0	OBB-A1184	19,-
C-Mount	0,5× (focus adjustable)				0	OBB-A1137	160,-
G-IVIOUTIL	1×				0	OBB-A1139	100,-

Model outfit			Mode	KERN		Order number	Price/piece
		OBE 111	OBE 112	OBE 113	OBE 114	-	ex works €
	HWF 10×/ø 18 mm	✓	11	11	44	OBB-A1403	45,-
Eyepieces	WF 16×/Ø 13 mm	0	00	00	00	OBB-A1354	45,-
(23,2 mm)	HWF 10×/ø 18 mm (with Pointer)	0	0	0	0	OBB-A1348	40,-
	HWF 10×/Ø 18 mm (reticule 0,1 mm) (non-adjustable)	0	0	0	0	OBB-A1349	60,-
	4×/0,10 W.D. 18,6 mm	✓	1	✓	✓	OBB-A1111	35,-
	10×/0,25 W.D. 6,5 mm	✓	✓	✓	✓	OBB-A1108	45,-
	40×/0,65 (spring-loaded) W.D. 0,47 mm	✓	✓	✓	✓	OBB-A1112	70,-
Achromatic	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	✓	✓	✓	✓	OBB-A1109	105,-
objectives	20×/0,40 (spring-loaded) W.D. 1,75 mm	0	0	0	0	OBB-A1110	85,-
	60×/0,85 (spring-loaded) W.D. 0,1 mm	0	0	0	0	OBB-A1113	105,-
	E-Plan 100×/0,80 (dry) (spring-loaded) W.D. 0,15 mm	0	0	0	0	OBB-A1442	170,-
	Plan 100×/1,0 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	OBB-A1441	175,-
Monocular tube	30° inclined/360° rotatable	✓				OBB-A1227	
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Diopter adjustment: One-sided		<b>✓</b>	~		OBB-A1123	
Trinocular tube	see binocular tube     Light distribution 20:80				✓	OBB-A1341	
Mechanical stage	Stage size W×D 125×115 mm     Travel 50×70 mm     Coaxial coarse and fine focusing knobs, scale: 2 μm	<b>✓</b>	~	~	~		
Condenser	Abbe N.A. 1,25 (aperture diaphragm)	✓	✓	✓	✓	OBB-A1101	85,-
Darkfield unit	Usable for 4× - 40× objectives	0	0	0	0	OBB-A1148	85,-
Polarising unit	Analyser/Polariser	0	0	0	0	OBB-A1276	265,-
III	3 W LED illumination system (transmitted) (non-rechargeable)	~	~		~		
Illumination	3 W LED illumination system (transmitted) (rechargeable)			✓			
	Blue	0	0	0	0	OBB-A1466	19,-
<b>Colour filters</b> for transmitted illumination	Green	0	0	0	0	OBB-A1467	19,-
	Yellow	0	0	0	0	OBB-A1468	19,-
	Grey	0	0	0	0	OBB-A1184	19,-
C.Mt	0,5× (focus adjustable)				0	OBB-A1137	160,-
C-Mount	1×				0	OBB-A1139	100,-





#### Note

Please request special conditions for a classroom set



Monocular version



Trinocular version



Butterfly tube

#### **EDUCATIONAL LINE**

# Elegant, dynamic and impressive – the new all-round compound microscope for schools, training and laboratories

#### **Features**

- The brand new OBE-12/13 range stands out through its exclusive, dynamic device, which is second to none in terms of sturdy construction and ergonomics. The clever storage compartment on the back will enables quick practical storage for your power cable. Thanks to the USB connection technology, it is also possible to supply power using an external powerbank
- The impressive, infinitely dimmable 3 W LED guarantees bright illumination of your sample
- A further highlight is the Butterfly tube
  which is integrated as standard and which
  enables you to achieve the ideal viewing
  angle. The height-adjustable and thereby
  focusable 1.25 Abbe condenser with aperture
  diaphragm is a further quality feature of the
  OBE range and guarantees the very best
  concentration of light
- Height adjustment of the fully-equipped mechnical stage is carried out using a coarse and fine focusing knob on both sides. The ergonomically designed coaxial drive enables you to work with the samples and move them rapidly
- A large selection of different eyepieces and objectives are available to you as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit liet
- Please find detailed information in the following model outfit list

#### Scope of application

 Training, haematology, sediment investigation, doctor's practise

#### Applications/Samples

 Translucent, thin, high-contrast, less complex samples (e.g. plant tissue, coloured cells/ parasites)

#### **Technical data**

- Finite optical system
- · Quadplex nosepiece
- $\bullet$  Butterfly 30° inclined
- Diopter adjustment: One-sided (for binocular and trinocular models)
- Overall dimensions W×D×H 360×150×320 mm
- Net weight approx. 4,6 kg

#### STANDARD

















Model		Standard configuration								
KERN	Tube	Illumination	ex works €							
OBE 121	Monocular	HWF 10×/Ø 18 mm	Achromatic		3W LED (transmitted)	425,-				
OBE 122	Binocular	HWF 10×/Ø 18 mm	Achromatic	4×/10×/40×	3W LED (transmitted)	520,-				
OBE 124	Trinocular	HWF 10×/Ø 18 mm	Achromatic		3W LED (transmitted)	630,-				
OBE 131	Monocular	HWF 10×/Ø 18 mm	Achromatic		3W LED (transmitted)	490,-				
OBE 132	Binocular	HWF 10×/Ø 18 mm	Achromatic	4×/10×/40×/100×	3W LED (transmitted)	580,-				
OBE 134	Trinocular	HWF 10×/Ø 18 mm	Achromatic		3W LED (transmitted)	690,-				

Model outfit				Model	KERN	ı		Order number	Price/piece excl. of VAT	
		OBE 121	OBE 122	OBE 124	OBE 131	OBE 132	OBE 134		ex works €	
	HWF 10×/ø 18 mm	✓	11	11	✓	11	11	OBB-A1403	45,-	
Eyepieces	WF 16×/ø 13 mm	0	00	00	0	00	00	OBB-A1354	45,-	
(23,2 mm)	HWF 10×/ø 18 mm (with pointer)	0	0	0	0	0	0	OBB-A1348	40,-	
	HWF 10×/Ø 18 mm (reticule 0,1 mm) (non-adjustable)	0	0	0	0	0	0	OBB-A1349	60,-	
	4×/0,10 W.D. 18,6 mm	✓	✓	1	1	1	✓	OBB-A1111	35,-	
	10×/0,25 W.D. 6,5 mm	✓	1	<b>✓</b>	1	✓	✓	OBB-A1108	45,-	
	40×/0,65 (spring-loaded) W.D. 0,47 mm	✓	✓	1	✓	1	✓	OBB-A1112	70,-	
Achromatic	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	0	0	0	✓	1	1	OBB-A1109	105,-	
objectives	20×/0,40 (spring-loaded) W.D. 1,75 mm	0	0	0	0	0	0	OBB-A1110	85,-	
	60×/0,85 (spring-loaded) W.D. 0,1 mm	0	0	0	0	0	0	OBB-A1113	105,-	
	E-Plan 100×/0,80 (dry) (spring-loaded) W.D. 0,15 mm	0	0	0	0	0	0	OBB-A1442	170,-	
	Plan 100×/1,0 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	0	0	OBB-A1441	175,-	
Monocular tube	30° inclined	<b>✓</b>			✓					
Binocular tube	Butterfly 30° inclined     Interpupillary distance 48 – 75 mm     Diopter adjustment: One-sided		<b>*</b>			1				
Trinocular tube	see binocular tube     Light distribution 20:80			~			1			
Mechanical stage	<ul> <li>Stage size W×D 125×115 mm</li> <li>Travel 50×70 mm</li> <li>Coaxial coarse and fine focusing knobs, scale: 2 μm</li> </ul>	~	~	~	~	✓	<b>✓</b>			
Condenser	Abbe N.A. 1,25 (aperture diaphragm)	✓	✓	✓	✓	✓	✓	OBB-A1101	85,-	
Darkfield unit	Usable for 4× - 40× objectives	0	0	0	0	0	0	OBB-A1148	85,-	
Illumination	3 W LED illumination system (transmitted)	✓	✓	1	✓	✓	✓			
	Blue	0	0	0	0	0	0	OBB-A1466	19,-	
Colour filters	Green	0	0	0	0	0	0	OBB-A1467	19,-	
for transmitted Illumination	Yellow	0	0	0	0	0	0	OBB-A1468	19,-	
	Grey	0	0	0	0	0	0	OBB-A1184	19,-	
	0,5× (focus adjustable)			0			0	OBB-A1137	160,-	
C-Mount	1×			0			0	OBB-A1139	100,-	





Trinocular version



Simple polarising attachment

#### LAB LINE

### The high-performance compound microscope for every laboratory with fixed, pre-centred Koehler illumination

#### **Features**

- The KERN OBF models are excellent, stable laboratory microscopes for all common routine applications. A central feature of this adaptable, robust microscope series is the stable mechanism which can be adjusted precisely
- · Depending on the application, there is a choice of models with strong, continuously dimmable 3 W LED or 20 W halogen illumination (Philips)
- The fixed, pre-centred and focusable 1,25 Abbe condenser with aperture diaphragm and field diaphragm gives you a simplified Koehler illumination, without having to move the centre
- The large mechanical stage and its specimen holder holds up to two samples at the same time and is quick and easy to focus using a coaxial coarse and fine focusing knob on both sides

- A large selection of eyepieces, objectives and colour filters as well as a darkfield condenser and a simple polarising unit are available to you as accessories
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model
- · Please find detailed information in the following model outfit list

#### Scope of application

· Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, oncology, entomology, vets, water analysis and breweries

#### Applications/Samples

· Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue)

#### **Technical data**

- · Finite optical system (DIN)
- · Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: One-sided
- · Overall dimensions W×D×H 395×200×380 mm
- · Net weight approx. 6,7 kg







								0		
360°	BINO	TRINO	ABBE	HAL	LED	230 V	1 DAY	DF	POLAR	SCAL
Mod	el								Stand	ard c

Model		Standard configuration									
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €					
OBF 121	Binocular	HWF 10×/Ø 18 mm	Achromatic		20 W Halogen (transmitted)	750,-					
OBF 122	Binocular	HWF 10×/Ø 18 mm	Plan		20 W Halogen (transmitted)	860,-					
OBF 123	Binocular	HWF 10×/Ø 18 mm	Plan	4×/10×/40×/100×	3 W LED (transmitted)	840,-					
OBF 131	Trinocular	HWF 10×/Ø 18 mm	Achromatic	4^/10^/40^/100^	20 W Halogen (transmitted)	860,-					
OBF 132	Trinocular	HWF 10×/Ø 18 mm	Plan		20 W Halogen (transmitted)	970,-					
OBF 133	Trinocular	HWF 10×/Ø 18 mm	Plan		3 W LED (transmitted)	950,-					

Model outfit				Model	KERN			Order number	Price/piece excl. of VAT	
		OBF 121	OBF 131	OBF 122	OBF 132	OBF 123	OBF 133		ex ci. of VAI ex works €	
	HWF 10×/ø 18 mm	11	11	11	11	11	11	OBB-A1403	45,-	
Eyepieces	WF 16×/ø 13 mm	00	00	00	00	00	00	OBB-A1354	45,-	
(23,2 mm)	HWF 10×/Ø 18 mm (with Pointer)	0	0	0	0	0	0	OBB-A1348	40,-	
	HWF 10×/Ø 18 mm (reticule 0,1 mm) (non-adjustable)	0	0	0	0	0	0	OBB-A1349	60,-	
	4×/0,10 W.D. 18,6 mm	<b>√</b>	✓					OBB-A1111	35,-	
	10×/0,25 W.D. 6,5 mm	1	1					OBB-A1108	45,-	
Achromatic	40×/0,65 (spring-loaded) W.D. 0,47 mm	1	1					OBB-A1112	70,-	
objectives	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm	1	1					OBB-A1109	105,-	
	20×/0,40 (spring-loaded) W.D. 1,75 mm	0	0					OBB-A1110	85,-	
	60×/0,85 (spring-loaded) W.D. 0,1 mm	0	0					OBB-A1113	105,-	
	4×/0,10 W.D. 14,5 mm			1	<b>✓</b>	✓	1	OBB-A1255	65,-	
	10×/0,25 W.D. 5,65 mm			1	<b>✓</b>	<b>√</b>	1	OBB-A1238	80,-	
	40×/0,65 (spring-loaded) W.D. 0,85 mm			1	1	1	1	OBB-A1256	115,-	
Plan objectives	100×/1,25 (oil) (spring-loaded) W.D. 0,07 mm			1	1	1	1	OBB-A 1239	150,-	
	20×/0,40 (spring-loaded) W.D. 1,5 mm			0	0	0	0	OBB-A 1249	115,-	
	60×/0,85 (spring-loaded) W.D. 0,07 mm			0	0	0	0	OBB-A1269	155,-	
	100×/1,0 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	0	0	OBB-A1441	175,-	
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     (for non-infinity system)     Diopter adjustment: One-sided		0	~	0	<b>✓</b>	0	OBB-A1129	,	
Trinocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 20:80 (for non-infinity system)     Diopter adjustment: One-sided	0	<b>✓</b>	0	~	0	~	OBB-A1345		
Mechanical stage	Stage size W×D 145×130 mm Travel 76×52 mm Coaxial coarse and fine focusing knobs, scale: 2 µm Two slide holder	~	<b>✓</b>	~	~	<b>✓</b>	<b>✓</b>			
Condenser	Abbe N.A. 1,25 precentered (aperture diaphragm)	1	✓	1	<b>✓</b>	✓	1	OBB-A1103	85,-	
Darkfield condenser	N.A. 0,85 - 0,91 (dry, paraboloid)	0	0	0	0	0	0	OBB-A1422	230,-	
	20 W Halogen spare bulb (transmitted)	✓	1	✓	✓			OBB-A1370	25,-	
Illumination	3 W LED illumination system (transmitted) (non-rechargeable)					1	1			
Polarising unit	Analyser/Polariser	0	0	0	0	0	0	OBB-A1277	265,-	
	Blue (built-in)	✓	✓	✓	✓	✓	✓			
Colour filters	Green	0	0	0	0	0	0	OBB-A1188	19,-	
for transmitted illumination	Yellow	0	0	0	0	0	0	OBB-A1165	19,-	
	Grey	0	0	0	0	0	0	OBB-A1183	19,-	
	0,47× (focus adjustable)				0		0	OBB-A1135	160,-	
	0,5× (focus adjustable)		0					OBB-A1515	185,-	
	, (,									
C-Mount	1×				0		0	OBB-A1142	100,-	





Trinocular version



Simple polarising attachment

#### LAB LINE

# The flexible laboratory assistant with infinity optical system and fixed, pre-centred Koehler illumination

#### **Features**

- The OBL series stands out through its infinity optical unit and is therefore ideally suited for all demanding transmitted illumination applications. The robust and ergonomic stand base guarantees safe and comfortable working
- Depending on the application, there is a choice of models with strong, continuously dimmable 3 W LED or 20 W halogen illumination (Philips)
- The fixed, pre-centred and focusable 1,25
   Abbe condenser with aperture diaphragm and field diaphragm gives you a simplified Koehler illumination, without having to move the centre
- The large mechanical stage and its specimen holder holds up to two samples at the same time and is quick and easy to focus using a coaxial coarse and fine focusing knob on both sides

- A large selection of eyepieces, objectives and colour filters as well as a darkfield condenser, a simple polarising unit, different phase contrast kits through to HBO and LED fluorescence units are available to you as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list
- Please find detailed information in the following model outfit list

#### Scope of application

 Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, oncology, entomology, vets, water analysis and breweries

#### Applications/Samples

 Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue)

#### Technical data

- · Infinity optical system
- Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- Diopter adjustment: One-sided
- Overall dimensions W×D×H 395×200×380 mm
- Net weight approx. 6,7 kg

STANDAR	D							
Ø	00		7	Ф	Ð	$\infty$	<b>—</b>	••
360°	BINO	TRINO	ABBE	HAL	LED	INFINITY	230 V	1 DAY

Model			Standard	d configuration		<b>Price</b> excl. of VAT	
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €	
OBL 125	Binocular	HWF 10×/Ø 20 mm	Infinity E-Plan		20 W Halogen (transmitted)	910,-	
OBL 127	Binocular	HWF 10×/Ø 20 mm	Infinity E-Plan	4× /10× /40× /100×	3 W LED (transmitted)	910,-	
OBL 135	Trinocular	HWF 10×/Ø 20 mm	Infinity E-Plan	4×/10×/40×/100×	20 W Halogen (transmitted)	1020,-	
OBL 137	Trinocular	HWF 10×/Ø 20 mm	Infinity E-Plan	]	3 W LED (transmitted)	1020	

Eyepieces (23,2 mm)	HWF 10×/Ø 20 mm  WF 16×/Ø 13 mm  HWF 10×/Ø 20 mm (with Pointer)  4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	OBL 125	OBL 135 ✓✓ ○○ ✓ ✓ ✓	OBL 127 ✓✓ ○○ ✓ ✓	OBL 137 ✓✓ ○○ ○	OBB-A1404 OBB-A1354 OBB-A1448 OBB-A1161	excl. of VAT ex works €  95,-  45,-  80,-
(23,2 mm)	WF 16×/Ø 13 mm  HWF 10×/Ø 20 mm (with Pointer)  4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	**  **  **  **  **  **  **  **  **  **	** ** ** ** ** ** ** ** ** ** ** ** **	00 0	00	OBB-A1354 OBB-A1448	45,-
(23,2 mm)	HWF 10×/Ø 20 mm (with Pointer)  4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	0 *  *  *  *  *  *  *  *  *  *  *  *  *	0 4	0 1	0	OBB-A1448	
Infinity	4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	<ul><li>✓</li><li>✓</li><li>✓</li><li>✓</li></ul>	✓ ✓ ✓	✓ ✓	✓		80,-
	10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	<ul><li>✓</li><li>✓</li><li>✓</li></ul>	<b>✓</b>	✓		OBB-A1161	
	40×/0,65 (spring-loaded) W.D. 0,58 mm 100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	<ul><li>✓</li><li>✓</li><li>○</li></ul>	✓		<b>✓</b>		90,-
	100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	<b>✓</b>		<b>✓</b>	l	OBB-A1159	150,-
	Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	0	✓		✓	OBB-A1160	265,-
	Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm			✓	✓	OBB-A1158	335,-
	, , , , ,	_	0	0	0	OBB-A1250	265,-
		0	0	0	0	OBB-A1270	360,-
	Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	0	OBB-A1437	445,-
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm (for infinity system)     Diopter adjustment: One-sided	<b>√</b>	0	~	0	OBB-A1130	
Trinocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm	0	<b>✓</b>	0		OBB-A1346	
Timoculai tube	<ul> <li>Light distribution 20:80 (for infinity system)</li> <li>Diopter adjustment: One-sided</li> </ul>				~	OBB-A1549	
Mechanical stage	• Stage size W×D 145×130 mm • Travel 76×52 mm • Coaxial coarse and fine focusing knobs, scale: 2 µm • Two slide holder	✓	<b>✓</b>	<b>✓</b>	<b>✓</b>		
Condenser	Abbe N.A. 1,25 precentered (aperture diaphragm)	✓	✓	✓	1	OBB-A1103	85,-
Darkfield condenser	N.A. 0,85 – 0,91 (dry, paraboloid)	0	0	0	0	OBB-A1422	230,-
	20 W Halogen spare bulb (transmitted)	✓	✓			OBB-A1370	25,-
Illumination	3 W LED illumination system (transmitted) (non-rechargeable)			✓	✓		
Polarising unit	Analyser/Polariser	0	0	0	0	OBB-A1277	265,-
	Single unit with ∞ PH-Plan objective 10×	0	0	0	0	OBB-A1215	310,-
Phase contrast	Single unit with ∞ PH-Plan objective 20×	0	0	0	0	OBB-A1217	350,-
units (including PH-condenser and	Single unit with ∞ PH-Plan objective 40×	0	0	0	0	OBB-A1219	420,-
PH-slides)	Single unit with ∞ PH-Plan objective 100×	0	0	0	0	OBB-A1213	510,-
	If required, there are several magnification levels, please conta	act our OP	TICS-Tean	n			
FI	100 W HBO Epi Fluorescence unit, three-hole slide (B/G) including centering objective	0	0	0	0	OBB-A1154	2350,-
Fluorescence unit	3 W LED Epi Fluorescence unit, three-hole slide (B/G) including centering objective	0	0	0	0	OBB-A1157	3110,-
	Blue (built-in)	✓	✓	✓	✓		
Colour filters	Green	0	0	0	0	OBB-A1188	19,-
for transmitted illumination	Yellow	0	0	0	0	OBB-A1165	19,-
	Grey	0	0	0	0	OBB-A1183	19,-
	0,47× (focus adjustable)	0	0	0		OBB-A1135	160,-
	0,5× (focus adjustable)				0	OBB-A1515	185,-
C-Mount		0	0	0		OBB-A1142	100,-
	1×				0	OBB-A1514	120,-





Mounted phase contrast condenser



Simple PH condenser with 40× PH slide

#### LAB LINE

# High-quality phase contrast microscope – specially pre-configured with a series of options for flexible expansion

#### **Features**

- We have developed this series specially for general applications with phase contrast method. In addition, the stable, modular construction system of the OBL series offers many more options
- A strong and continously adjustable 20 W halogen illumination unit (Philips) ensures the optimum lighting conditions
- A special fixed, pre-centred phase contrast condenser with aperture diaphragm as well as field diaphragm give you a simplified Koehler illumination and thereby a powerful phase-contrast display of your sample
- The large mechanical stage and its specimen holder holds up to two samples at the same time and is quick and easy to focus using a coaxial coarse and fine focusing knob on both sides

- A large selection of eyepieces, objectives and colour filters, a simple polarising unit as well as further phase contrast units are available to you as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list
- Please find detailed information in the following model outfit list

#### Scope of application

 Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, Sewage treatment plants, Oncology, entomology, vets, water analysis and breweries

#### Applications/Samples

 Specially for extremely translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue) with phase contrast

#### **Technical data**

- Infinity optical system
- Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- Diopter adjustment: One-sided
- Overall dimensions W×D×H 395×200×380 mm
- Net weight approx. 6,7 kg

Model		Standard configuration								
KERN	Tube	Eyepiece	Objective quality	Objectives	ex works €					
OBL 145	Binocular	HWF 10×/Ø 20 mm	Infinity E-Plan/ Plan	4×/PH10×/	20 W Halogen (transmitted)	1420,-				
OBL 155	Trinocular	HWF 10×/Ø 20 mm	Infinity E-Plan/ Plan	PH40×/100×	20 W Halogen (transmitted)	1530,-				

Model outfit		Mode	el KERN	Order number	Price/piece excl. of VAT
		OBL 145	OBL 155		ex works €
	HWF 10×/ø 20 mm	11	44	OBB-A1404	95,-
Eyepieces (23,2 mm)	WF 16×/Ø 13 mm	00	00	OBB-A1354	45,-
(==,=,	HWF 10×/ø 20 mm (with Pointer)	0	0	OBB-A1448	80,-
	4×/0,10 W.D. 12,1 mm	✓	✓	OBB-A1161	90,-
	10×/0,25 W.D. 2,1 mm	0	0	OBB-A1159	150,-
	40×/0,65 (spring-loaded) W.D. 0,58 mm	0	0	OBB-A1160	265,-
Infinity E-Plan objectives	100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm	✓	✓	OBB-A1158	335,-
E i iun objectives	Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm	0	0	OBB-A1250	265,-
	Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	0	0	OBB-A1270	360,-
	Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm	0	0	OBB-A1437	445,-
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm (for infinity system)     Diopter adjustment: One-sided	~	0	OBB-A1130	
Trinocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 20:80 (for infinity system)     Diopter adjustment: One-sided	0	<b>✓</b>	OBB-A1549	
Mechanical stage	<ul> <li>Stage size W×D 145×130 mm</li> <li>Travel 76×52 mm</li> <li>Coaxial coarse and fine focusing knobs, scale: 2 µm</li> <li>Two slide holder</li> </ul>		*		
PH condenser	Abbe N.A. 1,25 precentered, for bright field and phase contrast	✓	<b>✓</b>	OBB-A1398	150,-
	Infinity PH-Plan objective 10×	✓	✓	OBB-A1390	170,-
	Infinity PH-Plan objective 20×	0	0	OBB-A1391	200,-
	Infinity PH-Plan objective 40×	✓	✓	OBB-A1392	260,-
	Infinity PH-Plan objective 100×	0	0	OBB-A1393	345,-
Phase contrast units	PH slide 10×	✓	✓	OBB-A1399	75,-
	PH slide 20×	0	0	OBB-A1400	75,-
	PH slide 40×	✓	✓	OBB-A1401	75,-
	PH slide 100×	0	0	OBB-A1402	75,-
	Centering eyepiece	<b>✓</b>	✓		
Darkfield condenser	N.A. 0,85 – 0,91 (dry, paraboloid)	0	0	OBB-A1422	230,-
Illumination	20 W Halogen spare bulb (transmitted)	✓	✓	OBB-A1370	25,-
	Blue (built-in)	✓	✓		
Colour filters	Green	✓	✓	OBB-A1188	19,-
for transmitted illumination	Yellow	0	0	OBB-A1165	19,-
	Grey	0	0	OBB-A1183	19,-
	0,5× (focus adjustable)		0	OBB-A1515	185,-
C-Mount	1×		0	OBB-A1514	120,-
For further optional a	ccessories, please see the list of items for the OBL-12 and OBL-	13 series from page	21		





Easy connection to the PC, Laptop not included in delievery Please find the description of the software in chapter 9 (p. 97)

#### LAB LINE

# The digital all-round talent for your applications, documentation and live transmissions through an integrated camera

#### **Features**

- The KERN OBD is an excellent, stable laboratory microscope with integrated camera and infinity optical system, based on the OBL series
- A strong and continously adjustable 20 W halogen illumination unit (Philips) ensures the optimum lighting conditions
- The fixed, pre-centred and focusable 1,25
   Abbe condenser with aperture diaphragm and field diaphragm gives you a simplified Koehler illumination, without having to move the centre
- The large mechanical stage and its specimen holder holds up to two samples at the same time and is quick and easy to focus using a coaxial coarse and fine focusing knob on both sides

- A large selection of eyepieces, objectives and colour filters, a simple polarising unit as well as phase contrast units are available to you as accessories
- Multi-lingual software, USB cable, calibration scale as well as a protective dust cover, eye cups and multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

#### Scope of application

 Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, Sewage treatment plants, Oncology, entomology, vets, water analysis and breweries, if necessary, for training where there will be lots of people observing the screen/beamer at the same time

#### Applications/Samples

 Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue)

#### **Technical data**

- · Infinity optical system
- Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: One-sided
- Overall dimensions W×D×H 395×200×430 mm
- · Net weight approx. 7 kg

Model			Standard	d configuration		Price excl. of VAT	
KERN	RN Tube Eyepiece Objective quality Objectives Illumination						
OBD 127	Binocular/3MP digital/USB 2.0	HWF 10×/Ø 20 mm	Infinity E-Plan	4×/10×/40×/100×	20 W Halogen (transmitted)	990,- •	
OBD 128	Binocular/5MP digital/USB 2.0	HWF 10×/Ø 20 mm	Infinity E-Plan	4^/10^/40^/100^	20 W Halogen (transmitted)	1090,- 🔮	

ONLY WHILE STOCKS LAST

Price reduction

Eyepieces (23,2 mm)  H  Infinity E-Plan objectives  P P	HWF 10×/Ø 20 mm  WF 16×/Ø 13 mm  HWF 10×/Ø 20 mm (with Pointer)  4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  10×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm  Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable Interpupillary distance 50 – 75 mm	OBD 127	OBD 128  ✓✓  ○○  ○  ✓✓  ✓✓  ✓  ✓  ✓  ○  ○  ○  ○	OBB-A1404 OBB-A1354 OBB-A1448 OBB-A1161 OBB-A1159 OBB-A1160 OBB-A1158 OBB-A1250 OBB-A1270 OBB-A1437	excl. of VAT ex works  €  95,-  45,-  80,-  90,-  150,-  265,-  335,-  265,-  360,-  445,-
Eyepieces (23,2 mm)  H  Infinity E-Plan objectives  P P	WF 16×/Ø 13 mm  HWF 10×/Ø 20 mm (with Pointer)  A×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  40×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm  Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable	00 0 4 4 4 0	00 0 4 4 4 0	OBB-A1354  OBB-A1448  OBB-A1161  OBB-A1159  OBB-A1160  OBB-A1158  OBB-A1250  OBB-A1270	45,- 80,- 90,- 150,- 265,- 335,- 265,- 360,-
Infinity E-Plan objectives  P P P	HWF 10×/Ø 20 mm (with Pointer)  4×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  10×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm  Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable	0	· · · · · · · · · · · · · · · · · · ·	OBB-A1448  OBB-A1161  OBB-A1159  OBB-A1160  OBB-A1158  OBB-A1250  OBB-A1270	80,- 90,- 150,- 265,- 335,- 265,- 360,-
Infinity E-Plan objectives  P P P	1×/0,10 W.D. 12,1 mm  10×/0,25 W.D. 2,1 mm  10×/0,65 (spring-loaded) W.D. 0,58 mm  100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm  Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm  Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm  Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable	* * * * * * * * * * * * * * * * * * *	*     *     *     *     O	OBB-A1161  OBB-A1159  OBB-A1160  OBB-A1158  OBB-A1250  OBB-A1270	90,- 150,- 265,- 335,- 265,- 360,-
Infinity E-Plan objectives  P P	10×/0,25 W.D. 2,1 mm 40×/0,65 (spring-loaded) W.D. 0,58 mm 100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm	* * * * * * * * * * * * * * * * * * *	* * * * * * * * * * * * * * * * * * *	OBB-A1159  OBB-A1160  OBB-A1158  OBB-A1250  OBB-A1270	150,- 265,- 335,- 265,- 360,-
nfinity E-Plan objectives  P P P	100×/0,65 (spring-loaded) W.D. 0,58 mm 100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm	• • • • • • • • • • • • • • • • • • •	• • • • • • • • • • • • • • • • • • •	OBB-A1160 OBB-A1158 OBB-A1250 OBB-A1270	265,- 335,- 265,- 360,-
nfinity E-Plan objectives  P P	100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm   100×/0,40 (spring-loaded) W.D. 2,41 mm   100×/0,80 (spring-loaded) W.D. 0,33 mm   100×/1,15 (water) (spring-loaded) W.D. 0,18 mm   30° inclined/360° rotatable	0	0	OBB-A1158  OBB-A1250  OBB-A1270	335,- 265,- 360,-
E-Plan objectives P P P	Plan 20×/0,40 (spring-loaded) W.D. 2,41 mm Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm 30° inclined/360° rotatable	0	0	OBB-A1250 OBB-A1270	265,-
P P	Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm  Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable	0	0	OBB-A1270	360,-
P	Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm  30° inclined/360° rotatable				,
	30° inclined/360° rotatable	0	0	OBB-A1437	445,-
Trinocular (3MP)	Diopter adjustment: One-sided Light distribution 80:20 Built-in digital 3MP Camera with ½" CMOS USB port for PC without extra power supply With multilingual (DE, EN, FR, IT, ES) software "Microscope VIS" for Windows XP, Vista, 7, 8, 10	•		OBB-A1126	
Digital tube . Trinocular (5MP) .	30° inclined/360° rotatable Interpupillary distance 50 – 75 mm Diopter adjustment: One-sided Light distribution 80:20 Built-in digital 5 MP Camera with ½,5" CMOS USB port for PC without extra power supply With multilingual (DE, EN, FR, IT, ES) software "Microscope VIS" for Windows XP, Vista, 7, 8, 10		*	OBB-A1127	
	For calibrating the software measuring function, 0.01 mm division	<b>✓</b>	✓	OBB-A1224	70,-
Mechanical stage .	Stage size W×D 145×130 mm Travel 76×52 mm Two slide holder Coaxial coarse and fine focusing knobs, scale: 2 µm	<b>✓</b>	<b>*</b>		
Condenser A	Abbe N.A. 1,25 precentered (aperture diaphragm)	✓	✓	OBB-A1103	85,-
Darkfield condenser	N.A. 0,85 - 0,91 (dry, paraboloid)	0	0	OBB-A1422	230,-
Ilumination 2	20 W Halogen spare bulb (transmitted)	✓	✓	OBB-A1370	25,-
Polarising unit A	Analyser/Polariser	0	0	OBB-A1277	265,-
S	Single unit with ∞ PH-Plan objective 10×	0	0	OBB-A1215	310,-
Phase contrast	Single unit with ∞ PH-Plan objective 20×	0	0	OBB-A1217	350,-
unite (including	Single unit with ∞ PH-Plan objective 40×	0	0	OBB-A1219	420,-
N	Single unit with ∞ PH-Plan objective 100×	0	0	OBB-A1213	510,-
If	f required, there are several magnification levels, please conta	act our OPTICS-Tear	n		
	Blue (built-in)	✓	✓		
	Green	0	0	OBB-A1188	19,-
for transmitted   Y	/ellow	0	0	OBB-A1165	19,-
0	Grey	0	0	OBB-A1183	19,-





OBN-15: Mounted phase contrast condenser



Quintuple PH universal rotary condenser with  $10 \times /20 \times /40 \times /100 \times$ Infinity PH-Plan objectives (complete set,

for OBN-15 included)

#### **PROFESSIONAL LINE**

### Professionalism and versatility united in one microscope with Koehler illumination for demanding applications

#### **Features**

- · The OBN series stands out because of its unbeatable and consistently high quality and its ergonomic design. The range of modular components means that the OBN series can be individually customised for the professional user
- · Depending on the application, there is a choice of models with strong, continuously dimmable 3 W LED or 20 W halogen transmitted illumination (Philips)
- · In addition the halogen variant is available as a pre-configured phase contrast microscope, which, through the combination of a professional quintuple condenser wheel, phase contrast condenser and Infinity Plan phase contrast objectives makes it a high-quality, fully-equipped microscope for all applications related to this contrasting method
- · This series has a professional Koehler illumination unit with an adjustable field diaphragm as well as a height-adjustable 1,25 Abbe condenser which can be centred and which has an adjustable aperture diaphragm

- The extremely large mechanical stage with ergonomic, coaxial coarse and fine focusing knob on both sides enables you to adjust and focus your sample rapidly and accurately
- · A wide variety of modular systems, such as, for example, a swing-out condenser, various eyepieces, objectives, colour filters, phase contrast units, a darkfield condenser, a simple polarising unit, Butterfly tube, through to complete fluorescence units are available to you as accessories
- · This centring eyepiece for adjusting the phase contrast (OBN 158), a protective dust cover, eye cups as well as multi-lingual user instructions are included in the scope of delivery
- · A C-mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

· Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, Sewage treatment plants, Oncology, entomology, vets, water analysis and breweries

#### Applications/Samples

· Translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, bacteria, tissue)

#### **Technical data**

- · Infinity optical system
- Quintuple nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 390×200×400 mm
- · Net weight approx. 9 kg

#### STANDARD

















360°	TRINO	ABBE	HAL	ED LED	O PH	INFINITY	230 V	1 DAY
Mode	el							

Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €			
OBN 132	Trinocular	HWF 10×/Ø 20 mm	Infinity Plan	4×/10×/20×/	20 W Halogen (transmitted)	1370,-			
OBN 135	Trinocular	HWF 10×/Ø 20 mm	Infinity Plan	40×/100×	3 W LED (transmitted)	1390,-			
OBN 158	Trinocular	HWF 10×/Ø 20 mm	Infinity Plan	4×/PH10×/PH20×/ PH40×/PH100×	20 W Halogen (transmitted)	2540,-			

Model outfit			Model KERN	l	Order number	Price/piece excl. of VAT
	Thurs to 11 or		OBN 135	OBN 158		ex works €
Evepieces	HWF 10×/ø 20 mm	11	11	11	OBB-A1404	95,-
(23,2 mm)	WF 16×/ø 13 mm	00	00	00	OBB-A1354	45,-
	4×/0,10 W.D. 12,1 mm	✓	✓	✓	OBB-A1263	110,-
	10×/0,25 W.D. 4,64 mm	✓	✓	0	OBB-A1243	195,-
	20×/0,40 (spring-loaded) W.D. 2,41 mm		✓	0	OBB-A1250	265,-
Infinity	40×/0,66 (spring-loaded) W.D. 0,65 mm		✓	0	OBB-A1257	290,-
Plan achromatic objectives	100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm	✓	✓	0	OBB-A1240	315,-
	2,5×/0,07 W.D. 8,47 mm	0	0	0	OBB-A1247	265,-
	Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	0	0	0	OBB-A1270	360,-
	Plan 100×/1,15 (water) (spring-loaded) W.D. 0,18 mm	0	0	0	OBB-A1437	445,-
Trinocular tube	<ul> <li>Siedentopf 30° inclined/360° rotatable</li> <li>Interpupillary distance 50 - 75 mm</li> <li>Light distribution 100:0</li> <li>Diopter adjustment: Both-sided</li> </ul>	<b>✓</b>	<b>✓</b>	<b>√</b>		
irmocular tube	Butterfly 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 100:0     Diopter adjustment: Both-sided	0	0	0	OBB-A1382	520,-
Mechanical stage	Stage size W×D 175×145 mm     Travel 78×55 mm     Coaxial coarse and fine focusing knobs     Two slide holder		✓	~		
Condenser	Abbe N.A. 1,25 center-adjustable (aperture diaphragm)	✓	✓	0	OBB-A1102	85,-
Contract	Swing-out condenser N.A. 0,9/0,13 center-adjustable (aperture diaphragm)	0	0	0	OBB-A1104	220,-
Darkfield	N.A. 0,85 – 0,91 (dry, paraboloid)	0	0	0	OBB-A1421	265,-
condenser	N.A. 1,3 (oil, cardioid)	0	0	0	OBB-A1538	580,-
	20 W Halogen spare bulb (transmitted)	✓		0	OBB-A1370	25,-
Koehler illumination	3 W LED illumination system (transmitted) (non-rechargeable)		✓			
Polarising unit	Analyser/Polariser	0	0	0	OBB-A1283	265,-
	Quintuple hole turret with 10×/20×/40×/100× Infinity-PH-Plan objectives (complete set)	0	0	<b>✓</b>	OBB-A1237	1570,-
	Single unit with ∞ PH-Plan objective 10×	0	0		OBB-A1214	310,-
Phase contrast	Single unit with ∞ PH-Plan objective 20×	0	0		OBB-A1216	350,-
units	Single unit with ∞ PH-Plan objective 40×	0	0		OBB-A1218	420,-
	Single unit with ∞ PH-Plan objective 100×	0	0		OBB-A1212	510,-
	Centering eyepiece	0	0	✓		
	If required, there are several magnification levels, please cont	act our OPTIC	S-Team			
	1×	0	0	0	OBB-A1140	100,-
C-Mount	0,57× (focus adjustable)	0	0	0	OBB-A1136	160,-
	100 W HBO Epi Fluorescence unit 6-filter disc (UV/V/B/G) including centering objective	0	0	0	OBB-A1155	4340,-
Fluorescence unit	100 W HBO Epi Fluorescence unit, two-hole slide (B/G) including centering objective	0	0	0	OBB-A1153	3080,-
	3 W LED Epi Fluorescence unit (B/G) including centering objective	0	0	0	OBB-A1156	3110,-
	Blue	✓		✓		
Colour filters	Green	0	0	✓	OBB-A1188	19,-
for transmitted illumination	Yellow	0	0	0	OBB-A1165	19,-
	Grey	0	0	0	OBB-A1183	19,-





OBN 141/OBN 147



Illumination unit



Sextuple filter wheel OBN 148

#### **PROFESSIONAL LINE**

### The fluorescence microscope for the professional user

#### **Features**

- The fluorescence microscope in the OBN-14 series is based on the usual high quality and versatility of the OBN series. The outstanding, stable design in combination with high-quality optics set the standard in fluorescence microscopy in this class
- The powerful, dimmable 20W halogen illumination unit (Philips) and a 100W Epi fluorescence incident illumination unit on the OBN 147/OBN 148 models ensure perfect illumination and stimulation of your fluorescence samples
- · As an alternative, with the OBN 141 model we can offer you a fluorescence microscope with a 3W LED transmitted illumination unit and 3W LED Epi fluorescence incident illumination unit
- This series has a professional Koehler illumination unit with an adjustable field diaphragm as well as a height-adjustable 1,25 Abbe condenser which can be centred and which has an adjustable aperture diaphragm
- The extremely large mechanical stage with ergonomic, coaxial coarse and fine focusing knob on both sides enables you to adjust and focus your sample rapidly and accurately

- With the OBN 147/OBN 148 halogen variant you have a filter wheel which has up to 6 fittings. As standard this is fitted with a B/G or B/G/UV/V fluorescence filter. The OBN 141 LED variant is fitted with a B/G fluorescence filter with a changeover slider as standard. The changeover slider and the filter wheel mean that you can change the stimulation filter quickly
- · A large selection of eyepieces, objectives, colour filters, darkfield condensers as well as a Butterfly tube, polarising and phase contrast units can easily be integrated thanks to the modular construction system
- · The centring objective for adjusting the fluorescence, a protective dust cover, eye cups as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

· Haematology, urology, gynaecology, dermatology, pathology, microbiology and parasitology, immunology, Sewage treatment plants, Oncology, entomology, vets, water analysis and breweries

#### Applications/Samples

· Specially for translucent, thin, low-contrast, challenging samples (e.g. immunofluorescence, FISH, DAPI staining, etc.)

#### **Technical data**

- · Infinity optical system
- · Quintuple nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 530×220×490 mm
- · Net weight approx. 23 kg

STANDARD



























PTION			
0	•		inin
PH	DF	POLAR	SCALE

Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €			
OBN 141	Trinocular	HWF 10×/Ø 20 mm	Infinity Plan	4× /10× /20× /	LED + 3 W LED Epi Fluorescence (B/G)	4470,-			
OBN 147	Trinocular	WF 10×/Ø 20 mm	Infinity Plan	4×/10×/20×/ 40×/100×	Halogen + 100 W Epi Fluorescence (B/G)	4390,-			
OBN 148	Trinocular	HWF 10×/Ø 20 mm	Infinity Plan	40.7 100	Halogen + 100 W Epi Fluorescence (B/G/UV/V)	5630,-			

Model outfit			Model KERN	I	Order number	Price/piece				
				OBN 148		excl. of VAT ex works €				
	HWF 10×/ø 20 mm	44		11	OBB-A1404	95,-				
Eyepieces	WF 10×/Ø 20 mm		44		OBB-A1351	85,-				
(23,2 mm)	WF 16×/Ø 13 mm		00	00	OBB-A1354	45,-				
	WF 10×/Ø 20 mm (reticule 0,1 mm) (adjustable)	0	0	0	OBB-A1352	140,-				
	4×/0,10 W.D. 12,1 mm	<b>✓</b>	✓	✓	OBB-A1263	110,-				
	10×/0,25 W.D. 4,64 mm		✓	✓	OBB-A1243	195,-				
Infinity	20×/0,40 (spring-loaded) W.D. 2,41 mm	<b>✓</b>	<b>✓</b>	✓	OBB-A1250	265,-				
Plan achromatic	40×/0,66 (spring-loaded) W.D. 0,65 mm	<b>✓</b>	<b>✓</b>	✓	OBB-A1257	290,-				
objectives	100×/1,25 (oil) (spring-loaded) W.D. 0,19 mm	<b>✓</b>	✓	<b>✓</b>	OBB-A 1240	315,-				
	2,5×/0,07 W.D. 8,47 mm	0	0	0	OBB-A1247	265,-				
	Plan 60×/0,80 (spring-loaded) W.D. 0,33 mm	0	0	0	OBB-A1270	360,-				
Trinocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 100:0     Diopter adjustment: Both-sided     Butterfly 30° inclined/360° rotatable	*	~	~		,				
	Interpupillary distance 50 – 75 mm     Light distribution 100:0     Diopter adjustment: Both-sided	0	0	0	OBB-A1382	520,-				
Mechanical stage	<ul> <li>Stage size W×D 175×145 mm</li> <li>Travel 78×55 mm</li> <li>Coaxial coarse and fine focusing knobs</li> <li>Two slide holder</li> </ul>	*	<b>✓</b>	<b>✓</b>						
	Abbe N.A. 1,25 center-adjustable (aperture diaphragm)	✓	✓	✓	OBB-A1102	85,-				
Condenser	Swing-out condenser N.A. 0,9/0,13 center-adjustable (aperture diaphragm)		0	0	OBB-A1104	220,-				
Darkfield	N.A. 0,85 – 0,91 (dry, paraboloid)	0	0	0	OBB-A1421	265,-				
condenser	N.A. 1,3 (oil, cardioid)	0	0	0	OBB-A 1538	580,-				
Koehler illumination	20 W Halogen spare bulb (transmitted)	✓	✓	✓	OBB-A1370	25,-				
Polarising unit	Analyser/Polariser	0	0	0	OBB-A1283	265,-				
	Quintuple hole turret with 10×/20×/40×/100× Infinity-PH-Plan objectives (complete set)	0	0	0	OBB-A1237	1570,-				
	Single unit with ∞ PH-Plan objective 10×	0	0	0	OBB-A1214	310,-				
Phase contrast	Single unit with ∞ PH-Plan objective 20×	0	0	0	OBB-A1216	350,-				
units	Single unit with ∞ PH-Plan objective 40×	0	0	0	OBB-A1218	420,-				
	Single unit with ∞ PH-Plan objective 100×	0	0	0	OBB-A1212	510,-				
	If required, there are several magnification levels, please contact our OPTICS-Team									
	1×	0	0	0	OBB-A1140	100,-				
C-Mount	0,57× (focus adjustable)	0	0	0	OBB-A1136	160,-				
	100 W HBO Epi Fluorescence unit 6-filter disc (UV/V/B/G) including centering objective			✓						
Fluorescence unit	100 W HBO Epi Fluorescence unit, two-hole slide (B/G) including centering objective		<b>✓</b>							
	3 W LED Epi Fluorescence unit (B/G) including centering objective	<b>✓</b>								
	Blue	✓	✓	✓						
Colour filters	Green	0	0	0	OBB-A1188	19,-				
for transmitted illumination	Yellow	0	0	0	OBB-A1165	19,-				
	Grey	0	0	0	OBB-A1183	19,-				





N.A. 0,3 Abbe Condenser with phase contrast slide



Coaxial control knobs for x/y can be fitted either left or right

#### LAB LINE

### The inverted biological laboratory microscope – also with fluorescence

#### **Features**

- The OCM range stands out through its design which is ergonomic, robust and extremely stable. This design, with its large working distance, is particularly suitable for the monitoring and analysis of cell cultures, for example
- A strong and continuously adjustable 30W halogen illumination unit ensures the optimum illumination in the bright field of your samples. An additional Osram 100 W Epi fluorescence illumination unit is available to you as a fluorescence microscope (OCM 165) for perfect illumination and excitation of your fluorescence samples
- · A special Abbe N.A. 0.3 condenser with aperture diaphragm and large working distance of 72 mm guarantees the very best working practise in the bright field and with fluorescence applications
- · As standard, the OCM range is fitted with a trinocular eyepiece tube

- The mechanical stage including specimen holder (Ø 118 mm) means that you can work quickly and effectively. Further brackets for petri dishes are included with delivery or available as accessories
- · Further options such as, for example, a selection of eyepieces, objectives, specimen holders and other phase contrast units can be integrated as accessories
- · A dust cover as well as user instructions are included with the delivery
- · Please find detailed information in the following model outfit list

#### Scope of application

· Research and breeding of cell cultures and tissue cultures

#### Applications/Samples

• Particularly for viewing samples in culture vessels (flasks, petri dishes, microtitre plates), translucent, thin, low-contrast, challenging samples (e.g. living mammal cells, tissue, microorganisms if necessary, immunofluorescence, FISH, DAPI staining etc.)

#### **Technical data**

- · Infinity optical system
- Quintuple nosepiece
- · Siedentopf 45° inclined
- · Diopter adjustment: Both-sided

- · Overall dimensions W×D×H 304×599×530 mm
- · Net weight approx. 13,5 kg

#### **OCM 165**

- Overall dimensions W×D×H 304×782×530 mm
- · Net weight approx. 21 kg






















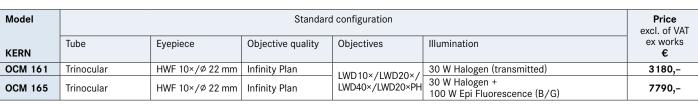

Model outfit		Mode	I KERN	Order number	Price/piece excl. of VAT
		OCM 161	OCM 165		ex works €
Eyepieces	HWF 10×/ø 22 mm (adjustable)	44	<b>**</b>	OBB-A1491	90,-
(30 mm)	HWF 10×/Ø 20 mm (reticule 0,1 mm) (adjustable)	0	0	OBB-A1523	140,-
Indinia.	4×/0,11 W.D. 12,1 mm	0	0	OBB-A1493	135,-
Infinity Plan achromatic objectives for long working	10×/0,25 W.D. 8,3 mm	<b>✓</b>	<b>✓</b>	OBB-A1494	160,-
	20×/0,40 W.D. 7,2 mm	<b>✓</b>	<b>✓</b>	OBB-A1495	215,-
distance	40×/0,60 W.D. 3,4 mm	✓	✓	OBB-A1496	280,-
Trinocular tube	45° inclined     Interpupillary distance 48-76 mm     Light distribution 100:0     Diopter adjustment: Both-sided	~	~		
Mechanical stage	Stage size W×D 210×241 mm Travel128×80 mm Coaxial coarse and fine focusing knobs The x/y control knobs can be fitted either left or right Suitable for attaching a 96-hole microtitre plate	~	<b>~</b>		
S	Drop specimen holder (Ø 110)	✓	✓	OBB-A1503	30,-
	Specimen holder for 35 mm culture dish	0	0	OBB-A1505	60,-
	Specimen holder for 54 mm culture dish	<b>✓</b>	<b>✓</b>	OBB-A1506	60,-
	Specimen holder for 65 mm culture dish	0	0	OBB-A1507	60,-
Condenser	Abbe N.A. 0,3 (aperture diaphragm), LWD 72 mm	✓	✓		
Illumination	30 W Halogen spare bulb (transmitted)	✓	<b>✓</b>	OBB-A1372	35,-
	Phase contrast slide (universal)	✓	✓	OBB-A1500	100,-
	Infinity PH-Plan objective 10×	0	0	OBB-A1497	215,-
Phase contrast units	Infinity PH-Plan objective 20×	✓	<b>✓</b>	OBB-A1498	280,-
	Infinity PH-Plan objective 40×	0	0	OBB-A1499	330,-
	Centering eyepiece	✓	<b>✓</b>	OBB-A1544	110,-
Fluorescence unit	100 W HBO Epi Fluorescence unit, two-hole slide (B/G)		✓		
	Blue	✓	✓	OBB-A1510	19,-
Colour filters	Green	✓	✓	OBB-A1511	19,-
for transmitted illumination	Yellow	0	0	OBB-A1512	19,-
	Grey	0	0	OBB-A1513	19,-
C Marriet	0,5×	0	0	OBB-A1515	185,-
C-Mount	1×	0	0	OBB-A1514	120,-



# TIPS FROM THE EXPERTS

In addition to our range of microscopes and refractometers, we also offer you a comprehensive range of professional measuring and testing technology as well as a comprehensive range of medical scales.

We will offer you reliable quality at fair prices and with the shortest delivery times. Our product specialists will give you professional advice, will work with you to find the right product and will provide comprehensive support after the purchase too.

This ensures a high level of investment security and a good feeling – KERN, the nice balance manufacturer from the wilds of Southern Germany.

For all details on our ranges, please see **www.kern-sohn.com** or request your own personal catalogues right now: Telephone +49 7433 9933-0

# # 1 Pure competence in balances

"KERN offers you a complete, carefully-designed range of laboratory balances, analytical balances, moisture analysers, industrial scales and test weights. And all this at an extremely attractive price with the same high level of quality.

Browse and find products in the KERN 2020 catalogue for Balances & Test service."



Stephan Ade, Head of Sales Balances and test service



# 2 Metallurgical microscopes







Illumination unit with filter disc



Stage and objectives

#### LAB LINE MET

### The metallurgical reflected light microscope for material testing and surface testing, as well as quality assurance in industry

#### **Features**

- The KERN OKM is an excellent metallurgical reflected light microscope, e.g. for surface quality testing of raw materials and finished products in industry
- The strong, continuously dimmable 30 W halogen reflected illumination unit (Philips) ensures excellent, high-contrast images
- · The illumination unit with an integrated 5-slot filter wheel for blue, green, yellow, grey and blank means that you can quickly change the colour filter for different contrast views
- · A large mechanical stage for reflected illumination applications is configured as standard. The coarse and fine focusing knob on both sides guarantees optimal adjustment and focusing of your sample

- A simple polarising unit (analyser and polariser) is included with delivery
- · A large selection of different eyepieces, objectives and a polarising unit are also
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

· Metallurgy, material testing, quality assurance

#### Applications/Samples

· Opaque and thick samples, workpieces (surfaces, fold lines, coatings)

#### **Technical data**

- · Infinity optical system
- · Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: One-sided
- · Overall dimensions W×D×H 440×200×460 mm
- · Net weight basic configuration approx. 8 kg



Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €			
OKM 172*	Binocular	HWF 10×/Ø 18 mm	Infinity Plan	5×/10×/ LWD 20×/	30 W Halogen (incident)	1790,-			
OKM 173	Trinocular	HWF 10×/Ø 18 mm	Infinity Plan	LWD40×	30 W Halogen (incident)	1890,-			

Model outfit		Mode	Model KERN		Price/piece excl. of VAT
		OKM 172	OKM 173	_	ex works €
	HWF 10×/ø 18 mm	✓	✓	OBB-A1403	45,-
	HWF 10×/Ø 18 mm (reticule 0,1 mm) (non-adjustable)	✓	✓	OBB-A1349	60,-
Eyepieces (23,2 mm)	WF 5×/ø 20 mm	0	0	OBB-A1355	55,-
(==,=,	WF 12,5×/ø 14 mm	0	0	OBB-A1353	75,-
	WF 16×/ø 13 mm	0	0	OBB-A1354	45,-
	5×/0,11 W.D. 12,10 mm	✓	✓	OBB-A1268	105,-
Infinity Plan achromatic	10×/0,25 W.D. 4,75 mm	<b>✓</b>	✓	OBB-A1244	195,-
objectives (no cover glass)	20×/0,40 (spring-loaded) W.D. 2,14 mm	0	0	OBB-A1251	265,-
(110 dover glads)	40×/0,65 (spring-loaded) W.D. 0,45 mm	0	0	OBB-A1258	290,-
Infinity	20×/0,40 W.D. 8,35 mm	✓	✓	OBB-A1252	290,-
Plan achromatic objectives (no cover glass) for long working	40×/0,65 W.D. 3,90 mm	✓	✓	OBB-A1259	405,-
	50×/0,70 (spring-loaded) W.D. 1,95 mm	0	0	OBB-A1266	450,-
distance	80×/0,80 (spring-loaded) W.D. 0,85 mm	0	0	OBB-A1271	500,-
Binocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Diopter adjustment: One-sided	✓	0	OBB-A1130	
Trinocular tube	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 80:20     Diopter adjustment: One-sided	0	~	OBB-A1346	
Mechanical stage	Stage size W×D 200×140 mm     Travel 76×52 mm     Coaxial coarse and fine focusing knobs	<b>✓</b>	<b>✓</b>		
Illumination	30 W Halogen spare bulb (incident)	✓	✓	OBB-A1372	35,-
Deflected	5-filter unit (Blue, Green, Yellow, Grey, Empty)	✓	✓		
Reflected illumination unit	Polarising unit (Incl. analyser and polariser slide)	✓	✓		
C.Marriet	1×		0	OBB-A1514	100,-
C-Mount	0,5× (focus adjustable)		0	OBB-A1515	160,-





Stage OKN



Stage OKO



Illumination unit

#### **PROFESSIONAL LINE MET**

### The fully-equipped reflected and transmitted light microscope for numerous applications in metallurgy

#### **Features**

- · The KERN OKN and OKO series are professional, versatile, metallurgical microscopes, which are used in testing metals and analysing surfaces
- · You can choose between two reflected illumination units: A 50 W halogen reflected illumination unit or a premium illumination unit with stunning 100 W reflected illumination for powerful performance
- A height-adjustable 1,25 Abbe condenser which can be centred as well as a field diaphragm are available for the transmitted light variants (KERN OKO series), for complete professional Koehler illumination
- The KERN OKO transmitted illumination variant is fitted with an open, mechanical stage, as standard. On the other hand, the KERN OKN reflected illumination variant has a closed, mechanical stage, as standard

- A simple polarising unit (analyser and polariser) is included with delivery
- · A large selection of accessories, such as, for example, a Butterfly tube, eyepieces and further objectives are available for longer working distances
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-Mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

· Metallurgy, material testing, quality assurance

#### Applications/Samples

· Opaque and thick samples, workpieces (surfaces, fold lines, coatings)

#### **Technical data**

- · Infinity optical system
- Quintuple nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 550×200×460 mm
- Net weight basic configuration approx. 14,5 kg





















Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €			
OKN 175*	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	5×/10×/LWD20×/ LWD40×	50 W Halogen (incident)	2230,-			
OKO 176*	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	5×/10×/LWD20×/ LWD40×/100×	50 W Halogen (incident) + 20 W (transmitted)	2770,-			
OKN 177	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	5×/10×/LWD20×/ LWD40×	100 W Halogen (incident)	2450,-			
ОКО 178	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	5×/10×/LWD20×/ LWD40×/100×	100 W Halogen (incident) + 20 W (transmitted)	2990,-			

Metallurgical microscopes

Model outfit			Model	KERN		Order number	Price/piece
		OKN 175	OKO 176	OKN 177	OKO 178	_	ex works €
	WF 10×/ø 18 mm	11	11	11	11	OBB-A1347	40,-
	WF 10×/Ø 18 mm (reticule 0,1 mm) (adjustable)	✓	✓	✓	✓	OBB-A1350	90,-
Eyepieces (23,2 mm)	WF 5×/ø 20 mm	00	00	00	00	OBB-A1355	55,-
,,	WF12,5×/ø 14 mm	00	00	00	00	OBB-A1353	75,-
	WF 16×/ø 13 mm	00	00	00	00	OBB-A1354	45,-
	5×/0,11 W.D. 6,73 mm	✓	✓	✓	✓	OBB-A1268	105,-
	10×/0,25 W.D. 4,19 mm	✓	✓	✓	✓	OBB-A1244	195,-
Infinity Plan objectives	20×/0,40 (spring-loaded) W.D. 2,14 mm	0	0	0	0	OBB-A1251	265,-
(no cover glass)	40×/0,65 (spring-loaded) W.D. 0,45 mm	0	0	0	0	OBB-A1258	290,-
	100×/1,25 (oil) (spring-loaded) W.D. 0,12 mm		✓		✓	OBB-A1241	315,-
indinia.	20×/0,40 W.D. 8,35 mm	✓	✓	✓	✓	OBB-A1252	290,-
Infinity Plan objectives	40×/0,65 W.D. 3,90 mm	<b>✓</b>	✓	✓	✓	OBB-A1259	405,-
(no cover glass) for long working	50×/0,70 (spring-loaded) W.D. 1,95 mm	0	0	0	0	OBB-A1266	450,-
distance	80×/0,80 (spring-loaded) W.D. 0,85 mm	0	0	0	0	OBB-A1271	500,-
	Siedentopf 30° inclined/360° rotatable     Interpupillary distance 50 – 75mm     Light distribution 100:0     Diopter adjustment: Both-sided	~	~	~	~		
Trinocular tube	Butterfly 30° inclined/360° rotatable     Interpupillary distance 50 – 75 mm     Light distribution 100:0     Diopter adjustment: Both-sided	0	0	0	0	OBB-A1382	520,-
Mechanical stage for reflection	Stage size W×D 200×140 mm Travel 78×55 mm Stage fast lowering unit Stage Up-Down moving range: max. 50 mm	~		~			
Mechanical stage for transmitted Illumination	Stage size W×D 175×145 mm     Travel 78×55 mm     Coaxial coarse and fine focusing knobs		✓		✓		
Stage plate	Plate for sample placement	✓	✓	✓	✓		
Glass plate	Glass plate		0		0	OBB-A1378	50,-
Reflected illumination unit	Polarising unit (Incl. analyser, polariser and blue filter slide)	✓	✓	✓	✓		
Condenser	Abbe N.A. 1,25 (aperture diaphragm)		✓		✓	OBB-A1380	100,-
	20 W Halogen spare bulb (transmitted)		✓		✓	OBB-A1370	25,-
Illumination	50 W Halogen spare bulb (incident)	✓	✓			OBB-A1207	19,-
	100 W Halogen spare bulb (incident)			✓	✓	OBB-A1377	35,-
Polarising unit	for transmitted illumination		✓		✓	OBB-A1470	135,-
	Blue		✓		✓		
Colour filters	Green		0		0	OBB-A1188	19,-
or transmitted Ilumination	Yellow		0		0	OBB-A1165	19,-
nunililativil	Grey		0		0	OBB-A1183	19,-
	1×	0	0	0	0	OBB-A1140	100,-
C-Mount	0,57× (focus adjustable)	0	0	0	0	OBB-A1136	160,-





Specimen stage and illumination unit



Analyser/Polariser

#### **LAB LINE MET**

#### The inverted metallurgical microscope for professional applications

#### **Features**

- The KERN OLM range is part of the range of inverted microscopes and stands out through its design which is ergonomic, robust and extremely stable. This range, with its large working distance is, for example, particularly suitable for surface quality testing of raw materials and finished products in industry
- · Strong and continuously adjustable 50W halogen illumination unit ensures the optimum illumination of the materials to be
- · As standard, the OLM range is fitted with a trinocular eyepiece tube
- · A simple polarising unit (analyser and polariser) is included with delivery

- · A large mechanical stage is included with delivery as standard. The coarse and fine focusing knob on both sides guarantees optimal adjustment and focusing
- Further options such as, for example, a large selection of objectives can be integrated as accessories
- · A dust cover as well as user instructions are included with the delivery
- · Please find detailed information in the following model outfit list

#### Scope of application

· Metallurgy, material testing, quality assurance

#### Applications/Samples

• Opaque and thick samples, workpieces (surfaces, fold lines, coatings)

#### **Technical data**

- · Infinity optical system
- Quintuple nosepiece
- · Siedentopf 30° inclined
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 271×379×747 mm
- · Net weight approx. 12,5 kg

















þ	
٧	1 DAY

Model		Standard configuration						
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €		
OLM 171	Trinocular	HWF 10×/Ø 22 mm	Infinity Plan	LWD5×/LWD10×/ LWD20×/LWD50×	50 W Halogen (incident)	3190,-		

Model outfit		Model KERN	Order number	Price/piece excl. of VAT	
		OLM 171		ex works €	
Eyepieces	HWF 10×/ø 22 mm (adjustable)	✓	OBB-A1491	90,-	
(30 mm)	HWF 10×/Ø 22 mm (reticule 0,1 mm) (adjustable)	✓	OBB-A1523	140,-	
	5×/0,13 W.D. 24,23 mm	✓	OBB-A1525	125,-	
Infinity	10×/0,25 W.D. 18,48 mm	✓	OBB-A1526	185,-	
Plan achromatic	20×/0,40 W.D. 8,35 mm	✓	OBB-A1527	245,-	
<b>objectives</b> for long working	50×/0,70 (spring-loaded) W.D. 1,95 mm	✓	OBB-A1528	320,-	
distance	80×/0,80 (spring-loaded) W.D. 0,85 mm	0	OBB-A1530	480,-	
	100×/0,90 (dry) W.D. 1,0 mm	0	OBB-A1531	720,-	
Trinocular tube	30° inclined     Interpupillary distance 48-76 mm     Light distribution 100:0     Diopter adjustment: Both-sided	<b>√</b>			
Mechanical stage	Stage size W×D 210×180 mm     Travel 50×50 mm     Coaxial coarse and fine focusing knobs	✓			
Illumination	50 W Halogen spare bulb (incident)	✓	OBB-A1207	19,-	
Reflected illumination unit	Polarising unit (Incl. analyser, polariser and colour filter slide)	<b>✓</b>			
	Blue	✓	OBB-A1510	19,-	
Colour filters	Green	0	OBB-A1511	19,-	
for transmitted illumination	Yellow	0	OBB-A1512	19,-	
	Grey	0	OBB-A1513	19,-	
0 M 1	0,5×	0	OBB-A1515	185,-	
C-Mount	1×	0	OBB-A1514	120,-	

<sup>✓ =</sup> Included with delivery


# 3 Polarising microscopes









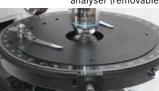

Center-adjustable and turnable polarisation stage



"Swing-Out" condenser

#### **PROFESSIONAL LINE POL**

### The flexible and powerful polarising microscope for all professional applications with reflected and transmitted light

#### **Features**

- · These devices are professional, fully-equipped polarising microscopes, which use the polarisation of light to analyse minerals, crystals and isotropic materials
- You can choose between a pure transmitted light variant (KERN OPM), a pure reflected light variant (KERN OPN) and a combi variant (KERN OPO). A complete Koehler illumination is integrated into all series as standard
- · As standard, the KERN OPM and OPO transmitted illumination variants have a height-adjustable 0,9/0,13 swing-out Abbe condenser which can be centred, for complete Koehler illumination
- A 360° revolving stage with 1° division, 6' fine division and locking function is integrated into all series as standard

- · As standard all series are fitted with a complete polarising unit with scale, a Bertrand lens, a  $\lambda$  +  $\frac{1}{4}\lambda$  Slip as well as a quartz wedge
- · A large selection of accessories such as, for example, a mechanical stage attachment as well as further objectives for a long working distance and filter units are also available
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-Mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

• Mineralogy, texture observations, material testing, observation of crystals

#### Applications/Samples

· More complex samples with polarising properties

#### **Technical data**

- · Infinity optical system
- Quintuple nosepiece
- Siedentopf 30° inclined
- · Diopter adjustment: One-sided
- · Overall dimensions W×D×H 500×200×500 mm
- · Net weight approx. 14,5 kg

#### STANDARD

















Model		Standard configuration							
KERN	Tube	Eyepiece	Objective quality	Objectives	Illumination	ex works €			
OPM 181	Trinocular	WF 10×/ø 20 mm	Infinity Plan	Non-stress 4×/10×/20×/40×	20 W Halogen (transmitted)	3430,-			
OPN 184	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	Non-stress 4×/10×/20×/40×	100 W Halogen (incident)	3640,-			
OPO 185	Trinocular	WF 10×/Ø 18 mm	Infinity Plan	Non-stress 4×/10×/20×/40×/60×	100 W Halogen (incident) + 20 W (transmitted)	4980,-			

Model outfit			Model KERN	N	Order number	Price/piece excl. of VAT
		OPM 181	OPN 184	OPO 185	_	ex works €
	WF 10×/18 mm		✓	✓	OBB-A1347	40,-
Eyepieces	WF 10×/18 mm (reticule 0,1 mm) (adjustable)		✓	✓	OBB-A1464	90,-
(23,2 mm)	WF 10×/20 mm	✓			OBB-A1351	85,-
	WF 10×/20 mm (reticule 0,1 mm) (adjustable)	<b>✓</b>			OBB-A1465	140,-
	4×/0,10 W.D. 12,1 mm	<b>✓</b>	✓	1	OBB-A1294	105,-
	10×/0,25 W.D. 4,64 mm	✓	✓	✓	OBB-A1289	200,-
Non-stress	20×/0,40 (spring-loaded) W.D. 2,41 mm	✓	✓	✓	OBB-A1290	265,-
Infinity Plan objectives	40×/0,65 (spring-loaded) W.D. 0,65 mm	✓		✓	OBB-A1292	310,-
	40×/0,65 (spring-loaded) (no cover glass) W.D. 3,9 mm	0	✓	0	OBB-A1288	390,-
	60×/0,80 (spring-loaded) W.D. 0,33 mm	0	0	✓	OBB-A1296	370,-
Indinia.	20×/0,40 W.D. 8,35 mm	0	0	0	OBB-A1291	345,-
Infinity Plan objectives	40×/0,65 W.D. 3,90 mm	0	0	0	OBB-A1293	390,-
(no cover glass) for long working	50×/0,70 (spring-loaded) W.D. 1,95 mm	0	0	0	OBB-A1295	415,-
distance	80×/0,80 (spring-loaded) W.D. 0,85 mm	0	0	0	OBB-A1297	500,-
Trinocular tube	<ul> <li>Siedentopf 30° inclined</li> <li>Interpupillary distance 50 – 75 mm</li> <li>Light distribution 100:0</li> <li>Diopter adjustment: One-sided</li> </ul>	<b>✓</b>	<b>✓</b>	<b>✓</b>		
Professional dedicated polarising trinocular head	icated polarising lo keep the reticular cross in the right-hand eyepiece in the		0	0	OBB-A1210	550,-
Analyser unit with scale	360° rotatable, lockable	✓	✓	✓		
Bertrand lens	Built-in, center-adjustable	✓	✓	✓	OBB-A1121	305,-
λ + ¼ λ Slip	λ Slip and ¼ λ Slip (combination)	✓	✓	✓	OBB-A1316	140,-
Quartz wedge	I - IV Class	✓	✓	✓	OBB-A1321	240,-
Revolving round stage	360° rotatable, center-adjustable, division 1°, Vernier division 6'	<b>✓</b>	<b>✓</b>	<b>✓</b>		
Polarising attached mechanical stage	Polarising attached mechanical stage	0	0	0	OBB-A1337	270,-
Swing-out condenser	N.A. 0,9/0,13 swing-out achromatic condenser (aperture diaphragm)	<b>✓</b>		<b>✓</b>	OBB-A1107	135,-
Polarising unit with scale (transmitted)	360° rotatable, lockable	✓		<b>✓</b>		
Koehler illumination	20 W Halogen spare bulb (transmitted)	✓		✓	OBB-A1370	25,-
Reflecting polarising unit replacement	12 V/50 W Halogen		0	0	OBB-A1207	19,-
bulb	12 V/100 W Halogen		✓	✓	OBB-A1377	35,-
	Blue	✓		✓	OBB-A1170	19,-
Colour filters for transmitted	Green	0		0	OBB-A1188	19,-
illumination	Yellow	0		0	OBB-A1165	19,-
	Grey	0		0	OBB-A1183	19,-
C Marriet	1×	0	0	0	OBB-A1140	100,-
C-Mount	0,57× (focus adjustable)	0	0	0	OBB-A1136	160,-





#### Cleaning sets for microscopes

#### **Features**

- This economical and fully equipped 7-piece cleaning set contains everything you need for the very best care of your microscope
- A silicon hand blower, dust brush, 60 ml of cleaning liquid, lint-free duster, optical cleaning cloths and cleaning swabs. You get all that in a high-quality KERN storage bag which you can also easily fix onto your belt
- You can use this set not only to gently clean your microscope, but also for example your camera, binoculars or all other optical surfaces

Model KERN	Description	Price excl. of VAT ex works €
OCS 901	7-piece cleaning sets for microscopes und other optical instruments	25,-

### Microscopy, measuring technology and testing services from a single source



Optical instruments catalogue



Balances & test service catalogue



Medical scales catalogue



SAUTER measuring equipment catalogue



DAkkS calibration service brochure



# TIPS FROM THE EXPERTS

In addition to our range of microscopes and refractometers, we also offer you a comprehensive range of professional measuring and testing technology as well as a comprehensive range of medical scales.

We will offer you reliable quality at fair prices and with the shortest delivery times. Our product specialists will give you professional advice, will work with you to find the right product and will provide comprehensive support after the purchase too.

This ensures a high level of investment security and a good feeling – KERN, the nice balance manufacturer from the wilds of Southern Germany.

For all details on our ranges, please see **www.kern-sohn.com** or request your own personal catalogues right now: Telephone +49 7433 9933-0

# # 1 Medical products from KERN – products you can rely on

For years an established name in hospitals, doctors' surgeries, rehabilitation clinics and nursing homes.

Wherever reliable quality is important, look no further than the complete KERN range of medical scales, from baby scales, to personal scales, chair scales, obesity scales through to hand grip dynamometers.

KERN 2020 catalogue for medical balances and devices – this is where you will find everything for day-to-day tasks in the medical environment!



Hendrik Neff, Product specialist Medical Scales



# **Stereomicroscopes**Stereo, Stereo-Zoom, Coaxial and Gem microscopes







With white stage plate



With black stage plate

#### **EDUCATIONAL LINE**

## The small robust model for school, training establishment or workshops

#### **Features**

- The KERN OSE-4 is an extremely robust, stable stereo microscope which is easy to use, it is ideal for all conventional applications in schools, workshops and training companies
- Depending on the model, optional LED reflected illumination as well as transmitted and reflected illumination ensure the very best illumination of your sample
- Despite its low price it has very good optical characteristics, which enable you to have sharp images over a large field of view
- A turnable objective with predefined magnifications is available to make your working procedures quicker and more efficient

- The eyepieces are fixed in the eyepiece tube, to stop them getting damaged or lost
- A special feature of this adaptable and yet robust microscope series is the stable mechanism of the microscope stand which can be adjusted precisely. It will also impress you with its functionality and ergonomic design
- A large selection of eyepieces as well as various additional external illumination units are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

#### Scope of application

 Training, in vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

• Samples with focus on three-dimensional impression (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- Optical system: Greenough optics
- Tube 45° inclined
- Interpupillary distance 55 75 mm
- Diopter adjustment: One-sided
- Overall dimensions W×D×H 230×130×330 mm
- Net weight approx. 2 kg

#### STANDARI



Model		Standard configuration					
	Tube Eyepiece Field of view Objective Stand Illumination				excl. of VAT ex works		
KERN			mm				€
OSE 416	Binocular	WF 10×/Ø 20 mm	Ø 20	1×/3×	Pillar style	0,21 W LED (incident); 0,21 W LED (transmitted)	255,-
OSE 417	Binocular	WF 10×/Ø 20 mm	Ø 20	2×/4×	Pillar style	0,21 W LED (incident); 0,21 W LED (transmitted)	255,-

Eyepiece		Specifications – Objectives								
	Magnification	1×	2×	3×	4×					
ME 5	Total magnification	5×	10×	15×	20×					
WF 5×	Field of view mm	ø 20	ø 10	Ø 6,7	Ø 5					
WE 40	Total magnification	10×	20×	30×	40×					
WF 10×	Field of view mm	ø 20	ø 10	Ø 6,7	Ø 5					
ME 45	Total magnification	15×	30×	45×	60×					
WF 15×	Field of view mm	ø 15	Ø 7,5	Ø 5	Ø 3,7					
WE 00	Total magnification	20×	40×	60×	80×					
WF 20×	Field of view mm	ø 10	Ø 6,5	Ø 4,3	Ø 3,2					
Working distance	Working distance		57 mm	57 mm	57 mm					

Model outfit		Mode	I KERN	Order number	<b>Price</b> /piece excl. of VAT
		OSE 416	OSE 417		ex works €
	WF 5×/ø 16,2 mm	00	00	OZB-A4101	35,-
Eyepieces	WF 10×/ø 20 mm	<b>44</b>	11	OZB-A4102	30,-
(30,5 mm)	WF 15×/ø 15 mm	00	00	OZB-A4103	30,-
	WF 20×/ø 10 mm	00	00	OZB-A4104	35,-
Stand	Pillar style, with 0,21 W LED illumination (transmission + incident)	✓	✓		
04	Frosted glass/ø 95 mm	✓	✓	OZB-A4805	19,-
Stage plate	Black-white/Ø 95 mm	✓	✓	OZB-A4806	19,-
External illumination	Please find the information about external illumination units in	the catalogue on p	age 90 and on our v	vebsite www.kern-sc	hn.com





Side view

The successor of the well-known series OSF-4

#### **EDUCATIONAL LINE**

### Stereo microscope with robust, ergonomic design, ideal for workshops, schools and training

#### **Features**

- · With its integrated handle as well as its stable arm curved stand, the KERN OSE OSE-42 has been specially developed for schools and workshops
- The incident and transmitted illumination unit included as standard can be optionally enabled for the very best illumination of your sample. Mobile use of the OSE 422 is also no problem through the use of rechargeable batteries
- Despite its low price it has very good optical characteristics, which enable you to have sharp images over a large field of view
- · An turnable objective with predefined magnifications is available to make your working procedures quicker and more efficient

- The eyepieces are fixed in the eyepiece tube, to stop them getting damaged or lost
- A special feature of this adaptable and yet robust microscope series is the stable mechanism of the microscope stand which can be adjusted precisely. It will also impress you with its functionality and ergonomic design
- A large selection of eyepieces as well as various additional external illumination units are available as accessories

#### Scope of application

· Training, in vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

· Samples with focus on three-dimensional impression (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Greenough optics
- · Brightness adjustable
- Tube 45° inclined
- Interpupillary distance 55 75 mm
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 200×180×300 mm
- · Net weight approx. 2 kg



Model	Standard configuration						Price
	Tube Eyepiece Field of view Objective Stand Illumination		excl. of VAT ex works				
KERN			mm				€
OSE 421	Binocular	WF 10×/Ø 20 mm	Ø 20	2×/4×	Arm curved	1 W LED (incident); 1 W LED (transmitted)	255,-
OSE 422	Binocular	WF 10×/Ø 20 mm	ø 20	2×/4×	Arm curved	1 W LED (incident); 1 W LED (transmitted) (Rechargeable batteries included)	265,-

Eyepiece	Specifications - Objectives							
	Magnification	2×	4×					
WF 5×	Total magnification	10×	20×					
WF5*	Field of view mm	ø 10	Ø 5					
WE 10	Total magnification	20×	40×					
WF 10×	Field of view mm	ø 10	Ø 5					
WE 15v	Total magnification	30×	60×					
WF 15×	Field of view mm	Ø 7,5	Ø 3,7					
WE 20	Total magnification	40×	80×					
WF 20×	Field of view mm	Ø 6,5	Ø 3,2					
Working distance		57 mm	57 mm					

Model outfit		Mode	el KERN	Order number	Price/piece excl. of VAT	
		OSF 421	OSF 422		ex works €	
	WF 5×/ø 16,2 mm	00	00	OZB-A4101	35,-	
Eyepieces	WF 10×/ø 20 mm	44	11	OZB-A4102	30,-	
(30,5 mm)	WF 15×/ø 15 mm	00	00	OZB-A4103	30,-	
	WF 20×/ø 10 mm	00	00	OZB-A4104	35,-	
Stand	Arm curved, with 1 W LED illumination (transmitted + incident)	<b>✓</b>	✓			
Ctoro minto	Frosted glass/Ø 59,5 mm	✓	✓	OZB-A4815	19,-	
Stage plate	Black-white/Ø 59,5 mm	✓	✓	OZB-A4816	19,-	
External illumination	Please find the information about external illumination units	in the catalogue on p	page 90 and on our v	vebsite www.kern-sc	hn.com	





Stage plate black



Stage plate white

#### **EDUCATIONAL LINE**

## The practical and robust product for schools, training centres, the workshop and laboratory

#### **Features**

- With its integrated handle as well as its stable arm curved stand, the KERN OSF-4G has been specially developed for schools and workshops
- The LED reflected and transmitted illumination included as standard guarantees the very best, continuously dimmable illumination of your sample
- As well as very good optical characteristics, its ergonomic working surface means that it offers the highest level of convenience in this class
- A turnable objective with three predefined magnifications is available to make your working procedures quicker and more effective

- The eyepieces are fixed in the eyepiece tube, to stop them getting damaged or lost
- The ergonomic shape and the stable mechanism which can be adjusted extremely accurately offer a high level of functionality and enable you to work quickly and efficiently with very little effort
- A large selection of eyepieces as well as various additional external illumination units are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

#### Scope of application

 Training, in vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

 Samples with focus on three-dimensional impression (depth, thickness), e.g. insects, seeds, circuit boards, components

#### Technical data

- Optical system: Greenough optics
- Brightness adjustable (separate)
- Tube 45° inclined
- Interpupillary distance 55 75 mm
- Diopter adjustment: One-sided
- Overall dimensions W×D×H 230×180×275 mm
- · Net weight approx. 2,5 kg

#### STANDARD



Model				Standard c	onfiguration		Price
	Tube Eyepiece		Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm				€
OSF 438	Binocular	WF 10×/Ø 20 mm	Ø 20	1×/2×/3×	Arm curved	1 W LED (incident); 0,35 W LED (transmitted)	315,- 🕛
OSF 439	Binocular	WF 10×/Ø 20 mm	Ø 20	1×/2×/4×	Arm curved	1 W LED (incident); 0,35 W LED (transmitted)	315,- 🕛

Eyepiece		S	pecifications - Objective	s	
	Magnification	1×	2×	3×	4×
WE 5	Total magnification	5×	10×	15×	20×
WF 5×	Field of view mm	ø 20	Ø 10	Ø 6,7	Ø 5
WE 40:	Total magnification	10×	20×	30×	40×
WF 10×	Field of view mm	ø 20	Ø 10	Ø 6,7	Ø 5
WE 45	Total magnification	15×	30×	45×	60×
WF 15×	Field of view mm	ø 15	Ø 7,5	Ø 5	Ø 3,7
WE 00:	Total magnification	20×	40×	60×	80×
WF 20×	Field of view mm	ø 10	Ø 6,5	Ø 4,3	Ø 3,2
Working distant	ce	57 mm	57 mm	57 mm	57 mm

Model outfit		Mode	I KERN	Order number	Price/piece excl. of VAT	
		OSF 438	OSF 439		ex works €	
	WF 5×/Ø 16,2 mm	00	00	OZB-A4101	35,-	
Eyepieces	WF 10×/ø 20 mm	<b>*</b>	11	OZB-A4102	30,-	
(30,5 mm)	WF 15×/ø 15 mm	00	00	OZB-A4103	30,-	
	WF 20×/ø 10 mm	00	00	OZB-A4104	35,-	
Stand	Arm curved, incl. handle, with LED illumination (0,35 W transmitted + 1 W incident)	✓	✓			
Ctoro minto	Frosted glass/Ø 59,5 mm	✓	✓	OZB-A4815	19,-	
Stage plate	Black-white / Ø 59,5 mm	✓	✓	OZB-A4816	19,-	
External illumination	Please find the information about external illumination un	nits in the catalogue on p	page 90 and on our	website www.kern-so	ohn.com	



#### LAB LINE

#### The affordable and flexible stereo zoom microscope for laboratories, inspection authorities and quality controls

#### **Features**

- The products in the KERN OZL-44 series are stereo zoom microscopes, which will impress you with their easy handling, flexibility as well as their stability and economical price
- The LED reflected and transmitted illumination included as standard guarantees the very best illumination of your sample
- As well as excellent optical characteristics and their large working surface, these models offer the highest level of comfort in this class - ideal for training companies, workshops as well as assembly and repair workstations, e.g. in the electronics industry
- · The zoom objective gives you continuous magnification of 7,5× - 36×
- The OZL-44 series is available as a binocular version. The eyepieces are fixed in the eyepiece tube, to stop them getting damaged or lost

- The pillar stand offers you the highest level of flexibility and the freedom to remove the microscope head and to integrate it into other modular systems, for example into a universal stand
- · A large selection of eyepieces, external illumination units as well as auxiliary objectives are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- · Please find detailed information in the following model outfit list

#### Scope of application


· In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

· Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Greenough optics
- Brightness adjustable (separate)
- Tube 45° inclined
- Magnification ratio: 4,8:1
- Interpupillary distance 55 75 mm
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 330×235×380 mm
- · Net weight approx. 5 kg

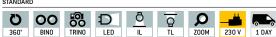

Model		Standard configuration									
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works				
KERN			mm	Zoom			€				
OZL 445	Binocular	WF 10×/Ø 20 mm	Ø 26,7 – 5,6	0,75× - 3,6×	Pillar style	1 W LED (incident); 0,35 W LED (transmitted)	420,-				

OZL 445		S	pecifications - Obje	ctives				
Eyepiece	Magnification	Standard		Auxiliary objectives				
		1,0×	0,5×	0,75×	1,5×	2,0×		
NF 5×	Total magnification	3,75× - 18×	1,875× - 9×	2,81× - 13,5×	5,625× - 27×	7,5×-36×		
WL 2V	Field of view mm	Ø 26 - 6	Ø 60 – 13	Ø 32 - 7	Ø 16-4	ø 12,5 - 3		
NE 40	Total magnification	7,5× - 36×	3,75× - 18×	5,625× - 27×	11,25× - 54×	15×-72×		
WF 10×	Field of view mm	Ø 26,7 – 5,6	ø 53,3 – 11,1	Ø 35,5 – 7,4	Ø 17,8 – 3,7	Ø13,3 - 2,8		
ME 45	Total magnification	11,25× - 54×	5,625× – 27×	8,44× - 40,5×	16,875× – 81×	22,5× - 108×		
WF 15×	Field of view mm	Ø 19-4,5	Ø 43 – 9,5	Ø 24 – 5,5	Ø 12 – 3	Ø 9,5 – 2		
WE 00:	Total magnification	15× - 72×	7,5×-36×	56,25× - 54×	22,5× - 108×	30× - 144×		
WF 20×	Field of view mm	Ø 12,5 - 3	Ø 28 – 6	Ø 16-3,5	Ø 8 – 2	Ø 6 – 1,5		
Working distan	ce	86 mm	178 mm	96 mm	42,5 mm	25,5 mm		
Maximum samp	ole height	100 mm	10 mm	60 mm	120 mm	135 mm		

Model outfit		Model KERN	Order number	Price/piece excl. of VAT						
		OZL 445		ex works €						
	WF 5×/Ø 16,2 mm	00	OZB-A4101	35,-						
Eyepieces	WF 10×/ø 20 mm	44	OZB-A4102	35,-						
(30,5 mm)	WF 15×/ø 15 mm	00	OZB-A4103	30,-						
	WF 20×/ø 10 mm	00	OZB-A4104	30,-						
	0,5×	0	OZB-A4201	75,-						
	0,75×	0	OZB-A4202	75,-						
Auxiliary objectives	1,5×	0	OZB-A4204	75,-						
	2,0×	0	OZB-A4205	80,-						
	Soldering protection lens	0	OZB-A4251	19,-						
Stand	Pillar style, with LED illumination (0,35 W transmitted + 1 W incident)	✓								
a	Frosted glass/Ø 95 mm	✓	OZB-A4805	19,-						
Stage plate	Black-white/ø 95 mm	✓	OZB-A4806	19,-						
External illumination	Please find the information about external illumination units in the catalogue on page 90 and on our website www.kern-sohn.com									







OZL 465 With ring illumination



OZL 467 With handle

#### LAB LINE

### The flexible, affordable all-rounder with zoom function for schools, training companies, inspection authorities and laboratories

#### **Features**

- The products in the KERN OZL-46 series are stereo zoom microscopes, which will impress you with their quality, easy handling, flexibility as well as their stability and economical price
- The LED reflected and transmitted illumination included as standard guarantees the very best illumination of your sample
- The highlight of the OZL 465/OZL 466 is the strong, continuously dimmable, integrated LED ring illumination in the objective housing, which ensures uniform, shadow-free illumination. An LED transmitted light variant is also included
- As well as excellent optical characteristics and their large working surface, these models offer the highest level of comfort in this class – ideal for training companies, workshops as well as assembly and repair workstations, e.g. in the electronics industry
- The zoom objective offers you continuous magnification from 7×-45×

- The KERN OZL-46 series is available as a binocular or trinocular version
- The pillar stand offers you the highest level of flexibility and the freedom to remove the microscope head and to integrate it into other modular systems, for example into a universal stand
- With its integrated handle as well as its stable arm curved stand, the KERN OZL 467/ OZL 468 has been specially developed for schools and workshops
- A large selection of eyepieces, external illumination units as well as auxiliary objectives are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list
- Please find detailed information in the following model outfit list

#### Scope of application

 In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

 Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- Optical system: Greenough optics
- Brightness adjustable (separate)
- $\bullet$  Tube 45° inclined
- Magnification ratio: 6,4:1
- Light distribution 50:50
- Interpupillary distance 55 75 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 300×240×420 mm
- Net weight approx. 4 kg

#### STANDARD



Model				Standard c	onfiguration		Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN		mm Zoom			€		
OZL 463	Binocular	HWF 10×/Ø 20 mm	Ø 28,6 – 4,4	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	500,-
OZL 464	Trinocular	HWF 10×/Ø 20 mm	Ø 28,6 – 4,4	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	570,-
OZL 465	Binocular	HWF 10×/Ø 20 mm	Ø 28,6 - 4,4	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	550,-
OZL 466	Trinocular	HWF 10×/Ø 20 mm	Ø 28,6 - 4,4	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	620,-
OZL 467	Binocular	HWF 10×/Ø 20 mm	Ø 28,6 – 4,4	0,7× - 4,5×	Arm curved	3 W LED (incident); 3 W LED (transmitted)	510,-
OZL 468	Trinocular	HWF 10×/Ø 20 mm	Ø 28,6 - 4,4	0,7×-4,5×	Arm curved	3 W LED (incident); 3 W LED (transmitted)	580,-

Eyepiece		Spec	fications - Objective	es					
	Magnification	Standard		Auxiliary objectives					
		1,0×	0,5×	0,75×	1,5×	2,0×			
HSWF 10×	Total magnification	7×-45×	3,5× - 22,5×	5,3×-33,8×	10,5×-67,5×	14× - 90×			
HOWF IUA	Field of view mm	Ø 28,6 – 4,4	Ø 57,1-8,9	Ø 38,1 – 5,9	Ø 19-3	Ø 14,3 - 2,2			
INA/E 45	Total magnification	10,5× - 67,5×	5,3× - 33,8×	7,9×-50,6×	15,5× - 101,3×	21× - 135×			
HWF 15×	Field of view mm	Ø 21,4 - 3,3	Ø 42,9 - 6,7	Ø 28,5 – 4,4	Ø 14,3 - 2,2	ø 10,7 - 1,7			
HSWF 20×	Total magnification	14× - 90×	7×-45×	10,5× - 67,5×	21×-135×	28× - 180×			
115WF 2U*	Field of view mm	ø 14,3 – 2,2	Ø 28,6 - 4,4	Ø 19,1-2,9	Ø 9,5 - 1,5	Ø 7,1 – 1,1			
INVE 05	Total magnification	17,5× - 122,5×	8,8×-56,3×	13,1×-91,9×	26,3× - 168,8×	35× - 225×			
HWF 25×	Field of view mm	ø 12,9 – 2,0	Ø 25,7 – 4,0	Ø 17,2 - 2,7	Ø 8,6 - 1,3	Ø 6,4 – 1,0			
Working distance		105 mm	177 mm	120 mm	47 mm	26 mm			
Maximum samp	le height	140 mm	35 mm	80 mm	165 mm	185 mm			

Model outfit				Model	KERN			Order number	Price/piece excl. of VAT	
				OZL 465			OZL 468		ex works €	
	HWF 10×/ø 20 mm	11	11	11	11	<b>4</b>	11	OZB-A4631	40,-	
Eyepieces	HSWF 15×/ø 15 mm	00	00	00	00	00	00	OZB-A4632	50,-	
(30,0 mm)	HWF 20×/ø 10 mm	00	00	00	00	00	00	OZB-A4633	50,-	
	HSWF 25×/ø 9 mm	00	00	00	00	00	00	OZB-A4634	55,-	
	0,5×	0	0			0	0	OZB-A4641	80,-	
	0,75×	0	0			0	0	OZB-A4644	80,-	
Auxiliary objectives	1,5×	0	0			0	0	OZB-A4642	80,-	
	2,0×	0	0			0	0	OZB-A4643	80,-	
	Soldering protection lens	0	0	0	0	0	0	OZB-A4645	35,-	
	1× (focus adjustable)		✓		✓		✓	OZB-A4809	50,-	
C-Mount	0,3× (focus adjustable)		0		0		0	OZB-A4810	85,-	
	0,5× (focus adjustable)		0		0		0	OZB-A4811	95,-	
Eyepiece camera adapter	1,0×; for fitting an eyepiece camera to the trinocular connection of the microscope		0		0		0	OZB-A4863	35,-	
	Pillar style, with 3 W-LED illumination (transmitted + incident)	1	✓							
Stand	Pillar style, with 3 W-LED illumination (transmitted)			✓	✓					
	Arm curved, incl. handle, with 3 W-LED illumination (transmitted + incident)					✓	<b>✓</b>			
Ring illumination	Integrated into the microscope head as incident illumination			<b>√</b>	<b>✓</b>					
Stage plate	Frosted glass/Ø 95 mm	✓	1	1	1	✓	<b>✓</b>	OZB-A4670	19,-	
Stage plate	Black-white/Ø 95 mm	1	1	✓	1	✓	✓	OZB-A4806	19,-	



#### LAB LINE

### Stereo zoom microscope with or without halogen illumination, for the laboratory, training centres, quality control or agriculture

#### **Features**

- The KERN OZL-45 stereo zoom microscope series will impress you with its excellent optical characteristics, easy operation and high level of ergonomic working comfort
- The Halogen incident and transmitted illumination included as standard guarantees the very best illumination of your sample
- The high-quality optics, together with a large working surface offers the highest level of comfort for your applications
- The zoom objective offers you continuous magnification from 7,5×-50×
- The KERN OZL-45 series is available as a binocular version
- The pillar stand offers you the highest level of flexibility and the freedom to remove the microscope head and to integrate it into other modular systems, for example into a universal stand

- A large selection of eyepieces, external illumination units as well as auxiliary objectives are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

#### Scope of application

 In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

• Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- Optical system: Greenough optics
- Brightness adjustable
- Tube 45° inclined
- Magnification ratio: 6,7:1
- Interpupillary distance 55 75 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 330×270×460 mm
- · Net weight approx. 5 kg

#### STANDARD



Model				Standard co	onfiguration		Price
	Tube	Eyepiece	Field of view	v Objective Stand Illumination		excl. of VAT ex works	
KERN	mm Zoom				€		
OZL 451	Binocular	HSWF 10×/Ø 23 mm	Ø 33 - 5	0,75×-5,0×	Pillar style	10 W Halogen (incident) 10 W Halogen (transmitted)	610,-

Eyepiece		Spec	ifications - Objectives					
	Magnification	Standard		Auxiliary objectives				
		1,0×	0,5×	0,75×	2,0×			
HWF 5×	Total magnification	3,75× - 25×	1,875× - 12,5×	2,813× - 18,75×	7,5×-50×			
HWF 3^	Field of view mm	Ø 31-4,6	Ø 61,3-9,2	Ø 41,3 - 6,1	Ø 16 – 2,5			
HSWF 10×	Total magnification	7,5× - 50×	3,75× - 25×	5,625× - 37,5×	15× - 100×			
	Field of view mm	Ø 33 - 5	Ø 65 – 10	Ø 44 – 6,7	Ø 16 – 2,5			
UMF 15v	Total magnification	11,25× - 75×	5,625× - 37,5×	8,438× - 56,25×	22,5× – 150×			
HWF 15×	Field of view mm	Ø 24 - 4,2	Ø 48 – 8,5	Ø 32 – 5,6	Ø 12 – 2			
LICWE 204	Total magnification	15× - 100×	7,5×-50×	11,25× – 75×	30× - 200×			
HSWF 20×	Field of view mm	Ø 20 - 3,5	Ø 40 – 7	Ø 26,7 – 4,7	Ø 10 – 1,8			
LIME SEV	Total magnification	18,75× – 125×	9,375× - 62,5×	14,063× - 93,75×	37,5× – 255×			
HWF 25×	Field of view mm	Ø 15,8 - 2,4	Ø 31,5 – 4,8	Ø 24,1-3,2	Ø 7,9 – 1,2			
Working distance	e	113 mm	177 mm	117 mm	35 mm			
Maximum samp	le height	120 mm	60 mm	90 mm	165 mm			

Model outfit		Model KERN	Order number	Price/piece excl. of VAT	
		OZL 451		ex works €	
	HWF 5×/ø 23,2 mm	00	OZB-A4112	40,-	
	HSWF 10×/ø 23 mm	<b>**</b>	OZB-A4118	40,-	
Eyepieces (30,0 mm)	HWF 15×/ø 15 mm	00	OZB-A4119	40,-	
(55,5)	HSWF 20×/ø 14,5 mm	00	OZB-A4120	60,-	
	HWF 25×/ø 11,7 mm	00	OZB-A4121	60,-	
	0,5×	0	OZB-A4209	120,-	
Auxiliary objectives	0,75×	0	OZB-A4210	120,-	
	2,0×	0	OZB-A4206	115,-	
Stand	Pillar style, with 12 V/10 W Halogen Illumination (transmitted + incident)	<b>✓</b>			
01	Frosted glass/Ø 95 mm	✓	OZB-A4805	19,-	
Stage plate	Black-white/Ø 95 mm	✓	OZB-A4806	19,-	
Illumination	10 W spare bulb (transmitted + incident)	<b>✓</b>	OZB-A4804	19,-	
Mechanical stage (Pre-assembling on request)	Stage size W×D 180×155 mm, Travel 75×55 mm, for transmitted and incident illumination	0	OZB-A4605	235,-	
External illumination	Please find the information about external illumination units in the catalogue on page 90 and on our website www.kern-sohn.com				





Dimmable, integrated LED ring illumination

#### LAB LINE

### The practical and flexible stereo zoom microscope with integrated LED ring illumination and large zoom range

#### **Features**

- The KERN OZL-456 stereo zoom microscope series will impress you with its excellent optical characteristics, easy operation and its integrated LED ring illumination unit
- The highlight of the KERN OZL-456 is the strong, continuoulsy dimmable, integrated LED ring illumination in the objective housing, which ensures uniform, shadow-free illumination. An LED transmitted light variant is also included
- With its built-in, top-quality optics and powerful, integrated LED illumination unit, this model is a special all-rounder for all areas of application
- The zoom objective offers you continuous magnification from 7,5×  $50\times$

- As standard, the KERN OZL-45R series is provided as a binocular version with 10× eyepieces with a field of view with a diameter of 23 mm
- The arm curved stand gives you a large working area as well as a precise adjustment mechanism
- A large selection of eyepieces as well as auxiliary objectives are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- Please find detailed information in the following model outfit list

#### Scope of application

 In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control

#### Applications/Samples

• Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- Optical system: Greenough optics
- · Incident illumination dimmable
- $\bullet$  Tube 45° inclined
- Magnification ratio: 6,7:1
- Interpupillary distance 55 75 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 320×275×420 mm
- Net weight approx. 4,5 kg

#### STANDARI



Model				Standard co	onfiguration		Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZL 456	Binocular	HSWF 10×/Ø 23 mm	Ø 33 – 5	0,75×-5,0×	Arm curved	1 W LED (incident); 0,21 W LED (transmitted)	690,-

Eyepiece		Specifications - Objectives						
	Magnification	Standard	Auxiliary objectives					
		1,0×	2,0×					
HWF 5×	Total magnification	3,75× - 25×	7,5×-50×					
IIWI J^	Field of view mm	Ø 31-4,6	Ø 16 – 2,5					
HSWF 10×	Total magnification	7,5×-50×	15× – 100×					
	Field of view mm	Ø 33-5	Ø 16-2,5					
HWF 15×	Total magnification	11,25× - 75×	22,5× - 150×					
HWF 15*	Field of view mm	Ø 24 - 4,2	Ø 12 - 2					
HCME 204	Total magnification	15× - 100×	30× - 200×					
HSWF 20×	Field of view mm	Ø 20 - 3,5	Ø 10 - 1,8					
LINE OF	Total magnification	18,75× – 125×	37,5×-255×					
HWF 25×	Field of view mm	Ø 15,8 – 2,4	Ø 7,9 – 1,2					
Working distant	ce	113 mm	35 mm					
Maximum samp	Maximum sample height		95 mm					

Model outfit		Model KERN	Order number	Price/piece excl. of VAT
		OZL 456		ex works €
	HWF 5×/ø 23,2 mm	00	OZB-A4112	40,-
	HSWF 10×/ø 23 mm	44	OZB-A4118	40,-
Eyepieces (30,0 mm)	HWF 15×/ø 15 mm	00	OZB-A4119	40,-
(22,2,	HSWF 20×/ø 14,5 mm	00	OZB-A4120	60,-
	HWF 25×/ø 11,7 mm	00	OZB-A4121	60,-
Stand	Arm curved, with LED illumination (0,21 W transmitted + 1 W incident)	<b>√</b>		
01	Frosted glass/Ø 95 mm	✓	OZB-A4805	19,-
Stage plate	Black-white/Ø 95 mm	✓	OZB-A4806	19,-
Mechanical stage (Pre-assembling on request)	Stage size W×D 180×155 mm, Travel 75×55 mm, for transmitted and incident illumination	0	OZB-A4605	235,-
External illumination	Please find the information about external illumination units	in the catalogue on page 90 and on o	our website www.kern-sc	hn.com



#### LAB LINE

### First-class optics and strong illumination combined with a high level of flexibility

#### **Features**

- The KERN OZM series is a range of excellent stereo zoom microscopes with above-average optical features
- · The ergonomic shape allows a simple, effortless working over a period of several hours
- The extraordinarily strong and continuously dimmable 3 W LED reflected and transmitted illumination ensures a flexible and particularly good level of illumination for your sample
- · With its large working distance, an extra large field of view and its brilliant resolution, the KERN OZM provides sharp, high-contrast, colour-true images
- The zoom objective gives you continuous magnification from 7,5×-45×
- · There is a choice of a binocular model as well as a trinocular model for connecting a camera for documentation purposes and for quality reports

- The pillar stand is particularly flexible due to its variable and sturdy adjustment mechanism and therefore enables ergonomic working procedures
- A large selection of eyepieces, (universal) stands, a darkfield kit, external illumination units as well as auxiliary objectives and more are available as accessories
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model
- · Please find detailed information in the following model outfit list

#### Scope of application


· In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control, electronics and semiconductor industry, assembly and repair

#### Applications/Samples

· Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Greenough optics
- Brightness adjustable (separate)
- Tube 45° inclined
- Magnification ratio: 6,4:1
- Light distribution OZM 543/544: 50:50
- Interpupillary distance 52 76 mm
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 330×285×440 mm
- · Net weight approx. 4,5 kg

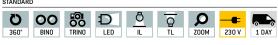

OPTION	
min	
SCALE	

Model				Standard c	onfiguration		Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZM 542	Binocular	HSWF 10×/Ø 23 mm	Ø 32,8 – 5,1	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	1070,-
OZM 544	Trinocular	HSWF 10×/Ø 23 mm	Ø 32,8 – 5,1	0,7× - 4,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	1350,-

Eyepiece		Specifications - Objectives									
	Magnification	Standard	Auxiliary objectives								
		1,0×	0,5×	0,7×	1,5×	2×					
HSWF 10×	Total magnification	7×-45×	3,5×-22,5×	4,9×-31,5×	10,5×-67,5×	14×-90×					
nswr Iu^	Field of view mm	Ø 32,8 - 5,1	ø 65,7 – 10,2	Ø 46,9 – 7,3	Ø 21,9 - 3,4	Ø 16,4 - 2,6					
SWF 15×	Total magnification	10,5× - 67,5×	5,3× - 33,8×	7,4× - 47,2×	15,8× - 101,3×	21× - 135×					
	Field of view mm	Ø 24,3 - 3,8	Ø 48,6 – 7,6	Ø 34,7 – 5,4	Ø 16,2 - 2,5	Ø 12,1 – 1,9					
014/E 00	Total magnification	14×-90×	7×-45×	9,8×-63×	21× - 135×	28× - 180×					
SWF 20×	Field of view mm	Ø 20 - 3,1	Ø 40 - 6,2	Ø 28,6 – 4,4	Ø 13,3 - 2,1	Ø 10 – 1,6					
OUE 00.	Total magnification	21× - 135×	10,5× - 67,5×	14,7× - 94,5×	31,5× - 202,5×	42× - 270×					
SWF 30×	Field of view mm	Ø 12,9 - 2	Ø 25,7 – 4	Ø 18,4 – 2,9	Ø 8,6 – 1,6	Ø 6,4 – 1					
Working distance		110 mm	195 mm	145 mm	50 mm	35 mm					
Maximum sample height		130 mm	30 mm	65 mm	160 mm	175 mm					

Model outfit		Model KERN		Order number	Price/piece excl. of VAT
		OZM 542	OZM 544		ex works €
	HSWF 10×/ø 23 mm	<b>√</b> √	44	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	00	OZB-A5504	60,-
	SWF 20×/Ø 14 mm	00	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/Ø 9 mm	00	00	OZB-A5506	85,-
,	HSWF 10×/ø 23 mm (reticule 0,1 mm)	0	0	OZB-A5512	95,-
	SWF 15×/Ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-
	SWF 20×/Ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-
	0,5×	0	0	OZB-A5612	135,-
	0,7×	0	0	OZB-A5613	135,-
Achromatic auxiliary objectives	1,5×	0	0	OZB-A5615	135,-
	2,0×	0	0	OZB-A5616	140,-
	Soldering protection lens	0	0	OZB-A5614	45,-
	0,3× (focus adjustable)		0	OZB-A5701	120,-
	0,5× (focus adjustable)		0	OZB-A5702	120,-
	1,0× (focus adjustable)		0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703		0	OZB-A5704	255,-
	for SLR cameras (Nikon)		0	OZB-A5706	255,-
	for SLR cameras (Olympus)		0	OZB-A5707	255,-
	for SLR cameras (Canon)		0	OZB-A5708	255,-
Darkfield unit	Darkfield unit	0	0	OZB-A4601	70,-
Object clamp	Object clamp	0	0	OBB-A6205	35,-
	Pillar style, without illumination				
Stand	Pillar style, with 3 W LED illumination (transmitted + incident)	<b>~</b>	<b>✓</b>		
	Please find more stands in the catalogue on page 86 and on ou	ır website www.ke	rn-sohn.com		
	Frosted glass/Ø 94,5 mm	✓	✓	OZB-A5192	19,-
Stage plate	Black-white/Ø 94,5 mm	✓	✓	OZB-A5191	19,-
	Glass/ø 94,5 mm	0	0	OZB-A5190	19,-
Mechanical stage	Stage size W×D 188×160 mm, Travel 76×65 mm, for transmitted and incident illumination	0	0	OZB-A5781	175,-
(Pre-assembling on request)	Stage size W×D 180×175 mm, Travel 100×86 mm, for incident illumination only	0	0	OZB-A5782	240,-
External illumination	Please find the information about external illumination units in	the catalogue on	page 90 and on ou	r website www.kern-so	hn.com
			✓ = Included v	Sels at a Processing	O = Ontio



#### LAB LINE

### Professional and powerful – thanks to its extremely large magnification range, strong illumination and first-class optics

#### **Features**

- The KERN OZP stereo zoom microscope stands out through its above-average magnification range and its robust shape which is also ergonomic, it enables effortless, simple working over a period of several hours
- The KERN OZP series is available as a strong, continuously adjustable 3 W LED reflected and transmitted light variant for the very best illumination of your sample or as a variant without illumination
- · With its large working distance, an extra large field of view and brilliant resolution, the KERN OZP provides sharp, high-contrast and colour-true images
- · The extremely large, continuously adjustable magnification range from 6 to 55 times magnification means that you can work quickly and effectively
- · There is a choice of a binocular model as well as a trinocular model for connecting a camera for documentation purposes and for quality reports

- The pillar stand is particularly flexible due to its variable and sturdy adjustment mechanism and therefore enables ergonomic working procedures
- A large selection of eyepieces, (universal) stands, a darkfield kit, external illumination units as well as auxiliary objectives and more are available as accessories
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model
- · Please find detailed information in the following model outfit list

#### Scope of application

 Zoology and botany, quality control, electronics and semiconductor industry, assembly and repair

#### Applications/Samples

· Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Greenough optics
- Brightness adjustable (separate)
- · Tube: 35° inclined
- Magnification ratio: 9,2:1
- Light distribution OZP 557/558: 50:50
- Interpupillary distance 52 76 mm
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 330×285×470 mm
- · Net weight approx. 4,5 kg



OPTION
www
SCALE

Model				Standard c	onfiguration		Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZP 556	Binocular	HSWF 10×/Ø 23 mm	Ø 38,3 – 4,2	0,6× - 5,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	1230,-
OZP 558	Trinocular	HSWF 10×/Ø 23 mm	Ø 38,3 – 4,2	0,6×-5,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	1530,-

Eyepiece		Specifications – Objectives									
	Magnification	Standard		Auxiliary objectives							
		1,0×	0,5×	0,7×	1,5×	2×					
HSWF 10×	Total magnification	6×-55×	3×-27,5×	4,2×-38,5×	9×-82,5×	12× - 110×					
nswr iu×	Field of view mm	Ø 38,3 - 4,2	Ø 76,7 – 8,4	Ø 54,8 – 6	Ø 25,6 – 2,8	Ø 19,2 - 2,1					
SWF 15×	Total magnification	9× - 82,5×	4,5×-41,25×	6,3× – 57,75×	13,5× - 123,75×	18× - 165×					
	Field of view mm	Ø 28,3 - 3,1	Ø 56,7 – 6,2	Ø 40,5 – 4,4	Ø 18,9 – 2,1	Ø 14,2 - 1,5					
2WE 204	Total magnification	12× - 110×	6×-55×	8,4× – 77×	18× - 165×	24× - 220×					
SWF 20×	Field of view mm	Ø 23,3 – 2,5	Ø 46,7 – 5,1	Ø 33,3 – 3,6	Ø 15,6 – 1,7	ø 11,7 - 1,3					
OWE 00	Total magnification	18× – 165×	9×-82,5×	12,6× - 115,5×	27× - 247,5×	36×-330×					
SWF 30×	Field of view mm	ø 15 – 1,6	Ø 30 – 3,3	Ø 21,4 - 2,3	Ø 10 – 1,1	Ø 7,5 - 0,8					
Working distance		108 mm	195 mm	145 mm	50 mm	35 mm					
Maximum sample height		110 mm	10 mm	45 mm	140 mm	150 mm					

Model outfit		Mode	el KERN	Order number	Price/piece excl. of VAT				
		OZP 556	OZP 558		ex works €				
	HSWF 10×/ø 23 mm	44	44	OZB-A5503	60,-				
	SWF 15×/ø 17 mm	00	00	OZB-A5504	60,-				
	SWF 20×/ø 14 mm	00	00	OZB-A5505	70,-				
<b>Eyepieces</b> (30,0 mm)	SWF 30×/ø 9 mm	00	00	OZB-A5506	85,-				
,	HSWF 10×/Ø 23 mm (reticule 0,1 mm)	0	0	OZB-A5512	95,-				
	SWF 15×/ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-				
	SWF 20×/Ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-				
	0,5×	0	0	OZB-A5612	135,-				
	0,7×	0	0	OZB-A5613	135,-				
Achromatic auxiliary objectives	1,5×	0	0	OZB-A5615	135,-				
	2,0×	0	0	OZB-A5616	140,-				
	Soldering protection lens	0	0	OZB-A5614	45,-				
	0,3× (focus adjustable)		0	OZB-A5701	120,-				
	0,5× (focus adjustable)		0	OZB-A5702	120,-				
	1,0× (focus adjustable)		0	OZB-A5703	95,-				
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703		0	OZB-A5704	255,-				
	for SLR cameras (Nikon)		0	OZB-A5706	255,-				
	for SLR cameras (Olympus)		0	OZB-A5707	255,-				
	for SLR cameras (Canon)		0	OZB-A5708	255,-				
Darkfield unit	Darkfield unit	0	0	OZB-A4601	70,-				
Object clamp	Object clamp	0	0	OBB-A6205	35,-				
	Pillar style, without illumination								
Stand	Pillar style, with 3 W LED illumination (transmitted + incident)	<b>~</b>	<b>✓</b>						
	Please find more stands in the catalogue on page 86 and on our website www.kern-sohn.com								
	Frosted glass/Ø 94,5 mm		✓	OZB-A5192	19,-				
Stage plate	Black-white/Ø 94,5 mm	✓	✓	OZB-A5191	19,-				
	Glass/Ø 94,5 mm		0	OZB-A5190	19,-				
Mechanical stage	Stage size W×D 188×160 mm, Travel 76×65 mm, for incident and transmitted illumination	0	0	OZB-A5781	175,-				
(Pre-assembling on request)	Stage size W×D 180×175 mm, Travel 100×86 mm, for incident illumination only	0	0	OZB-A5782	240,-				
External illumination	Please find the information about external illumination units in	the catalogue on p	page 90 and on ou	r website www.kern-so	hn.com				
			✓ = Included v	24. 1. 2	0 = 0				







**OZR 563** Without illumination

#### **PROFESSIONAL LINE**

### Professional stereo zoom microscope with parallel optics for excellent images, depth of field, contrast and fatigue-free working

#### **Features**

- The KERN OZR series is a special, high-quality stereo zoom microscope with parallel optics for demanding analyses
- · The KERN OZR series is available as a strong, continuously adjustable 3 W LED reflected and transmitted light variant for the very best illumination of your sample or as a variant without illumination
- · The parallel optical system is a high-quality optical system and provides excellent images with the best contrast, colour and depth of field with fatigue-free working. Refocusing is also only necessary in very few cases when magnifying the zoom
- · The continuously adjustable magnification range from 8 to 50 times magnification means that you can work quickly and effectively
- · As standard, the models of the KERN OZR series are trinocular and are therefore equipped for connecting a camera for documentation purposes and for quality reports

- The pillar stand is particularly flexible due to its variable and sturdy adjustment mechanism and therefore enables ergonomic working procedures
- A large selection of eyepieces, (universal) stands, a darkfield kit, external illumination units as well as auxiliary objectives and more are available as accessories
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model
- · Please find detailed information in the following model outfit list

#### Scope of application

· In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control, electronics and semiconductor industry, assembly and repair

#### Applications/Samples

· Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Parallel optics
- Brightness adjustable (separate)
- Tube 45° inclined
- · Magnification ratio: 6,25:1
- Light distribution 50:50
- Interpupillary distance 52 76 mm
- · Diopter adjustment: Both-sided
- · Overall dimensions W×D×H 305×300×540 mm
- · Net weight approx. 5,5 kg

#### STANDARD



OPTION
min
SCALE

Model		Standard configuration					Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZR 563	Trinocular	HWF 10×/Ø 22 mm	Ø 27,5 – 4,4	0,8× - 5×	Pillar style	-	1550,-
OZR 564	Trinocular	HWF 10×/Ø 22 mm	Ø 27,5 - 4,4	0,8×-5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	1750,-

Eyepiece		Speci	fications - Objectives				
	Magnification	Standard Plan		Achromatic objectives			
		1,0×	0,5×	0,7×	1,5× (Auxiliary)		
LIME 100	Total magnification	8×-50×	4× - 25×	5,6×-35×	12×-75×		
HWF 10×	Field of view mm	Ø 27,5 – 4,4	Ø 55 – 8,8	ø 39,3 - 6,3	ø 18,33 – 2,93		
SWF 15×	Total magnification	12×-75×	6×-37,5×	8,4× - 5,5×	18× - 112,5×		
	Field of view mm	Ø 21,25 - 3,4	Ø 42,5 - 6,8	Ø 30,36 – 4,86	Ø 14,17 – 2,27		
NA/E 00	Total magnification	16× – 100×	8×-50×	11,2×-70×	24× - 150×		
SWF 20×	Field of view mm	Ø 17,5 - 2,8	Ø 35 – 5,6	Ø 25 - 4	Ø 11,67 – 1,87		
NA/E 00	Total magnification	24× – 150×	12× - 75×	16,8× - 105×	36× – 225×		
SWF 30×	Field of view mm	Ø 11,25 – 1,8	Ø 22,5 - 3,6	Ø 16,1 – 2,57	Ø 7,5 – 1,2		
Norking distand	ce	91 mm	186 mm	135 mm	40 mm		
Maximum samp	le height	100 mm	30 mm	80 mm	125 mm		

Model outfit		Mode	I KERN	Order number	Price/piece
		OZR 563	OZR 564		ex works €
	HWF 10×/Ø 22 mm	<b>√</b> √	44	OZB-A5502	60,-
	SWF 15×/Ø 17 mm	00	00	OZB-A5504	60,-
Eyepieces (30,0 mm)	SWF 20×/Ø 14 mm	00	00	OZB-A5505	70,-
	SWF 30×/ø 9 mm	00	00	OZB-A5506	85,-
,00,0,	HWF 10×/Ø 22 mm (reticule 0,1 mm)	0	0	OZB-A5511	110,-
	SWF 15×/Ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-
	SWF 20×/Ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-
Plan achromatic objective	1,0×	✓	✓	OZB-A5603	490,-
	0,5×	0	0	OZB-A5601	205,-
Achromatic objectives	0,7×	0	0	OZB-A5602	205,-
-2,00000	1,5× Only in combination with OZB-A5603	0	0	OZB-A5604	265,-
Trinocular Division 100:0		✓	✓	OZB-A5401	575,-
beamsplitter Division 50:50		0	0	OZB-A5402	575,-
C-Mount	0,3× (focus adjustable)	0	0	OZB-A5701	120,-
	0,5× (focus adjustable)	0	0	OZB-A5702	120,-
	1,0× (focus adjustable)	0	0	OZB-A5703	95,-
	1,0× (with micrometer) only in combination with OZB-A5703	0	0	OZB-A5704	255,-
	for SLR cameras (Nikon)	0	0	OZB-A5706	255,-
	for SLR cameras (Olympus)	0	0	OZB-A5707	255,-
	for SLR cameras (Canon)	0	0	OZB-A5708	255,-
Darkfield unit	Darkfield unit		0	OZB-A4601	70,-
Object clamp	Object clamp	0	0	OBB-A6205	35,-
	Pillar style, without illumination	✓			
Stand	Pillar style, with 3 W LED illumination (transmitted + incident)		✓		
	Frosted glass/Ø 94,5 mm		✓	OZB-A5192	19,-
Stage plate	Black-white/Ø 94,5 mm	✓	✓	OZB-A5191	19,-
	Glass/Ø 94,5 mm		0	OZB-A5190	19,-
Mechanical stage Pre-assembling on	Stage size W×D 188×160 mm, Travel 76×65 mm, for incident and transmitted illumination	0	0	OZB-A5781	175,-
Pre-assembling on request)	Stage size W×D 180×175 mm, Travel 100×86 mm, for incident illumination only	0	0	OZB-A5782	240,-
External Illumination	Please find the information about external illumination units in	the catalogue on p	age 90 and on our  ✓ = Included w		ohn.com







OZS 573
Without illumination

#### **PROFESSIONAL LINE**

### High-quality parallel optics for the best images, depth of field and contrast – with extra large zoom range

#### **Features**

- The devices in the KERN OZS series are special, high-quality stereo zoom microscopes with parallel optics and an above-average zoom range for demanding analyses
- The KERN OZS series is available as a strong, continuously adjustable 3 W LED reflected and transmitted light variant for the very best illumination of your sample or as a variant without illumination
- The parallel optics is a high-quality optical system and provides excellent images with the best contrast, colour and depth of field with fatigue-free working
- The exceptionally large, adjustable magnification range from 8 to 80 times magnification gives you continuous zoom for the views you need
- As standard, the models of the KERN OZS series are trinocular and therefore equipped for connecting a camera for documentation purposes and for quality reports

- These models are particularly flexible due to the modular construction, the pillar stand as well as the powerful, longlife LED reflected and transmitted illumination which can be integrated as an option
- A large selection of eyepieces, (universal) stands, a darkfield kit, external illumination units as well as auxiliary objectives and more are available as accessories
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- Please find detailed information in the following model outfit list

#### Scope of application

 In vitro fertilisation, detection of parasites, zoology and botany, tissue preparation, section, quality control, electronics and semiconductor industry, assembly and repair

#### Applications/Samples

 Samples with focus on three-dimesnional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

#### **Technical data**

- · Optical system: Parallel optics
- Brightness adjustable (separate)
- Tube 45° inclined
- Magnification ratio: 10:1
- Light distribution 50:50
- Interpupillary distance 52 76 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 305×300×540 mm
- Net weight approx. 5,5 kg

## STANDARD No. Standard Standard D</t

Model		Standard configuration					Price
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZS 573*	Trinocular	HWF 10×/Ø 22 mm	Ø 27,5 – 2,75	0,8× - 8×	Pillar style	-	1950,-
OZS 574	Trinocular	HWF 10×/Ø 22 mm	Ø 27,5 – 2,75	0,8×-8×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	2980,-

Eyepiece		Specifications - Objectives							
	Magnification	Standard Plan		Achromatic objectives					
		1,0×	0,5×	0,7×	1,5× (Auxiliary)				
INVE 10v	Total magnification	8×-80×	4×-40×	5,6× - 56×	12× - 120×				
HWF 10×	Field of view mm	ø 27,5 – 2,75	Ø 55 – 5,5	Ø 39,3 - 3,93	ø 18,33 – 1,83				
SWF 15×	Total magnification	12× - 120×	6×-60×	8,4× - 84×	18× – 180×				
	Field of view mm	Ø 21,25 - 2,13	Ø 42,5 - 4,25	Ø 30,36 – 3,04	Ø 14,17 - 1,42				
DWE 00	Total magnification	16× - 160×	8×-80×	11,2× - 112×	24× - 240×				
SWF 20×	Field of view mm	ø 17,5 – 1,75	Ø 35 – 3,5	Ø 25 - 2,5	Ø 11,67 – 1,17				
014/E 00	Total magnification	24× - 240×	12× – 120×	16,8× - 168×	36×-360×				
SWF 30×	Field of view mm	ø 11,25 - 1,13	Ø 22,5 – 2,25	ø 16,1 – 1,61	Ø 7,5 – 0,75				
Working distand	ce	91 mm	186 mm	135 mm	40 mm				
Maximum samp	Maximum sample height		30 mm	80 mm	125 mm				

Model outfit		Mode	I KERN	Order number	Price/piece
		OZS 573	OZS 574		ex works €
	HWF 10×/ø 22 mm	<b>√</b> √	44	OZB-A5502	60,-
	SWF 15×/Ø 17 mm	00	00	OZB-A5504	60,-
	SWF 20×/Ø 14 mm	00	00	OZB-A5505	70,-
E <b>yepieces</b> (30,0 mm)	SWF 30×/Ø 9 mm	00	00	OZB-A5506	85,-
	HWF 10×/Ø 22 mm (reticule 0,1 mm)	0	0	OZB-A5511	110,-
	SWF 15×/ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-
	SWF 20×/ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-
Plan achromatic objective	1,0×	<b>✓</b>	✓	OZB-A5603	490,-
	0,5×	0	0	OZB-A5601	205,-
Achromatic objectives	0,7×	0	0	OZB-A5602	205,-
,	1,5× Only in combination with OZB-A5603	0	0	OZB-A5604	265,-
Trinocular Division 100:0		✓	✓	OZB-A5401	575,-
peamsplitter	Division 50:50	0	0	OZB-A5402	575,-
	0,3× (focus adjustable)	0	0	OZB-A5701	120,-
	0,5× (focus adjustable)	0	0	OZB-A5702	120,-
	1,0× (focus adjustable)	0	0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703	0	0	OZB-A5704	255,-
	for SLR cameras (Nikon)	0	0	OZB-A5706	255,-
	for SLR cameras (Olympus)	0	0	OZB-A5707	255,-
	for SLR cameras (Canon)	0	0	OZB-A5708	255,-
Darkfield unit	Darkfield unit		0	OZB-A4601	70,-
Object clamp	Object clamp	0	0	OBB-A6205	35,-
	Pillar style, without illumination	✓			
Stand	Pillar style, with 3 W LED illumination (transmitted + incident)		✓		
	Frosted glass / Ø 94,5 mm		✓	OZB-A5192	19,-
Stage plate	Black-white/Ø 94,5 mm	✓	✓	OZB-A5191	19,-
	Glass/Ø 94,5 mm		0	OZB-A5190	19,-
Mechanical stage (Pre-assembling on request)	Stage size W×D 188×160 mm, Travel 76×65 mm, for incident and transmitted illumination	0	0	OZB-A5781	175,-
	Stage size W×D 180×175 mm, Travel 100×86 mm, for incident illumination only	0	0	OZB-A5782	240,-
External Ilumination	Please find the information about external illumination units in	the catalogue on p	age 90 and on our	website www.kern-so	hn.com
			✓ = Included w	ith dolivory	O = Op





Plug in for power supply

#### **PROFESSIONAL LINE**

### The coaxial with parallel optics for excellent contrast and depth of field

#### **Features**

- · The KERN OZC has been developed specially to meet requirements for high contrast and depth of field. These devices are absolutely essential for the LCD/LED electronics industry
- The coaxial 2 W LED reflected illumination which is integrated into the objective guarantees selective depth of focus, so that even low-lying sections can be recorded (e.g. the bottom of a drilled hole)
- The parallel optics is a high-quality optical system and provides excellent images with the best contrast, colour and depth of field with fatigue-free working. Refocusing is also only necessary in very few cases when magnifying the zoom
- · The large, adjustable magnification range from 18 to 65 times gives you continuous zoom when you are working

- · As standard, the KERN OZC is trinocular and is therefore equipped for connecting a camera for documentation purposes and for quality reports
- The arm curved stand ensures precise adjustment and focusing of your sample. The stand base is particularly heavy and therefore offers a high level of stability and an extremely secure footing
- · A large selection of eyepieces and a mechanical stage extension are available as accessories
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- A C-mount adapter is required to connect a camera. You can select this adapter from the following model outfit list
- · Please find detailed information in the following model outfit list

#### Scope of application

· LCD/LED electronics, semiconductor technology

#### Applications/Samples

• Samples with focus on three-dimesnional impression (depth, thickness), zoom for variable magnification, e.g. LCD/LED electronics, circuit boards, ICs

#### **Technical data**

- · Optical system: Parallel optics
- · Brightness adjustable
- Tube 45° inclined
- · Magnification ratio: 3,6:1
- Light distribution 50:50
- Interpupillary distance 52 76 mm
- · Diopter adjustment: Both-sided Overall dimensions W×D×H
- 305×180×405 mm
- · Net weight approx. 6,6 kg.



















with	
SCALE	

Model		Standard configuration			Price		
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZC 583	Trinocular	HSWF 10×/Ø 23 mm	Ø 12,78 - 3,5	1,8× - 6,5×	Arm curved	2 W LED (coaxial incident)	2990,-

Stereomicroscopes

Eyepiece	Specifications - Objectives			
	Magnification	Standard		
		1,0×		
HWF 10×	Total magnification	18×-65×		
HWF 10*	Field of view mm	Ø 12,78 - 3,5		
SWF 15×	Total magnification	27×-97,5×		
SWF 15*	Field of view mm	Ø 9,5 – 2,6		
SWF 20×	Total magnification	36× - 130×		
3WF 20^	Field of view mm	Ø 7,78- 2,2		
SWF 30×	Total magnification	54× - 195×		
3WF 3U^	Field of view mm	Ø 5 – 1,4		
Working distance		92 mm		
Maximum sample he	ight	35 mm		

Model outfit		Model KERN	Order number	Price/piece excl. of VAT
		OZC 583		ex works €
	HSWF 10×/ø 23 mm	<b>√</b> √	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	OZB-A5504	60,-
	SWF 20×/ø 14 mm	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/ø 9 mm	00	OZB-A5506	85,-
(,,	HSWF 10×/ø 23 mm (reticule 0,1 mm)	0	OZB-A5512	95,-
	SWF 15×/ø 17 mm (reticule 0,05 mm)	0	OZB-A5513	110,-
	SWF 20×/ø 14 mm (reticule 0,05 mm)	0	OZB-A5514	110,-
	0,3× (focus adjustable)	0	OZB-A5701	120,-
	0,5× (focus adjustable)	0	OZB-A5702	120,-
	1,0× (focus adjustable)	0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703	0	OZB-A5704	255,-
	for SLR cameras (Nikon)	0	OZB-A5706	255,-
	for SLR cameras (Olympus)	0	OZB-A5707	255,-
	for SLR cameras (Canon)	0	OZB-A5708	255,-
Stand	Arm curved, without illumination	✓		
External illumination	Please find the information about external illumination units in the o	catalogue on page 90 and on	our website www.kern-sc	hn.com





Side view

#### LAB LINE

#### The specialist for jewellers and the gem industry

#### **Features**

- The KERN OZG series has been specially developed for jewellers and mineral observations in the gem industry. Precious stones and gems can be checked and handled with this stereo zoom microscope
- You have a choice of a strong halogen transmitted illumination unit as well as halogen reflected and transmitted illumination variants, each with an additional frontal illumination
- · As well as very good optical characteristics, this model forms an ideal package with its dark field unit with object clamp which is included in the scope of delivery
- The KERN OZG 493 is fitted with a pole stand which has both integrated bright halogen light units with incident and transmitted illumination, as well as additional front lighting
- · A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery
- · Please find detailed information in the following model outfit list

#### Scope of application

· Jewellers and gem industry

#### Applications/Samples

• Samples with focus on three-dimesnional impression (depth, thickness), zoom for variable magnification, special stand for processing workpieces e.g. gems, components, precious stones

#### **Technical data**

- Optical system: Greenough optics
- Brightness adjustable
- Tube 45° inclined
- Interpupillary distance 55 75 mm
- Diopter adjustment: Both-sided
- Magnification ratio: 5,1:1
- · Overall dimensions W×D×H 310×170×350 mm
- · Net weight approx. 5 kg

#### STANDARD

















Ī	
	1 DAY

Model	Standard configuration						
	Tube	Eyepiece	Field of view	Objective	Stand	Illumination	excl. of VAT ex works
KERN			mm	Zoom			€
OZG 493	Binocular	WF 10×/Ø 20 mm	Ø 26,7 – 5,6	0,7× - 3,6×	Pillar style	10 W Halogen (incident) 10 W Halogen (transmitted) 10 W Fluorescence (front illumination)	495,-

OZG 493	Specifications - Objectives					
Eyepiece	Magnification	Standard				
		1,0×				
WF 5×	Total magnification	3,75× - 18×				
WF 3^	Field of view mm	Ø 26 - 6				
WF 10×	Total magnification	7,5×-36×				
WF IU×	Field of view mm	Ø 26,7 – 5,6				
WF 15×	Total magnification	11,25× - 54×				
WF 15*	Field of view mm	Ø 19 – 4,5				
WF 20×	Total magnification	15×-72×				
WF 2U^	Field of view mm	Ø 12,5 - 3				
Working distance	86 mm					

Model outfit		Model KERN	Order number	Price/piece excl. of VAT ex works €
		OZG 493		
	WF 5×/Ø 16,2 mm	00	OZB-A4101	35,-
Eyepieces	WF 10×/Ø 20 mm	44	OZB-A4102	35,-
(30,5 mm)	WF 15×/ø 15 mm	00	OZB-A4103	30,-
	WF 20×/Ø 10 mm	00	OZB-A4104	30,-
Darkfield unit	Darkfield unit	✓	OZB-A4601	70,-
Object clamp	Object clamp (steel wire)	✓	OZB-A4604	30,-
Stand	Pillar style, with 12 V/10 W Halogen (transmitted + incident) and 10 W Fluorescent illumination (front)	✓		
Chara alaka	Frosted glass/Ø 95 mm	✓	OZB-A4805	19,-
Stage plate	Black-white / Ø 95 mm ✓	OZB-A4806	19,-	
Illumination	10 W spare bulb (transmitted + incident)	✓	OZB-A4804	19,-

# 5 Digital microscope sets





Our all-round compound microscope as a comprehensive digital solution for schools, training and laboratories

#### **Features**

- Laboratory microscopes from the OBE and OBF range are now also available to you as a comprehensive digital solution for your live investigations. Optionally available with an mounted tablet or C-mount camera. Naturally, the appropriate C-mount adapter is included with the delivery
- The mounted KERN ODC 241 tablet-camera has been specially developed for simple, convenient and direct investigation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- The mounted C-mount camera is available in different versions and can be used anywhere

- For detailed information on the individual components, see the relevant product description of the individual item
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery

#### **Technical data**

- Finite optical system
- · Quadplex nosepiece
- Siedentopf 30° inclined/360° rotatable
- · Diopter adjustment: One-sided
- Eyepiece: HWF 10×/Ø 18 mm

#### OBE-1

 Overall dimensions W×D×H 320×180×365 mm

- Net weight approx. 5,5 kg
- · Objective quality: Achromatic
- Objectives OBE 104: 4×/10×/40×
- Objectives OBE 114: 4×/10×/40×/100×
- Illumination: 3 W LED (transmitted)

#### OBF-1

- Overall dimensions W×D×H 395×200×380 mm
- Net weight approx. 7,7 kg
- · Objective quality OBF 131: Achromatic
- Objective quality OBF 132/133: Plan
- Objectives: 4×/10×/40×/100×
- Illumination OBF 131/132:
   20 W Halogen (transmitted)
- · Illumination OBF 133: 3 W LED (transmitted)

Model			Standard configuration (camer	ra)		Price excl. of VAT
KERN	Included camera	Resolution camera	Interface	Sensor	Details microscope, camera	ex works €
OBE 104C825	ODC 825	5 MP	USB 2.0 (6,8 – 55 FPS)	CMOS 1/2,5"		1085,-
OBE 114C825	000 823	JIVIF	030 2.0 (0,0 - 33 1 + 3)	GIVIO3 1/2,3	KERN Optics catalogue Page 13, 92	1145,-
OBE 104C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"		1230,-
OBE 114C832	000 832	3 IVIP	USB 3.0 (14,2 = 101,2 FP3)	CIVIOS 1/2,5		1290,-
OBE 104T241	ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue	1590,-
OBE 114T241	ODC 241	3 IVIP	WLAN, 03B 2.0, HDWI, 3D (15-30 FF3)	CIVIOS 1/2,5	Page 13, 96	1650,-
OBF 131C825				CMOS 1/2,5"	KERN Optics catalogue Page 18, 92	1355,-
OBF 132C825	ODC 825	5 MP	USB 2.0 (6,8 - 55 FPS)			1465,-
OBF 133C825					1 480 10, 72	1445,-
OBF 131C832					KEBN Oution antalanus	1500,-
OBF 132C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 18, 92	1610,-
OBF 133C832					1 480 10, 72	1590,-
OBF 131T241					KEBN Outing antalanus	1860,-
OBF 132T241	ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 18, 96	1970,-
OBF 133T241						1950,-



### The digital laboratory assistant with infinity optical unit and fixed, Köhler illumination

#### **Features**

- Laboratory microscopes from the OBL and OBN range are now also available to you as a comprehensive digital solution for your live investigations. Optionally available with an mounted tablet or C-mount camera. Naturally, the appropriate C-mount adapter is included with the delivery.
- The mounted KERN ODC 241 tablet-camera has been specially developed for simple, convenient and direct investigation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- The mounted C-mount camera is available in different versions and can be used anywhere
- For detailed information on the individual components, see the relevant product description of the individual item

 A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery

#### Technical data

- Infinity optical system
- Siedentopf 30° inclined/360° rotatable
- Eyepiece: HWF 10×/Ø 20 mm

#### OBL-1

- Quadplex nosepiece
- · Diopter adjustment: One-sided
- Overall dimensions W×D×H 395×200×380 mm
- Net weight approx. 7,7 kg
- Objective quality: Infinity E-Plan
- Objectives: 4×/10×/40×/100×
- Illumination OBL 135:
   20 W Halogen (transmitted)
- Illumination OBL 137: 3 W LED (transmitted)

#### OBN-1

- Quintuple nosepiece
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 390×200×400 mm
- Net weight approx. 10 kg
- Objective quality: Infinity Plan
- Objectives: 4×/10×/20×/40×/100×
- Illumination OBN 132:
   20 W Halogen (transmitted)
- Illumination OBN 135: 3 W LED (transmitted)

Model		Standard configuration (camera)					
	Included camera	Resolution camera	Interface	Sensor	Details microscope, camera	ex works €	
OBL 135C825 OBL 137C825	ODC 825	5 MP	USB 2.0 (6,8 - 55 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 20, 92	1515,- 1540,-	
OBL 135C832 OBL 137C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 20, 92	1660,- 1685,-	
OBL 135T241 OBL 137T241	ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 20, 96	2020,- 2045,-	
OBN 132C825 OBN 135C825	ODC 825	5 MP	USB 2.0 (6,8 - 55 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 26, 92	1865,- 1885,-	
OBN 132C832 OBN 135C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 26, 92	2010,- 2030,-	
OBN 132T241 OBN 135T241	ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 26, 96	2370,- 2390,-	









OZL 466 with camera



OZL 468 with camera



OZL 464 with tablet



OZL 466 with tablet



OZL 468 with tablet

### The flexible, affordable all-rounder with zoom function as a digital solution for schools, training companies, inspection authorities and laboratories

#### **Features**

- The flexible, cost-effective OZL-46 range is now also available to you as a comprehensive digital solution for your live investigations.
   Optionally available with an mounted tablet or C-mount camera. Naturally, the appropriate C-mount adapter is included with the delivery.
- The mounted KERN ODC 241 tablet-camera has been specially developed for simple, convenient and direct investigation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- The mounted C-mount camera is available in different versions and can be used anywhere

- For detailed information on the individual components, see the relevant product description of the individual item
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery

#### Technical data

- Optical system: Greenough optics
- Brightness adjustable (separate)
- Tube 45° inclined
- Magnification ratio: 6,4:1
- Light distribution 50:50
- Interpupillary distance 55 75 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H 300×240×420 mm
- Net weight approx. 5 kg
- Eyepiece: HWF 10×/Ø 20 mm
- Field of view: Ø 28,6 4,4 mm
- Objective: 0,7× 4,5×
- Stand OZL 464/466: Pillar style
- Stand OZL 468: Arm curved
- Illumination: 3 W LED (incident + transmitted)

Model		Standard configuration (camera)						
KERN	Included camera	Resolution camera	Interface	Sensor	Details microscope, camera	ex works €		
OZL 464C825						1000,-		
OZL 466C825	ODC 825	5 MP	USB 2.0 (6,8 - 55 FPS)	CMOS 1/2,5"	KERN Optics catalogue	1050,-		
OZL 468C825						1010,-		
OZL 464C832					Page 54, 92	1145,-		
OZL 466C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"		1195,-		
OZL 468C832						1155,-		
OZL 464T241					VEDN Onting antalanus	1505,-		
OZL 466T241	ODC 241	ODC 241 5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 54, 96	1555,-		
OZL 468T241					rage 54, 90	1515,-		









OZP-5 with camera



OZP-5 with tablet

### First-class optics as well as strong illumination combined with a high level of flexibility and digital tools

#### **Features**

- Stereomicroscopes from the OZM and OZP range are now also available to you as a comprehensive digital solution for your live investigations. Optionally available with an mounted tablet or C-mount camera. Naturally, the appropriate C-mount adapter is included with the delivery
- The mounted KERN ODC 241 tablet-camera has been specially developed for simple, convenient and direct investigation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- The mounted C-mount camera is available in different versions and can be used anywhere
- For detailed information on the individual components, see the relevant product description of the individual item
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of the delivery

#### **Technical data**

- Optical system: Greenough optics
- Brightness adjustable (separate)
- Light distribution: 50:50
- Diopter adjustment: Both-sided
- Net weight approx. 5,5 kg
- Eyepiece: HSWF 10×/Ø 23 mm
- · Stand: Pillar style
- Illumination: 3 W LED (incident + transmitted)

#### OZM-5


- Tube 45° inclined
- Magnification ratio: 6,4:1
- Interpupillary distance 52 76 mm
- Overall dimensions W×D×H 330×285×440 mm
- Field of view: Ø 32,8 5,1 mm
- Objective: 0,7× 4,5×

#### OZP-5

- Tube: 35° inclined
- Magnification ratio: 9,2:1
- Interpupillary distance 52 76 mm
- Overall dimensions W×D×H 330×285×470 mm
- Field of view: Ø 38,3 4,2 mm
- Objective: 0,6×-5,5×

Model			Standard configuration (camera)			<b>Price</b> excl. of VAT
	Included camera	Resolution camera	Interface	Sensor	Details microscope, camera	ex works €
OZM 544C825	ODC 825	5 MP	USB 2.0 (6,8 - 55 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 60, 92	1805,-
OZM 544C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 60, 92	1950,-
OZP 558C825	ODC 825	5 MP	USB 2.0 (6,8 – 55 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 62, 92	1985,-
OZP 558C832	ODC 832	5 MP	USB 3.0 (14,2 - 101,2 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 62, 92	2130,-
OZP 558T241	ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD (15 - 30 FPS)	CMOS 1/2,5"	KERN Optics catalogue Page 62, 96	2490,-

# 6 Stereo microscope sets

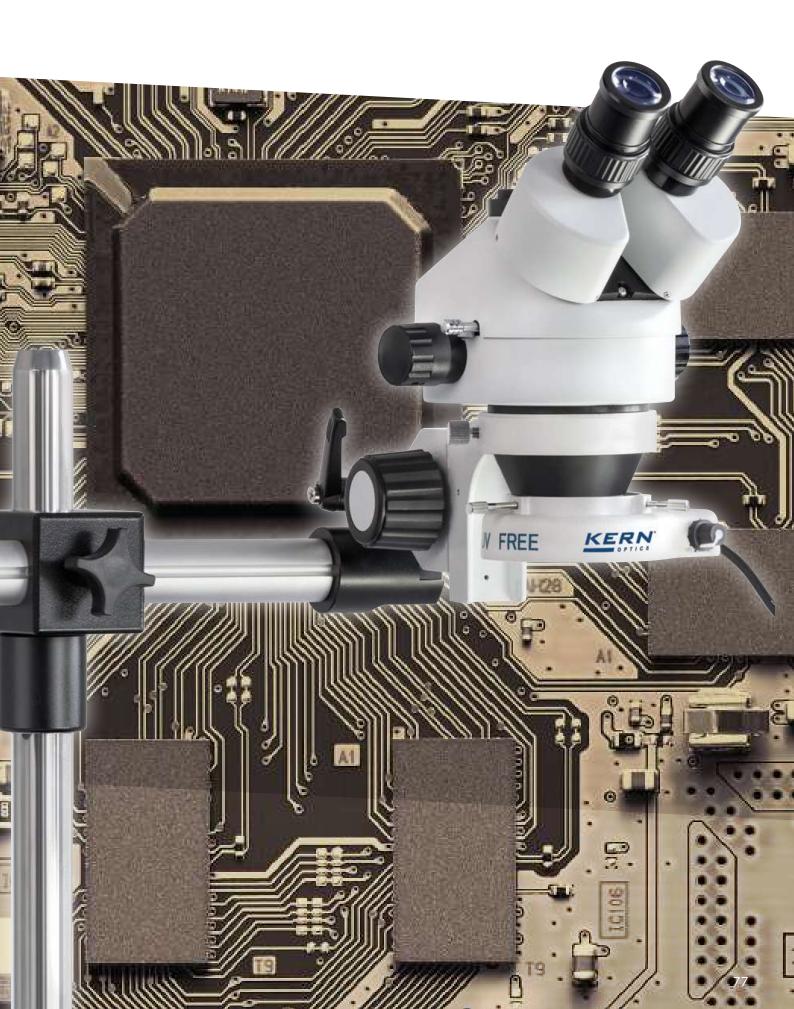



### Predefined stereo microscope sets with PREMIUM universal stand and illumination for your functional workplace

#### **Features**

- Sets which have already been defined, consisting of a stereo microscope head (p. 81), a universal stand (p. 86/87), a holder (p. 88), a ring illumination (p. 90) and a dust cover (p. 88) from our range
- Simple convenient affordable
- This saves you spending time on configuration and being spoilt for choice in the combination of different components. In this way you get an in expensive and highly-flexible solution for your microscope workplace

Model	Microscope h	nead	Stand	Holder	Illumination	<b>Price</b> excl. of VAT
KERN	Tube	Objective Zoom				ex works €
OZM 912	Binocular (OZM 546)	0,7×-4,5×	Telescopic arm with plate (OZB-A5201)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1440,-
OZM 913	Trinocular (OZM 547)	0,7×-4,5×	Telescopic arm with plate (OZB-A5201)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1740,-
OZM 932	Binocular (OZM 546)	0,7×-4,5×	ball-beared double telescopic arm with plate (OZB-A5203)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1730,-
OZM 933	Trinocular (OZM 547)	0,7×-4,5×	ball-beared double telescopic arm with plate (OZB-A5203)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	2030,-
OZM 952	Binocular (OZM 546)	0,7×-4,5×	Jointed arm with clamp (OZB-A5212)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1390,-
OZM 953	Trinocular (OZM 547)	0,7×-4,5×	Jointed arm with clamp (OZB-A5212)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1690,-
OZM 982	Binocular (OZM 546)	0,7×-4,5×	Spring loaded arm with clamp (OZB-A6302)	With coarse focusing knob Adjustable torque of the hand wheels	4,5 W LED ring light (OBB-A6102)	1560,-
OZM 983	Trinocular (OZM 547)	0,7×-4,5×	Spring loaded arm with clamp (OZB-A6302)	With coarse focusing knob Adjustable torque of the hand wheels	4,5 W LED ring light (OBB-A6102)	1860,-





OZL 961/963

### Predefined stereo microscope sets with ECO universal stand and illumination for your functional workplace

OZM 902/903

#### **Features**

- Sets which have already been defined (except OSE 409), consisting of a stereo microscope head (p. 81), a universal stand (p. 86/87), a holder (p. 88), a ring illumination (p. 90) and a dust cover (p. 88) from our range
- Simple convenient affordable
- This saves you spending time on configuration and being spoilt for choice in the combination of different components. In this way you get an in expensive and highly-flexible solution for your microscope workplace

Model	Microscope h	nead	Stand	Holder	Illumination	<b>Price</b> excl. of VAT
KERN	Tube	Objective Zoom				ex works €
OSE 409	Binocular (WF 10×/ Ø 20 mm)	1x (WD: 230 mm)	Swivel arm with block pedestal	With coarse focusing knob Adjustable torque of the hand wheels	3W LED goose neck (integrated)	355,-
OZL 961	Binocular (OZL 461)	0,7×-4,5×	Telescopic arm with plate	With coarse focusing knob Adjustable torque of the hand wheels	4,5 W LED ring light (OBB-A6102)	880,-
OZL 963	Trinocular (OZL 462)	0,7×-4,5×	Telescopic arm with plate	With coarse focusing knob Adjustable torque of the hand wheels	4,5 W LED ring light (OBB-A6102)	950,-
OZM 902	Binocular (OZM 546)	0,7×-4,5×	Telescopic arm with plate (OZB-A1201)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1290,-
OZM 903	Trinocular (OZM 547)	0,7×-4,5×	Telescopic arm with plate (OZB-A1201)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1570,-
OZM 922	Binocular (OZM 546)	0,7×-4,5×	ball-beared double telescopic arm with plate (OZB-A1203)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1340,-
OZM 923	Trinocular (OZM 547)	0,7×-4,5×	ball-beared double telescopic arm with plate (OZB-A1203)	With coarse focusing knob Adjustable torque of the hand wheels (OZB-A5301)	4,5 W LED ring light (OBB-A6102)	1620,-

# 6 Stereo microscope modular system

You can find sample diagrams of showing the configuration of a modular system like this on pages 83, 84 and 85 below.





Head of the microscope series OSF-5 (OSF 512, 514, 516)



Head of the microscope series OZL-46 (OZL 461, 462)



Head of the microscope series OZM-5 (OZM 546, 547)



Head of the microscope series OZP-5 (OZP 551, 552)



Head of the microscope series OZO-5 (OZO 556, 557)

### Individuality, variety and flexible working through our modular construction system ► Stereo microscope heads

#### **Features**

- To enable the highest level of flexibility for your special requirements and applications, we have a large selection of stereo microscope heads, universal stands and external illumination units, which are easy to combine
- Through the different properties of the stereo microscope heads, as well as the flexibility of the universal stands and the professional fixing of our brackets, we can configure your ideal microscope to suit your needs
- There are various microscope heads available from our product range for this purpose, both as binocular or trinocular versions
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the corresponding model outfit lists on the following pages

#### Technical data

- Optical system: Greenough optics
- Further technical data and model features is located in the tables below on the following pages
- OSF-5: p. 82
- OZL-46: p. 82
- OZM-5: p. 83
- OZP-5: p. 84
- OZO-5: p. 85

Model KERN	Tube	Tube angle	Eyepieces (included)	Interpupilary distance	Objective Zoom	Magnification ratio	Diopter adjustment	Price excl. of VAT ex works €
OSF 512	Binocular	45°	HSWF 10×/Ø 23 mm	52-76 mm	1×/2×	_	One-sided (-6/6)	320,-
OSF 514	Binocular	45°	HSWF 10×/Ø 23 mm	52-76 mm	1×/3×	-	One-sided (-6/6)	320,-
OSF 516	Binocular	45°	HSWF 10×/Ø 23 mm	52-76 mm	2×/4×	-	One-sided (-6/6)	320,-
OZL 461	Binocular	45°	HWF 10×/Ø 20 mm	55-75 mm	0,7×-4,5×	6,4:1	Both-sided (-5/5)	390,-
OZL 462	Trinocular	45°	HWF 10×/Ø 20 mm	52-76 mm	0,7×-4,5×	6,4:1	Both-sided (-5/5)	460,-
OZM 546	Binocular	45°	HSWF 10×/Ø 23 mm	52-76 mm	0,7×-4,5×	6,4:1	Both-sided (-6/6)	740,-
OZM 547	Trinocular	45°	HSWF 10×/Ø 23 mm	52-76 mm	0,7×-4,5×	6,4:1	Both-sided (-6/6)	1040,-
OZP 551	Binocular	35°	HSWF 10×/Ø 23 mm	52-76 mm	0,6×-5,5×	9,2:1	Both-sided (-6/6)	920,-
OZP 552	Trinocular	35°	HSWF 10×/Ø 23 mm	52-76 mm	0,6×-5,5×	9,2:1	Both-sided (-6/6)	1220,-
OZO 556	Binocular	35°	HSWF 10×/Ø 23 mm	52-76 mm	0,8×-7×	8,8:1	Both-sided (-6/6)	1350,-
OZO 557	Trinocular	35°	HSWF 10×/ø 23 mm	52-76 mm	0,8×-7×	8,8:1	Both-sided (-6/6)	1690,-

#### Fittings and accessories for the heads for the OSF-5 microscope range (OSF 512, OSF 514, OSF 516)

Eyepiece		Sį	pecifications - Objective	es .	
	Magnification	1×	2×	3×	4×
LICWE 10v	Total magnification	10×	20×	30×	40×
HSWF 10×	Field of view mm	ø 23	Ø 11,5	Ø 7,67	ø 5,75
	Total magnification	15×	30×	45×	60×
SWF 15×	Field of view mm	Ø 17	Ø 8,5	Ø 5,67	Ø 4,25
014/5 00	Total magnification	20×	40×	60×	80×
SWF 20×	Field of view mm	Ø 14	Ø 7	Ø 4,67	Ø 3,5
014/5-00	Total magnification	30×	60×	90×	120×
SWF 30×	Field of view mm	ø 9	Ø 4,5	Ø 3	Ø 2,25
Working distance	Working distance		105 mm	105 mm	105 mm

Model outfit			Model KERN		Order number	Price/piece excl. of VAT
		OSF 512	OSF 514	OSF 516		ex works €
	HSWF 10×/ø 23 mm	44	44	44	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	00	00	OZB-A5504	60,-
	SWF 20×/Ø 14 mm	00	00	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/Ø 9 mm	00	00	00	OZB-A5506	85,-
	HSWF 10×/Ø 23 mm (reticule 0,1 mm)	0	0	0	OZB-A5512	95,-
	SWF 15×/Ø 17 mm (reticule 0,05 mm)	0	0	0	OZB-A5513	110,-
	SWF 20×/Ø 14 mm (reticule 0,05 mm)	0	0	0	OZB-A5514	110,-

✓ = Included with delivery

O = Option

#### Fittings and accessories for the heads for the OZL-46 microscope range (OZL 461, OZL 462)

Eyepiece		Speci	ifications – Objective	es		
	Magnification	Standard		Auxiliar	y objectives	
		1,0×	0,5×	0,75×	1,5×	2,0×
HSWF 10×	Total magnification	7×-45×	3,5× - 22,5×	5,3×-33,8×	10,5× - 67,5×	14× - 90×
	Field of view mm	Ø 28,6 - 4,4	Ø 57,1-8,9	Ø 38,1 – 5,9	Ø 19-3	ø 14,3 - 2,2
HWF 15×	Total magnification	10,5× - 67,5×	5,3×-33,8×	7,9×-50,6×	15,5× - 101,3×	21× - 135×
HWF 13*	Field of view mm	Ø 21,4 - 3,3	Ø 42,9 – 6,7	Ø 28,5 - 4,4	Ø 14,3 - 2,2	ø 10,7 – 1,7
LICATE 20%	Total magnification	14×-90×	7×-45×	10,5× - 67,5×	21× - 135×	28× - 180×
HSWF 20×	Field of view mm	Ø 14,3 - 2,2	Ø 28,6 - 4,4	Ø 19,1-2,9	Ø 9,5 – 1,5	Ø 7,1 – 1,1
LIME 25~	Total magnification	17,5× - 122,5×	8,8×-56,3×	13,1×-91,9×	26,3× - 168,8×	35× - 225×
HWF 25×	Field of view mm	Ø 12,9 - 2,0	Ø 25,7 – 4,0	ø 17,2 - 2,7	Ø 8,6 – 1,3	Ø 6,4 – 1,0
Working distance		105 mm	177 mm	95 mm	47 mm	26 mm

Model outfit		Model KERN			Price/piece excl. of VAT
		OZL 461	OZL 462		ex works €
	HWF 10×/ø 20 mm	44	11	OZB-A4631	40,-
Eyepieces	HSWF 15×/Ø 15 mm	00	00	OZB-A4632	50,-
(30,0 mm)	HWF 20×/ø 10 mm	00	00	OZB-A4633	50,-
	HSWF 25×/ø 9 mm	00	00	OZB-A4634	55,-
	0,5×	0	0	OZB-A4641	80,-
Auvilianu ahiaatiwaa	0,75×	0	0	OZB-A4644	80,-
Auxiliary objectives	1,5×	0	0	OZB-A4642	80,-
	2,0×	0	0	OZB-A4643	80,-
	1× (focus adjustable)		✓	OZB-A4809	50,-
C-Mount	0,3× (focus adjustable)		0	OZB-A4810	85,-
	0,5× (focus adjustable)		0	OZB-A4811	95,-

✓ = Included with delivery

O = Option

Eyepiece		Spe	cifications - Obje	ctives					
	Magnification	Standard	Standard			Auxiliary objectives			
		1,0×	0,37×	0,5×	0,7×	1,5×	2×		
HSWF 10×	Total magnification	7× - 45×	2,59×-16,65×	3,5×-22,5×	4,9× - 31,5×	10,5× - 67,5×	14×-90×		
	Field of view mm	Ø 32,8 - 5,1	Ø 88,8 – 13,8	ø 65,7 – 10,2	ø 46,9 – 7,3	Ø 21,9 - 3,4	Ø 16,4 – 2,6		
	Total magnification	10,5× - 67,5×	3,89× - 25×	5,3× - 33,8×	7,4× - 47,2×	15,8×-101,3×	21× - 135×		
SWF 15×	Field of view mm	Ø 24,3 - 3,8	Ø 65,6 – 10,2	Ø 48,6 - 7,6	Ø 34,7 - 5,4	Ø 16,2 - 2,5	Ø 12,1 – 1,9		
CWE 20v	Total magnification	14×-90×	5,18×-33,3×	7×-45×	9,8×-63×	21× - 135×	28× - 180×		
SWF 20×	Field of view mm	Ø 20 – 3,1	Ø 54,1 – 8,4	Ø 40 - 6,2	Ø 28,6 - 4,4	Ø 13,3 - 2,1	Ø 10 – 1,6		
014/5 00	Total magnification	21× - 135×	7,77× - 50×	10,5× - 67,5×	14,7× - 94,5×	31,5×-202,5×	42× - 270×		
SWF 30×	Field of view mm	Ø 12,9 - 2	Ø 34,7 – 5,4	Ø 25,7 - 4	Ø 18,4 – 2,9	Ø 8,6 – 1,6	Ø 6,4 – 1		
Working distant	ce	110 mm	275 mm	195 mm	145 mm	50 mm	35 mm		

Model outfit		Mode	el KERN	Order number	Price/piece excl. of VAT
		OZM 546	OZM 547		ex works €
	HSWF 10×/ø 23 mm	44	11	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	00	OZB-A5504	60,-
	SWF 20×/ø 14 mm	00	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/ø 9 mm	00	00	OZB-A5506	85,-
(,,	HSWF 10×/ø 23 mm (reticule 0,1 mm)	0	0	OZB-A5512	95,-
	SWF 15×/Ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-
	SWF 20×/ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-
	0,37× only in combination with a universal stand	0	0	OZB-A5611	150,-
	0,5×	0	0	OZB-A5612	135,-
Achromatic	0,7×	0	0	OZB-A5613	135,-
auxiliary objectives	1,5×	0	0	OZB-A5615	135,-
	2,0×	0	0	OZB-A5616	140,-
	Soldering protection lens	0	0	OZB-A5614	45,-
	0,3× (focus adjustable)		0	OZB-A5701	120,-
	0,5× (focus adjustable)		0	OZB-A5702	120,-
	1,0× (focus adjustable)		0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703		0	OZB-A5704	255,-
	for SLR cameras (Nikon)		0	OZB-A5706	255,-
	for SLR cameras (Olympus)		0	OZB-A5707	255,-
	for SLR cameras (Canon)		0	OZB-A5708	255,-

✓ = Included with delivery O = Option

#### Functionality of our stereo microscope modular system

# **Step 1:**Select a microscope head (from page 81), a universal stand (page 86/87), a bracket (page 88) and a ring illumination unit (page 90), in order to generate a customised model.



#### Fittings and accessories for the heads for the OZP-5 microscope range (OZP 551, OZP 552)

Eyepiece	Specifications – Objectives							
	Magnification	Standard		Aı	uxiliary objecti	xiliary objectives		
		1,0×	0,37×	0,5×	0,7×	1,5×	2×	
HSWF 10×	Total magnification	6× - 55×	2,96× - 25,9×	3× - 27,5×	4,2× - 38,5×	9×-82,5×	12× - 110×	
HSWF 10^	Field of view mm	Ø 38,3 - 4,2	Ø 74,3 - 8,5	Ø 76,7 - 8,4	Ø 54,8 – 6	Ø 25,6 – 2,8	Ø 19,2 - 2,1	
OME 45	Total magnification	9×-82,5×	4,44×-38,9×	4,5× - 41,25×	6,3× - 57,75×	13,5×-123,75×	18× – 165×	
SWF 15×	Field of view mm	Ø 28,3 - 3,1	ø 57,4 - 6,6	Ø 56,7 – 6,2	Ø 40,5 – 4,4	Ø 18,9 – 2,1	Ø 14,2 – 1,5	
SWF 20×	Total magnification	12× - 110×	5,92×-51,8×	6×-55×	8,4× - 77×	18× – 165×	24× - 220×	
SWF 20*	Field of view mm	Ø 23,3 - 2,5	Ø 47,3 – 5,4	Ø 46,7 – 5,1	Ø 33,3 – 3,6	Ø 15,6 – 1,7	Ø 11,7 – 1,3	
SME 20%	Total magnification	18× - 165×	8,88× - 77,7×	9×-82,5×	12,6×-115,5×	27×-247,5×	36×-330×	
SWF 30×	Field of view mm	Ø 15 – 1,6	Ø 30,4 - 3,5	Ø 30 - 3,3	Ø 21,4 – 2,3	Ø 10 – 1,1	Ø 7,5 - 0,8	
Working distance		108 mm	275 mm	195 mm	145 mm	50 mm	35 mm	

Model outfit		Mode	el KERN	Order number	Price/piece excl. of VAT
		OZP 551	OZP 552		ex works €
	HSWF 10×/ø 23 mm	44	11	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	00	OZB-A5504	60,-
	SWF 20×/ø 14 mm	00	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/ø 9 mm	00	00	OZB-A5506	85,-
(22,2)	HSWF 10×/ø 23 mm (reticule 0,1 mm)	0	0	OZB-A5512	95,-
	SWF 15×/ø 17 mm (reticule 0,05 mm)	0	0	OZB-A5513	110,-
	SWF 20×/Ø 14 mm (reticule 0,05 mm)	0	0	OZB-A5514	110,-
	0,37× only in combination with a universal stand	0	0	OZB-A5611	150,-
	0,5×	0	0	OZB-A5612	135,-
Achromatic	0,7×	0	0	OZB-A5613	135,-
auxiliary objectives	1,5×	0	0	OZB-A5615	135,-
	2,0×	0	0	OZB-A5616	140,-
	Soldering protection lens	0	0	OZB-A5614	45,-
	0,3× (focus adjustable)		0	OZB-A5701	120,-
	0,5× (focus adjustable)		0	OZB-A5702	120,-
	1,0× (focus adjustable)		0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703		0	OZB-A5704	255,-
	for SLR cameras (Nikon)		0	OZB-A5706	255,-
	for SLR cameras (Olympus)		0	OZB-A5707	255,-
	for SLR cameras (Canon)		0	OZB-A5708	255,-

✓ = Included with delivery O = Option

#### Functionality of our stereo microscope modular system

#### Step 2:

Further illumination units (page 90) and a suitable protective dust cover (page 88) give you the opportunity to adapt the configuration, expansion and field of application of your ideal microscope individually to suit your own requirements

Fiber illumination Polarising ring light Dust cover







Fittings and accessories for the heads for the OZO-5 microscope range (OZO 556, OZO 557)

Eyepiece	Specifications - Objectives								
	Magnification	Standard		Auxiliary objectives					
		1,0×	0,37×	0,5×	0,7×	1,5×	2×		
HSWF 10×	Total magnification	8×-70×	2,96× - 25,9×	4×-35×	5,6×-49×	12× - 105×	16× - 140×		
HSWF 10^	Field of view mm	Ø 28,75 - 3,3	Ø 74,3 - 8,5	Ø 57,5 - 6,6	Ø 41,1 - 4,7	Ø 19,2 - 2,2	Ø 14,4 – 1,6		
CME 15	Total magnification	12× - 105×	4,44× - 38,9×	6×-52,5×	8,4×-73,5×	18× - 157,5×	24×-210×		
SWF 15×	Field of view mm	Ø 21,25 - 2,4	Ø 57,4 – 6,6	Ø 42,5 – 4,9	Ø 30,4 - 3,5	Ø 14,2 – 1,6	Ø 10,6 – 1,2		
SWF 20×	Total magnification	16× - 140×	5,92×-51,8×	8×-70×	11,2×-98×	24× - 210×	32× - 280×		
SWF 20^	Field of view mm	ø 17,5 - 2	Ø 47,3 - 5,4	Ø 35 - 4	Ø 25 - 2,9	ø 11,7 - 1,3	Ø 8,75 – 1		
SME 20	Total magnification	24× - 210×	8,88×-77,7×	12× - 105×	16,8× - 147×	36×-315×	48×-420×		
SWF 30×	Field of view mm	ø 11,25 - 1,3	Ø 30,4 - 3,5	Ø 22,5 - 2,6	Ø 16,1 – 1,8	Ø 7,5 - 0,9	Ø 5,625 - 0,6		
Working distance		108 mm	275 mm	195 mm	145 mm	50 mm	35 mm		

Model outfit	Model KERN				Order number	Price/piece excl. of VAT	
		OZO 551	OZO 552	OZO 553	OZO 554		ex works €
	HSWF 10×/ø 23 mm	44	11	44	44	OZB-A5503	60,-
	SWF 15×/ø 17 mm	00	00	00	00	OZB-A5504	60,-
	SWF 20×/ø 14 mm	00	00	00	00	OZB-A5505	70,-
Eyepieces (30,0 mm)	SWF 30×/ø 9 mm	00	00	00	00	OZB-A5506	85,-
(,,	HSWF 10×/ø 23 mm (reticule 0,1 mm)	0	0	0	0	OZB-A5512	95,-
	SWF 15×/ø 17 mm (reticule 0,05 mm)	0	0	0	0	OZB-A5513	110,-
	SWF 20×/ø 14 mm (reticule 0,05 mm)	0	0	0	0	OZB-A5514	110,-
	0,37× only in combination with a universal stand	0	0	0	0	OZB-A5611	150,-
	0,5×	0	0	0	0	OZB-A5612	135,-
Achromatic	0,7×	0	0	0	0	OZB-A5613	135,-
auxiliary objectives	1,5×	0	0	0	0	OZB-A5615	135,-
	2,0×	0	0	0	0	OZB-A5616	140,-
	Soldering protection lens	0	0	0	0	OZB-A5614	45,-
	0,3× (focus adjustable)			0	0	OZB-A5701	120,-
	0,5× (focus adjustable)			0	0	OZB-A5702	120,-
	1,0× (focus adjustable)			0	0	OZB-A5703	95,-
C-Mount	1,0× (with micrometer) only in combination with OZB-A5703			0	0	OZB-A5704	255,-
	for SLR cameras (Nikon)			0	0	OZB-A5706	255,-
	for SLR cameras (Olympus)			0	0	OZB-A5707	255,-
	for SLR cameras (Canon)			0	0	OZB-A5708	255,-

✓ = Included with delivery

#### Functionality of our stereo microscope modular system

#### Step 3:

When using a trinocular microscope configuration, select the microscope camera (from page 92) which meets your requirements. To find the appropriate C-mount adapter, which is essential to correctly connect the camera, please see the fitting lists of the selected microscope head (from page 82).







O = Option



### Individuality, variety and flexible working through our modular construction system ► PREMIUM universal stands

#### **Features**

- With our universal stands and basic stands, as well as microscope heads and external illumination units, you can configure your microscope to your own specifications and adapt it to your application
- Thanks to the versatile, adjustable universal stands it is possible to work in the very best way in all areas with the most varied of samples
- Large universal stands are available as stand base variants as well as with the option of a clamp for the edge or the centre of a bench. Depending on the model, you have the choice of a telescopic arm stand, a jointed arm stand or a telescopic double arm universal stand with ball bearings

#### **Technical data**

• Column height: 515 mm

#### OZB-A5201/OZB-A5211/OZB-A5221

• Length telescopic arm: 614 mm

#### OZB-A5202/OZB-A5212/OZB-A5222

Length jointed arm: 553 mm

#### OZB-A5203/OZB-A5213/OZB-A5223

· Length double arm: 545 mm

Model KERN	Description	Price excl. of VAT ex works €
OZB-A5201	Telescopic arm - Plate - excl. holder	460,-
OZB-A5211	Telescopic arm - Clamp Edge of bench (Range: max. 62 mm) - excl. holder	450,-
OZB-A5221	Telescopic arm - Clamp Centre of bench (hole required) - excl. holder	490,-
OZB-A5202 *	Jointed arm - Plate - excl. holder	490,-
OZB-A5212	Jointed arm - Clamp Edge of bench (Range: max. 62 mm) - excl. holder	480,-
OZB-A5222	Jointed arm - Clamp Centre of bench (hole required) - excl. holder	520,-
OZB-A5203	Telescopic double arm with ball bearings - Plate - excl. holder	740,-
OZB-A5213	Telescopic double arm with ball bearings – Clamp Edge of bench (Range: max. 62 mm) – excl. holder	730,-
OZB-A5223	Telescopic double arm with ball bearings - Clamp Centre of bench (hole required) - excl. holder	760,-















OZB-A6301

### Individuality, variety and flexible working through our modular construction system ► ECO universal stands

#### **Features**

- With our universal stands and basic stands, as well as microscope heads and external illumination units, you can configure your microscope to your own specifications and adapt it to your application
- Thanks to the versatile, adjustable universal stands it is possible to work in the very best way in all areas with the most varied of samples
- Small universal stands are available as stand base variants as well as with the option of a clamp for the edge of a bench. Depending on the model, you have the choice of a telescopic arm stand or a telescopic double arm universal stand with ball bearings
- The spring loaded universal stands including bench clamp will make your daily work with your stereo microscope easier. Now including coarse adjustment knob for easy, flexible focussing

#### **Technical data**

#### OZB-A1201/OZB-A1211

- Column height: 430 mm
- · Length telescopic arm: 385 mm

#### OZB-A1203/OZB-A1213

- · Column height: 430 mm
- · Length telescopic arm: 480 mm

#### OZB-A6302

- Height spring loaded arm: 525 mm
- Length spring loaded arm: 620 mm

#### OZB-A6301

· Column height: 300 mm

#### OZB-A6303

· Height spring loaded arm: 400 mm

· Length spring loaded arm: 850 mm

Model KERN	Description	Price excl. of VAT ex works €
OZB-A1201	Telescopic arm - Plate - excl. holder	330,-
OZB-A1211	Telescopic arm - Clamp Edge of bench (Range: max. 40 mm) - excl. holder	320,-
OZB-A1203	Jointed arm - Plate - excl. holder	380,-
OZB-A1213	Jointed arm - Clamp Edge of bench (Range: max. 40 mm) - excl. holder	360,-
OZB-A6302	Spring loaded arm (Pneumatic spring) - Clamp (Range: max. 50 mm) - with holder (Coarse focusing knob)	680,-
OZB-A6303 🔤	Spring loaded arm (Compression spring) - Clamp (Range: max. 50 mm) - with holder (Coarse focusing knob)	420,-
OZB-A6301	Pillar style stand with "C"-shape base – excl. holder	265,-



OZB-A5301



OZB-A5306

### Individuality, variety and flexible working through our modular construction system ► Holders

#### **Features**

- There are two microscope head holders available for these flexible, modular systems.
   These brackets are suitable for all stereo microscope stands and universal stands (excluding spring loaded arm), to make focusing possible
- The first variant available is a holder with adjustable handwheel as well as adjustment of the torque for your configuration
- For professional applications you have the choice of a mount with coarse and fine focusing knob for the very best focusing operation
- Diameter of the connector for the microscope head: 76 mm
- Diameter of the connector for the stand:
   25 mm

Model KERN	Description	Price excl. of VAT ex works €
OZB-A5301	Holder with adjustable torque of the hand wheel. Suitable for all universal stands (except of spring loaded arm) and for all basic stands as possible accessories.	145,-
OZB-A5306	Holder with coaxial coarse and fine focusing knob and adjustable torque of the hand wheel. Suitable for all universal stands (except of spring loaded arm) and for all basic stands as possible accessories.	245,-

#### Stereo microscope modular system - Dust covers KERN OBB-C



### Individuality, variety and flexible working through our modular construction system ► Dust covers

#### Features

- When working with microscopes, we offer dust covers to give greater ease of use.
   By using these, you can easily avoid the time-consuming cleaning work which is necessary with routine use of your microscope
- Depending on the size of your microscope set or your microscope configuration you can select between three different models
- Please find detailed information in the following model outfit list

Model KERN	Description	Suitable for	Price excl. of VAT ex works €
OBB-A1387	Size 1: 485×440 mm	Stereo microscope heads	30,-
OBB-A1388	Size 2: 600×600 mm	Stereo microscope heads in combination with basic stands	30,-
OBB-A1389	Size 3: 650×750 mm	Stereo microscope sets, stereo microscope heads in combination with universal stands	30,-

### 7 External light sources for stereomicroscopes

Ring illumination and cold light sources



### Professional illumination units guarantee outstanding, uniform and strong illumination

These illumination units are also available with UK mains plug. For more information on this, visit our online shop or give us a ring







OBB-A6102



OZB-A4571 OZB-A4572

#### OZB-A7101

#### **Features**

- Choose your favourite external illumination here to achieve maximum flexibility and greatest possible ease of use in stereo microscopy
- These professional illumination units provide a quality of light at a high, constant intensity at all times
- Regardless of whether your choice is space-saving ring lights or cold light sources using optical fiber, our range is all you can wish for
- With the OZB-A7101 polarisation ring illumination unit, you also have an excellent component which has been specially optimised for observing shiny surfaces
- Naturally, these external illumination units also fit your standard stereo microscope
- Exception: The ring illumination units cannot be used in combination with the following ranges: OSE-1, OSF-4G, OZL-45R, OZC-5 and OZG-4

Model KERN	Illuminance	Inner Ø mm	Colour temperature	Brightness adjustable	Illumination by segments	Polarising filter	Price excl. of VAT ex works €
OZB-A4571	3W-LED	60	7000 - 11000	✓			90,-
OZB-A4572	3W-LED	60	6500 - 7000	✓	✓		140,-
OBB-A6102	4,5W-LED	63	approx. 7600	✓			135,-
OZB-A7101	4,5W-LED	62	6500 - 7000	✓		✓	590,-

✓ = Included with delivery

O = Option

#### Fiber illumination KERN OZB-IF









Application example

#### **Features**

 With the OZB-A4516 20 W-LED goose neck illumination unit with focusable light beam, you can adjust the illumination to suit your needs. Spot or scattered radiation means that you can achieve the very best illumination of your sample.

Model KERN	Description	Length	Illuminance	Colour temperature	Brightness adjustable	Price excl. of VAT ex works €
OZB-A4515	Dual fiber LED unit	300	6W	5600 - 6300	✓	150,-
OZB-A4516	Dual fibre LED cold light source	540	20W	6400	✓	620,-

✓ = Included with delivery

O = Option

# 8 Microscope cameras & Software



### Specialists in microscopy for measurement, counting, documentation, archiving and image processing

#### **Features**

- A large selection of microscope cameras is available for your individual applications
- The universal microscope cameras can be used anywhere and can be connected to the microscope as well as to a laptop or PC using the USB cable (USB 2.0 or USB 3.0, see table)
- The power supply is through the USB cable, which means that no additional power supply is required
- Your daily work is made significantly easier with the very best synchronisation, a high frame rate as well as stable image performance together with our camera software microscope VIS KERN OXM 901 which we deliver with the product
- For details about our software please refer to the "Camera software microscope VIS KERN OXM 901" product group in the catalogue (page 97) or see www.kern-sohn.com
- These universal cameras can also be connected to all microscopes available on the market offering the appropriate C-mount adapter for the particular microscope

#### **Accessories**

 Object micrometer, for calibrating the software measuring function, division 0.01 mm, KERN ODC-A2403, € 20,-

#### C-mount cameras - USB 2.0/3.0 KERN ODC-82 · ODC-83



#### **Features**

- Through the proven CMOS technology, in connection with the USB 2.0 or USB 3.0 the images are shown quickly and clearly
- These cameras are also ideal for more demanding applications, such as, for example, darkfield, phase contrast and for fluorescence applications
- As well as the camera, the delivery includes our multi-lingual camera software Microscope VIS Basic KERN OXM 901, a USB cable (length: 2 m) various eyepiece adapters and an object micrometer to calibrate the software
- Please order the appropriate C-mount adapter to fit your KERN microscope now

STANDAR	D		
•	$\stackrel{\bullet}{\longleftrightarrow}$		
USB 2.0	USB 3.0	SOFTWARE	1 DAY

Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Colour/ Monochrome	Supported operating system	Price excl. of VAT ex works €
ODC 824	3,1 MP	USB 2.0	11,5 - 45	CMOS	1/2"	colour	Win XP, Vista, 7, 8, 10	370,-
ODC 825	5,1 MP	USB 2.0	6,8 – 55	CMOS	1/2,5"	colour	Win XP, Vista, 7, 8, 10	385,-
ODC 831	3,1 MP	USB 3.0	27,3 - 53,3	CMOS	1/3"	colour	Win XP, Vista, 7, 8, 10	500,-
ODC 832	5,1 MP	USB 3.0	14,2 - 101,2	CMOS	1/2,5"	colour	Win XP, Vista, 7, 8, 10	530,-

#### C-mount camera - High resolution KERN ODC-84





#### **Features**

- The high-resolution, professional ODC-84 range offers you an impressive 20 megapixel resolution which will give you bright detailed views of your sample. By using the integrated USB 3.0 interface, live images are transferred to the KERN OXM 902 for processing and documentation
- Power supply is through the USB interface so that there is no requirement for an external power source.
- As well as the camera, the delivery includes our multi-lingual camera software Microscope VIS Pro KERN OXM 902, a USB cable (length: 2 m), various eyepiece adapters and an object micrometre to calibrate the software
- Please order the appropriate C-mount adapter (only 1,0× possible) to fit your KERN microscope now



Can only be used in combination with compound microscopes

# USB 2.0 SOFTWARE 1 DAY

Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Colour/ Monochrome	Supported operating system	Price excl. of VAT ex works €
ODC 841	20 MP	USB 3.0	15 – 60	CMOS	1"	colour	Win XP, Vista, 7, 8, 10	1090,-



#### **Features**

- The ODC 851 HDMI microscope camera has been specially developed for direct HDMI connection to your HDMI compatible display device. The images can be stored straight onto the SD card which is delivered with the product or can be transferred to your PC or laptop for further processing using the USB 2.0 cable in combination with the OXM 901 software.
- The HDMI autofocus camera ODC 852 offers you a perfect, effective solution for modern microscopy. The autofocus function automatically detects and adjusts the focus level so that you always have a razor-sharp image. Ideal for all applications in connection with a KERN stereo microscope.
- · Realtime images can be transferred to an HDMI-compatible display device using the HDMI connection and they can also be stored on the SD card which was delivered

- with the product. As an alternative, data can also be transferred using the WLAN module (ODC 852) or USB 2.0 cable to a PC or laptop in combination with the KERN OXM 902 software which is included with the delivery
- Power supply is from an external 12 V power unit
- Scope of delivery ODC 851: Camera, USB mouse, USB 2.0 cable (length: 2 m), HDMI cable (length: 2 m), SD card (16 GB) and camera software Microscope VIS Basic KERN OXM 901
- · Scope of delivery ODC 852: Camera, USB mouse, USB 2.0 cable (length: 2 m), HDMI cable (length: 2 m), SD card (16 GB), WiFi adapter and camera software Microscope VIS Pro KERN OXM 902
- · Please order the appropriate C-mount adapter to fit your KERN microscope now















6	
/	1 DAY

Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Colour/ Monochrome	Supported operating system	Price excl. of VAT ex works €
ODC 851	2 MP	HDMI, USB 2.0, SD	30 - 60	CMOS	1/2,8"	colour	Win XP, Vista, 7, 8, 10	990,-
ODC 852* NEW	5 MP	HDMI, USB 2.0, SD, WLAN	25 - 60	CMOS	1/1,8"	colour	Win XP, Vista, 7, 8, 10	1060,-

<sup>\*</sup>can only be used in combination with stereo microscopes

New model

#### C-mount camera - High resolution KERN ODC-84







#### The cooled camera for professional fluorescence examinations

#### **Features**

- The ODC 861 camera with Peltier cooling technology has been specially designed for fluorescent applications. It is able to significantly compensate for image noise associated with weak lighting. Due to its high resolution and light-sensitive Sony CMOS colour sensor it proves first-class images. The practical, sturdy storage box serves as protection and for transportation of this premium camera
- · Realtime images can be transferred straight to a PC or laptop using the integrated USB 3.0 interface. As an alternative, 2 USB 2.0 interfaces are available, to operate the camera with the KERN OXM 902 software which is included with the delivery
- · Power supply is from an external 12 V power unit
- · Please order the appropriate C-mount adapter (only 1,0× possible) to fit your KERN microscope now



Can only be used in combination with compound microscopes







Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Colour/ Monochrome	Supported operating system	Price excl. of VAT ex works €
ODC 861	20 MP	USB 3.0	5 – 30	CMOS	1"	colour	Win XP, Vista, 7, 8, 10	1990,-



ODC-87, ODC-88



Eyepiece camera fixed into the tube

#### **Features**

- With the KERN eyepiece cameras you can convert your standard microscope to a digital microscope, by replacing one eyepiece of your non-digital microscope with an eyepiece camera and connect this to your computer via USB
- · The universal eyepiece can be connected to the microscope as well as to a laptop or PC using the USB cable (2.0 or 3.0, see table)
- The power supply is through the USB cable, which means that no additional power supply is required
- · Your daily work is made significantly easier with the very best synchronisation, a high frame rate as well as stable image performance together with our software
- As well as the camera, the delivery includes a simplified version of our multi-lingual camera software Microscope VIS KERN OXM 901, a USB cable (length: 1,5 m) and an object micrometer to calibrate the software











Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Colour/ Monochrome	Supported operating system	Price excl. of VAT ex works €
ODC 872	1,3 MP	USB 2.0	7,5 - 12,5	CMOS	1/3"	colour	Win XP, Vista, 7, 8, 10	185,-
ODC 874	3 MP	USB 2.0	3 - 7,5	CMOS	1/2,7"	colour	Win XP, Vista, 7, 8, 10	220,-
ODC 881	5 MP	USB 3.0	15 – 30	CMOS	1/2,5"	colour	Win XP, Vista, 7, 8, 10	295,-

USB microscope - USB 2.0 KERN ODC-89

### The digital USB microscope for rapid testing or for hobby use



#### **Features**

- The USB hand-held microscope is designed for rapid and simple observations. Ideally suited for coins, plants, insects and skin samples for all hobby scientists, children and students
- With the USB microscope you can easily adjust the magnification to suit all conventional samples. The zoom range can be adjusted to a magnification of 10× as well as 200×
- The eight LEDs fitted in the ring shape ensure strong and effective illumination of your sample. Use the adjustment wheel on the cable to control the illumination setting
- · As well as the camera, you will also find a simplified version of our multi-lingual camera software Microscope VIS KERN OXM 901 included with delivery
- · There are two stands available for you to use as a column
- · Cable length: 1,4 m

#### Stand with focus wheel:

· Work area: 150×80mm • Focus range: 60 mm

Overall dimensions: 150×80×135 mm

**ODC 895** 










Model KERN	Resolution	Interface	FPS	Sensor	Sensor size	Supported operating system	Magnifica- tion levels	Focusing stand	Illumination	Price excl. of VAT ex works €
ODC 895	2 MP	USB 2.0	15 – 30	CMOS	1/3,2"	Win XP, Vista, 7, 8, 10	10×, 200×	Focus wheel	8× LED	230,-







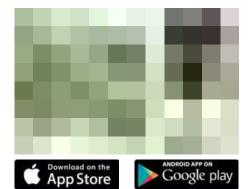

### Innovative hand-held microscope for mobile applications with immediate display of the image on a smartphone or tablet

#### **Features**

- The digital WLAN hand-held microscope is designed for rapid and simple surface observations. Ideally suited for coins, bank notes, stamps, circuit boards, plants, insects, gems and skin samples for industrial use, for all hobby scientists, children and students
- The KERN ODC 910 WLAN microscope has been specially developed for direct connection to your WLAN-enabled smartphone or tablet with iOS or android
- · During the live transfer to your smartphone or tablet you can take photos and videos of the sample you are investigating, and these can also be stored on your device. For larger videos you can also insert a mini SD card directly into the microscope
- · With the WLAN microscope you can easily adjust the magnification to suit all conventional samples. The focus can be adjusted to a magnification of 10× as well as 200×

- The six LEDs fitted in a ring shape ensure strong and effective illumination of your sample. Use the adjustment wheel on the microscope to control the illumination setting.
- You can download the app for the ODC 910 WLAN microscope from the Apple App Store or the Android Google Play Store free of charge and this app enables you to directly transfer images and videos from the microscope to your smartphone or tablet through a simple connection
- · The scope of delivery includes the WLAN microscope with integrated rechargeable battery pack, a flexible column which is easy to adjust and which has a swan neck so that you can achieve the ideal height setting, as well as a mains adapter

#### STANDARD










	odel ERN	Resolution	Interface	FPS	Sensor	Sensor size	Supported operating system	Magnifica- tion levels	Focusing stand	Illumination	Price excl. of VAT ex works €
0	DC 910	2 MP	WLAN, SD	15 – 30	CMOS	1/4"	Android, iOS	10×, 200×	Goose neck	6× LED	170,-





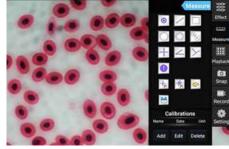

Integrated software with measuring function

Digital microscopy brought up to date - tablet with integrated camera for optimal observation and digital documentation of the sample

#### **Features**

- · A 2-in-1 solution in digital microscopy as a universal system for trinocular microscopes with C-mount adapter. The ODC 241 microscope-tablet-camera consists of a large Android tablet in combination with a 5-MP camera
- The KERN ODC 241 tablet-camera has been specially developed for simple and direct observation of the sample on the screen. Ideally suited for school pupils and students in education or for demonstration purposes in the laboratory
- · As well a live transfer of the image to the Android table, the integrated 5-MP camera also means that images and videos can be created for the documentation.
- · Simple measuring functions such as, for example, functions for measuring distance, surfaces and angles as well as a manual counting function are also available


- · Automatic white balance and automatic contrast adjustment can be performed quickly and easily, which enables efficient working procedures
- · A range of additional functions are provided through the integrated interfaces, such as, for example
  - Data storage on a USB stick or SD card
  - Connection to a USB mouse
  - Transfer of the live image to an external screen using HDMI
- Transfer of stored data to external receivers using WLAN
- · The delivery includes the tablet camera with pre-installed software as well as the mains adapter

#### **Technical data**


- 9.7" LCD-Touchscreen
- · Screen resolution: 2048×1536 pixels
- · CPU: Quad Core Cortex-A17; 1,8 GHz
- · Overall dimensions W×D×H 238×51×206 mm
- Net weight approx. 0,65 kg



Model KERN	Resolution Camera	Interface	FPS	Sensor	Sensor size	Supported operating system	Price excl. of VAT ex works €
ODC 241	5 MP	WLAN, USB 2.0, HDMI, SD	15 – 30	CMOS	1/2,5"	Android 5.1	890,-

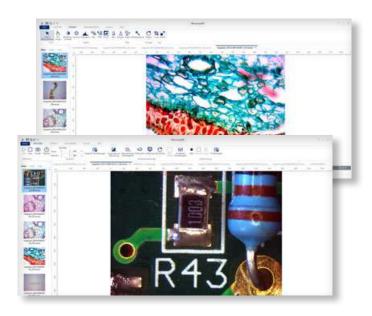



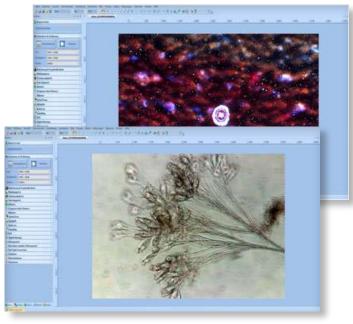

The digital specialist for measurement, counting and archiving – free of charge with all KERN microsope cameras

#### Features

#### OXM 901\*

- The camera software KERN microscope
   VIS Basic is a multi-lingual software, which we have developed specially for all available Kern microscope cameras
- As well as the streaming function for the object to be viewed, the software also offers you an image snapshot function, as well as a video function
- Various measuring functions such as, for example, functions for measuring distance, surfaces and angles and a manual counting function are available. In addition there are extensive image processing and documentation functions available, and of course an exporting function to Microsoft Office applications Word® and Excel®
- With the display settings you can display different measurements, grid sizes, scales and rulers for the very best measuring procedure
- Automatic white balance and automatic contrast adjustment can be performed quickly and easily, which enables efficient working procedures

#### OXM 902

- With the camera software KERN microscope VIS Pro in essence, all functions of the Basic variant are supported, however, in addition many other features are integrated which can be used for image analysis at a more professional level.
- The following highlights are included:
- Image Stitching
- Image Stacking
- Expanded measuring functions
- Auto counting function
- With this software it is possible operate all available KERN microscope cameras

#### **Technical data**

- Can be used for Microsoft Windows XP, Windows Vista, Windows 7, Windows 8, 8.1 and Windows 10
- Depending on the language setting of your Windows operating system the software KERN microscope VIS Basic & Pro will be identified and installed in the current language. This can be changed manually at any time
- The software is available in the following languages:
- **OXM 901**: German, English, Spanish, Italian, French, Portuguese, Polish
- OXM 902: German, English, Spanish, Italian, French, Portuguese, Polish, Russian, Turkish, Chinese, Japanese, Korean
- As well as the software CD, a USB cable and an object micrometer are included with all KERN cameras as well as all digital microscopes Please refer to the documentation for the software in the download area of our website www.kern-sohn.com for further details

<sup>\*</sup>Cannot be used in combination with the following cameras: ODC 841, ODC 852, ODC 861



### Refractometers

9 Analogue refractometers – type: hand-held 100

10 Digital refractometers – type: hand-held

106

11 Abbe refractometers – type: desktop

112





**Daniel Junger** Manager KERN Optics

Tel. +49 7433 9933-155 Fax +49 7433 9933-29195 daniel.junger@kern-sohn.com







Now also available with calibration certificate, see page 114!

#### Refractive index measurement for laboratories and the industry

#### **Features**

- The KERN ORA refractometers are universal, maintenance free analogue handheld refractometers
- The handy and robust design allows the easy, efficient and sustainable use in everyday life
- Manually calculated conversions and errors of the user are avoided by multiple selectable scales
- These scales are especially developed, exactly calculated and checked. They are also characterized by their thin and clear lines
- The optical system and the prism cover are made of special material which allows a low-tolerance measuring
- All ORA models are equipped with an eyepiece for easy and smooth setting for many different diopter strengths

- The models marked with "ATC" have an automatic temperature compensation which enables accurate measurement at different ambient temperatures (10 °C/30 °C)
- $\bullet$  The follwoing accessory-parts are included:
  - Storage box
  - Calibration liquid
  - Calibration block (if required)
  - Pipette
- Small screwdriver
- Cleaning tissue
- Further accessories are available optionally

#### **Technical data**

- Die-cast housing of copper-aluminium alloy, chrome coated
- Measurement temperature without ATC: 20 °C
- Measurement temperature range with ATC: 10 °C/30 °C
- Dimensions of the box: 205×75×55 mm (depending on the model)
- Product length: approx. 130 200 mm (depending on the model)
- Net weight approx. 135 600 g (depending on the model)





#### Scope of application: Sugar

The following models are particularly suitable for the measurement of the "BRIX" value. They are used to determine the sugar content in food, especially in fruit, vegetables, juice and soft drinks. In the same ideal way these refractometers serve for monitoring processes in the industry (coolant monitoring, oils, lubricants and fats).

The main scope of applications is:

- Industry: Monitoring of lubricants for process and quality control
- Food industry: Beverages, fruits and sweets
- Agriculture: Determination of the degree of ripeness of fruits for quality control in harvesting
- · Restaurants and large-scale catering establishment

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 10BB	Brix	0 - 10 %	0,1 %		80,-
ORA 10BA	Brix	0 - 10 %	0,1 %	✓	85,-
ORA 18BB	Brix	0 – 18 %	0,1 %		80,-
ORA 20BB	Brix	0 – 20 %	0,1 %		80,-
ORA 20BA	Brix	0 – 20 %	0,1 %	✓	85,-
ORA 32BB	Brix	0 - 32 %	0,2 %		80,-
ORA 32BA	Brix	0 - 32 %	0,2 %	✓	85,-
ORA 62BB	Brix	28 - 62 %	0,2 %		80,-
ORA 62BA	Brix	28 - 62 %	0,2 %	✓	85,-
ORA 82BB	Brix	45 - 82 %	0,5 %		80,-
ORA 80BB	Brix	0 - 80 %	0,5 %		80,-



#### Scope of application: Honey

The following models are particularly suitable for the measurement of the "BRIX" value, as well as the water content in honey and "degrees Baumé" to determine the relative density of liquids.

- Beekeeping
- Honey production

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 3HB	Brix Baumé Water content	58 – 92 % 38 – 43 °Bé 12 – 27 %	0,5 % 0,5 °Bé 1 %		85,-
ORA 3HA	Brix Baumé Water content	58 – 92 % 38 – 43 °Bé 12 – 27 %	0,5 % 0,5 °Bé 1 %	~	90,-
ORA 6HB	Water content	12 - 30 %	0,1 %		95,-
ORA 6HA	Water content	12 - 30 %	0,1 %	✓	100,-



#### Scope of application: Salt

The following models are particularly suitable for the measurement and concentration control of the mass fraction of natrium chloride in water as well as of the content of NaCl (salt) in water. This is often used in the preparation and the cooking of sauces, bases for pastries, the production of brines (e.g. for white cheese) and the preparation of seafood and marinades for meat.

The main scope of applications is:

- Food industry
- · Restaurants and large-scale catering establishment
- · Aquaristic: Fishkeepers/Fishfarmers in sea and sweetwater

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 1SB	Salinity specific gravity	0 – 100 ‰ 1,000 – 1,070 sg	1 ‰ 0,001 sg		80,-
ORA 1SA	Salinity specific gravity	0 – 100 ‰ 1,000 – 1,070 sg	1 ‰ 0,001 sg	✓	85,-
ORA 2SB	Salt (NaCl)	0 - 28 %	0,2 %		80,-
ORA 2SA	Salt (NaCl)	0 – 28 %	0,2 %	✓	85,-
ORA 3SB	Salt (NaCl) Brix	0 - 28 % 0 - 32 %	0,2 % 0,2 %		80,-
ORA 3SA	Salt (NaCl) Brix	0 - 28 % 0 - 32 %	0,2 % 0,2 %	✓	85,-



#### Scope of application: Wine

The following models are particularly suitable for the measurement of the content of sugar in fruits. It indicates the expected °Alcohol of the fruit. The degree of ripeness of fruit (fruit-sugar) can also be determined, such as e.g. grapes.

- · Agriculture: Wine-growing and fruit-growing
- · Wine-production
- Must and alcohol production

°Oe = Degree Oechsle, °KMW = Klosterneuburger Must balance

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 1WB	Oechsle KMW (Babo) Brix	0 – 140 °Oe 0 – 25 °KMW 0 – 32 %	1 °Oe 0,25 °KMW 0,2 %		80,-
ORA 1WA	Oechsle KMW (Babo) Brix	0 – 140 °Oe 0 – 25 °KMW 0 – 32 %	1 °Oe 0,25 °KMW 0,2 %	~	85,-
ORA 3WB	Oechsle Brix	30 - 140 °Oe 0 - 32 %	1 °Oe 0,2 %		80,-
ORA 3WA	Oechsle Brix	30 - 140 °Oe 0 - 32 %	1 °Oe 0,2 %	✓	85,-
ORA 7WB	Oechsle KMW (Babo) Brix	30 - 140 °Oe 0 - 25 °KMW 0 - 32 %	1 °Oe 0,2 °KMW 0,2 %		80,-
ORA 7WA	Oechsle KMW (Babo) Brix	30 - 140 °Oe 0 - 25 °KMW 0 - 32 %	1 °Oe 0,2 °KMW 0,2 %	<b>✓</b>	85,-
ORA 1AB	Percentage by volume Percentage by volume	0 - 50 % (v/v) 50 - 80 % (v/v)	1 % (v/v) 2,5 % (v/v)		80,-
ORA 2AB	Percentage by mass Percentage by mass	0 - 50 % (w/w) 50 - 80 % (w/w)	1 % (w/w) 2,5 % (w/w)		80,-



#### Scope of application: Urine

The following models are particularly suitable for the measurement of the specific gravity (sg) in urine, the quantitiy of serum (serumproteine) in urine (doping control among athletes), and the refractive index.

The main scope of applications is:

- Hospitals
- Doctor's surgeries/Physicians
- Medical training institutions
- Nursing homes
- · Sports medicine (doping test)
- Veterinary

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 2PB	Serum protein Urine (spec. gravity) Refractive index	0 – 12 g/dl 1,000 – 1,050 sgU 1,3330 – 1,3600 nD	0,2 g/dl 0,002 sgU 0,0005 nD		80,-
ORA 2PA	Serum protein Urine (spec. gravity) Refractive index	0 – 12 g/dl 1,000 – 1,050 sgU 1,3330 – 1,3600 nD	0,2 g/dl 0,002 sgU 0,0005 nD	~	85,-
ORA 5PB	Serum protein Urine (s. g. dog) Urine (s. g. cat)	2 – 14 g/dl 1,000 – 1,060 sgU 1,000 – 1,060 sgU	0,1 g/dl 0,001 sgU 0,001 sgU		80,-



#### Scope of application: Industry/Automotive

The following models are particularly suitable for the measurement and determination of AdBlue, glycol concentration (ethylene (EG) and propylene (PG)), battery fluid (BF), urea, the freezing point of fountain solution (CW). Furthermore these models are suitable for the measurement of thermal exchange systems.

- · Automotive industry: Car-workshops and producers
- · Chemical industry
- · Solar industry: Antifreeze monitoring
- · Geothermal industry: Brine-concentration-measurement for ground heat
- Forestry/Lumbermen

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 4FB	EG (G11/12) PG (G13) CW BF	-50 - 0 °C -50 - 0 °C -40 - 0 °C 1,10 - 1,40 kg/l	1 °C 1 °C 5 °C 0,01 kg/l		80,-
ORA 4FA	EG (G11/12) PG (G13) CW BF	-50 - 0 °C -50 - 0 °C -40 - 0 °C 1,10 - 1,40 kg/l	1 °C 1 °C 5 °C 0,01 kg/l	~	85,-
ORA 1UB	Urea	0 – 40 %	0,2 %		80,-
ORA 1UA	Urea	0 - 40 %	0,2 %	✓	85,-
ORA 4UB	Urea EG (G11/12) PG (G13) CW BF	30 - 35 % -50 - 0 ° C -50 - 0 ° C -40 - 0 ° C 1,10 - 1,40 kg/l	0,2 % 1 °C 1 °C 5 °C 0,01 kg/l		80,-
ORA 4UA	Urea EG (G11/12) PG (G13) CW BF	30 - 35 % -50 - 0 ° C -50 - 0 ° C -40 - 0 ° C 1,10 - 1,40 kg/l	0,2 % 1 °C 1 °C 5 °C 0,01 kg/l	<b>*</b>	85,-




#### Scope of application: Expert applications

The following models have a special large measuring range for the refractive index and large divided scales for the measurement and clear reading of Brix values.

The main scope of applications is:

• Universal application, especially when extra large measuring ranges are required

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 80BE	Brix	0 - 50 % 50 - 80 %	0,5 % 0,5 %		140,-
ORA 90BE	Brix	0 - 42 % 42 - 71 % 71 - 90 %	0,2 % 0,2 % 0,2 %		295,- 🕛
ORA 1RE	Refractive index	1,333 – 1,405 nD 1,405 – 1,468 nD 1,468 – 1,517 nD	0,005 nD 0,005 nD 0,005 nD		295,- 🕛
ORA 4RR	Refractive index	1,440 - 1,520 nD	0,001 nD		85,-

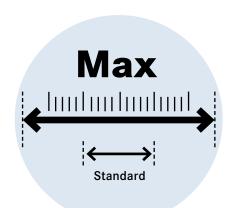









#### Scope of application: Gemmology/Jewellery

The Gem models have a special refracting-index range for jewellery. For this refractometer there is a nice leather bag in the scope of delivery included.

The main scope of applications is:

- Jewellers
- · Training/Education
- Jewellery industry

Model KERN	Scales	Measuring range	Division	ATC	Price excl. of VAT ex works €
ORA 1GG	Refractive index	1,30 - 1,81 nD	0,01 nD		225,-





ORA 1GG



#### Accessory parts: Analogue refractometer - ORA



Model KERN	Description	Price excl. of VAT ex works €
ORA-A1101	Prism coverplate with integrated LED illumination	19,-
ORA-A2103	Leather bag for analog refractometers	19,-
ORA-A2107	Leather bag for Gem refractometers (Spare part)	30,-
ORA-A1001	Calibration liquid – distilled water Volume: 2,5 ml	19,-
ORA-A1002	Contact liquid – Clove oil (for Calibration value 19,6%) Volume: 2,5 ml	19,-
ORA-A1003	Calibration liquid – saturated salt solution Volume: 2,5 ml	19,-
ORA-A1004	Contact liquid – Clove oil (for Calibration value 78,8%) Volume: 2,5 ml	19,-
ORA-A 1005	Calibration block for models ORA 82BB, ORA 3HA, ORA 3HB, ORA 6HA, ORA 6HB , ORA 4RR	19,-
ORA-A1007	Contact liquid – Diiodomethane "Standard" (Refractive index: 1,74 nD) Volume: 2,5 ml	19,-
ORA-A3001	Contact liquid – Diiodomethane "Pro" (Refractive index: 1,79 nD) Volume: 2 ml	35,-
ORA-A1008	Calibration block for model ORA 1GG	19,-
ORA-A2001	Prism coverplate (spare part)	19,-

Rela	Relationship overview – refractometer calibration (analogue)						
Model refractometer	Calibration value	Calibration liquid	Article number liquid	Calibration block	Article number calibration block		
ORA 10BA; ORA 10BB; ORA 18BB; ORA 1WA; ORA 1WB; ORA 20BA; ORA 20BB; ORA 32BA; ORA 32BB; ORA 3SA; ORA 3SB; ORA 3WA; ORA 3WB; ORA 7WA; ORA 7WB; ORA 80BB; ORA 80BE	0 % Brix	distilled water	ORA-A1001	-	-		
ORA 1UA; ORA 1UB	0 % Urea	distilled water	- - - ORA-A1001	-	-		
ORA 4FA; ORA 4FB; ORA 4UA; ORA 4UB	0 °C EG/PG/CW	distilled water		-			
ORA 1SA; ORA 1SB	0 ‰ Salinity	distilled water		-			
ORA 2SA; ORA 2SB	0 % Salt (NaCl)	distilled water		-			
ORA 2AB	0 % Vol (weight)	distilled water		-			
ORA 2PA; ORA 2PB; ORA 5PB	1,000 sg Urine	distilled water		-			
ORA 62BA; ORA 62BB	29,6 % Brix	saturated salt solution	ORA-A1003	-	-		
ORA 3HA; ORA 3HB; ORA 82BB	78,8 % Brix	Clove oil CAS 8000-34-8	ORA-A1004	yes	ORA-A1005		
ORA 4RR	1,4875 nD	Clove oil CAS 8000-34-8	ORA-A1004	yes	ORA-A1005		
ORA 6HA; ORA 6HB	19,6 % Water content	Clove oil CAS 8000-34-8	ORA-A1002	yes	ORA-A1005		
ORA 1GG	1,515 nD	Diiodomethane CAS 90-11-9	ORA-A1007	yes	ORA-A1008		





Transport and storage case



Rear view, screw-on battery compartment cover

#### Digital measurement of refraction index for universal application • ECO refractometer

#### **Features**

- The KERN ORF refractometers are accurate and universal maintenance free digital handheld refractometers
- The large display is easy to read. Mistakes in reading are avoided
- The typical and practical design is suitable for a quick and convenient everyday use and is characterized by its easy-using and robustness
- The large, easy-to-read display with integrated temperature display supports the user to reliably determine the measurement
- The integrated automatic temperature compensation (ATC at ORF 45BE), avoids the manual conversion of the measurement. This allows a quick and efficient usage of the instrument

- Due to the fact that the refractometer has been calibrated at the factory, this guarantees that it can be used immediately for accurately measuring your sample.
- The follwoing accessory-parts are included:
- Calibration liquid
- Pipette
- Storage box
- 2 × AAA batteries
- Leather bag
- Small screwdriver
- Cleaning tissue

#### **Technical data**

- Measurement temperature: 10 °C 40 °C
- Overall dimensions W×D×H 145×67×40 mm
- Net weight approx. 200 g
- Power supply: 2 × AAA (1,5 V)
- Lifetime of the battery: approx. 5.000 measurements
- ATC (Automatic Temperature Compensation), does not apply to the refraction index scale
- Minimum sample volume: 2-3 drops
- Automatic energy management (AUTO-OFF after 90 seconds)



Now also available with calibration certificate, see page 114!

STANDAR





Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 45BE	Brix	0 – 45 %	± 0,2 %	0,1 %	290,-
ORF 45RE	Refractive index	1,3330 - 1,4098 nD	± 0,0003 nD	0,0001 nD	290,-





Transport and storage case



Rear view, screw-on battery compartment cover



IP65: Protected against dust and water splashes

### Digital refractive index measurement for laboratories and the industry for multi-application ► PREMIUM refractometer

#### **Features**

- The KERN ORF refractometers are accurate and universal maintenance free digital handheld refractometers
- The large display is easy to read. Mistakes in reading are avoided
- The typical and practical design is suitable for a quick and convenient everyday use and is characterized by its easy-using and robustness
- The PREMIUM refractometers from the KERN ORF range are protected to international IP65 protection class, against dust and water splashes. After use, you can rinse the refractometer under running water
- · The large, easy-to-read TFT colour display with integrated temperature display supports the user to reliably determine the measurement
- · A large selection of models is available with single or multiple scales. This allows the use in various applications

- · The instrument comes with an optimized software that can show a result in different scales
- The integrated automatic temperature compensation (ATC), avoids the manual conversion of the measurement. This allows a quick and efficient usage of the instrument
- Due to the fact that the refractometer has been calibrated at the factory, this guarantees that it can be used immediately for accurately measuring your sample.
- The follwoing accessory-parts are included:
  - Calibration liquid
- Pipette
- Storage box
- 2 × AAA batteries
- Leather bag
- Small screwdriver
- Cleaning tissue

#### **Technical data**

- Measurement temperature: 5 °C 40 °C
- Overall dimensions W×D×H 145×67×40 mm
- Net weight approx. 200 g
- Power supply: 2 × AAA (1,5 V)
- Lifetime of the battery: approx. 3.750 measurements
- · ATC (Automatic Temperature Compensation), does not apply to the refraction index scale
- Minimum sample volume: 2-3 drops
- · Automatic energy management (AUTO-OFF after 90 seconds)



Now also available with calibration certificate, see page 114!









The following models are particularly suitable for the measurement of the "BRIX" value. They are used to determine the sugar content in food, especially in fruit, vegetables, juice and sweet or soft drinks. In the same ideal way, these refractometers serve in monitoring processes in the industry (coolant monitoring, oils, lubricants and fats). Alternativly, the dispaly can be switched to show the refractive index.

The main scope of applications is:

- Industry: Monitoring of lubricants in machines and quality control
- Food industry: Beverages, fruits and sweets
- · Agriculture: Determination of the degree of ripeness of fruit for quality control in harvesting
- Restaurants and large-scale catering establishment

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 45BM	Brix Refractive index	0 - 45 % 1,3330 - 1,4098 nD	± 0,2 % ± 0,0003 nD	0,1 % 0,0001 nD	340,-
ORF 92BM	Brix Refractive index	58 - 92 % 1,4370 - 1,5233 nD	± 0,2 % ± 0,0003 nD	0,1 % 0,0001 nD	370,-
ORF 85BM	Brix Refractive index	0 - 85 % 1,3330 - 1,5100 nD	± 0,2 % ± 0,0003 nD	0,1 % 0,0001 nD	395,-



#### Scope of application: Honey

The following models are particularly suitable for the measurement of the "BRIX" value, the water content in honey according to the International Honey Commission (IHC2002) and "degrees Baumé" to determine the relative density of liquids. Alternatively the display can be switched to show the refractive index.

- · Beekeeping
- Honey production

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 92HM	Brix Baumé Water content Refractive index	58 - 92 % 38 - 43 °Bé 13 - 25 % 1,4370 - 1,5233 nD	± 0,2 % ± 0,2 °Bé ± 0,2 % ± 0,0003 nD	0,1 % 0,1 °Bé 0,1 % 0,0001 nD	370,-



#### Scope of application: Salt

The following models are particularly suitable to determin the concentration of NaCl (salt) in water. This is often used for the preparation and for the cooking of sauces, bases for pastries, the production of brines (e.g. for white cheese) and the preparation of seafood and marinades for meat. Alternatively the display can be switched to show the refractive index.

The main scope of applications is:

- Food industry
- · Restaurants, and large-scale catering establishment, canteens



Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 3SM	Brix Salt (NaCl) Refractive index	0 - 45 % 0 - 28 % 1,3330 - 1,4100 nD	± 0,2 % ± 0,2 % ± 0,0003 nD	0,1 % 0,1 % 0,0001 nD	340,-

#### Scope of application: Wine

The following models are particularly suitable for the measurement of the sugar content in fruit. It indicates the expected °Alcohol of the fruit. The degree of ripeness of fruit (fruit-sugar) can also be determined, such as e.g. grapes.

The main scope of applications is:

- · Agriculture: Wine-growing (viticulture) and fruit-growing
- · Wine-production
- Must and alcohol production



°Oe = Degree Oechsle, °KMW = Klosterneuburger Most Waage

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 2WM	Mass SW Vol. AP Oechsle KMW (Babo)	0 - 35 % 0 - 22 % 0 - 150 °Oe 0 - 25 °KMW	± 0,2 % ± 0,2 % ± 1 °Oe ± 0,2 °KMW	0,1 % 0,1 % 1 °Oe 0,1 °KMW	340,-

The following models are particularly suitable for the measurement of the specific gravity (sg) in urine, the quantitiy of serum (serumproteine) in urine (doping control among athletes), and the refractive index.

The main scope of applications is:

- Hospitals
- Doctor's surgeries/Physicians
- Medical training institutions
- Nursing homes
- · Sports medicine (doping test)

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 1PM	Serum protein Urine (spec. gravity) Refractive index	0 – 12 g/dl 1,000 – 1,050 sgU 1,3330 – 1,3900 nD	± 0,1 g/dl ± 0,001 sgU ± 0,0003 nD	0,1 g/dl 0,001 sgU 0,001 nD	340,-



#### Scope of application: Industry/Automotive

The following models are particularly suitable for the measurement and determination of AdBlue, glycol concentration (ethylene (EG) and propylene (PG)), battery fluid (BF), urea, the freezing point of fountain solution (CW). Furthermore these models are suitable for the measurement of thermal exchange systems.

The main scope of applications is:

- · Automotive industry: Car-workshops and producers
- · Chemical industry
- · Solar industry: Antifreeze monitoring
- Geothermal industry: Brine-concentration-measurement for ground heat
- Forestry/Lumbermen

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORF 2UM	EG PG BF CW	-50 - 0 °C -50 - 0 °C 1.00 - 1.50 kg/l -40 - 0 °C	± 0,5 °C ± 0,5 °C ± 0,01 kg/l ± 0,5 °C	0,1 °C 0,1 °C 0,01 kg/l 0,1 °C	340,-
ORF 5UM	EG PG Urea CW	-50 - 0 °C -50 - 0 °C 0 - 40 % -40 - 0 °C	± 0,5 °C ± 0,5 °C ± 0,2 % ± 0,5 °C	0,1 °C 0,1 °C 0,1 % 0,1 °C	340,-
ORF 6US	Urea Refractive index	0 - 40 % 1,3330 - 1,4100 nD	± 0,2 % ± 0,0003 nD	0,1 % 0,0001 nD	340,-




11

#### Scope of application: Expert applications

The following model has a special large measuring range for the refractive index.

The main scope of applications is:

• Universal measuring instrument, especially for applications with extra large measuring ranges

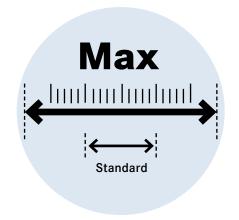

	odel ERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
OF	RF 1RS	Refractive index	1,3330 - 1,5400 nD	± 0,0005 nD	0,0001 nD	430,-

#### Accessory parts: Digital refractometer - ORF

Model KERN	Description	Price excl. of VAT ex works €
ORF-A1005	Prism cover for digital refractometers	25,-
ORA-A1001	Calibration liquid – distilled water Volume: 2,5 ml	19,-
ORA-A1006	Calibration liquid - Triethyl citrate Volume: 2,5 ml	19,-
ORD-A2104	Leather bag for digital refractometer (Spare part)	19,-



Calibration liquid/ Contact liquid

Relationship overview – refractometer calibration (digital)							
Model refractometer	Calibration value	Calibration liquid	Article number liquid	Calibration block	Article number calibration block		
ORF 45BM; ORF 85BM; ORF 3SM	0 % Brix	distilled water	ORA-A1001	-	-		
ORF 2WM	0 °KMW	distilled water	ORA-A1001	-	-		
ORF 1PM; ORF 1RS	1,3330 nD	distilled water	ORA-A1001	-	-		
ORF 2UM; ORF 5UM	0 °C EG/PG/CW	distilled water	ORA-A1001	-	-		
ORF 6US	0 % Urea	distilled water	ORA-A1001	-	-		
ORF 92BM; ORF 92HM	60 % Brix	Triethyl citrate CAS 77-93-0	ORA-A1006	-	-		



Now also available with calibration certificate, see page 114!

## Refractive index measurement for pharmacy, laboratories and industry

#### **Features**

- The KERN ORT refractometers are universal analog Abbe refractometers
- The handy and robust design allows the easy, efficient and sustainable use in everyday life
- The integrated scale allows the use in different applications and provides the best possible security to read the measurement results accurately
- The scope of delivery includes:
- Calibration solution
- Calibration block
- Pipette
- Small screwdriver
- Cleaning tissue
- Digital thermometer
- Accessories are available as options

#### **Technical data**

- Measurement temperature: 20  $^{\circ}\text{C}$
- Overall dimensions W×D×H 180×90×240 mm
- Net weight approx. 1950 g

12

TANDAN








#### Scope of application: Industry/Pharmacy/Laboratory

The following model is a simple yet highly reliable Abbe refractometer with a digital thermometer. Liquid, solid and pasty samples can be evaluated. This refractometer is robust, accurate and easy to use. Optionally a solide aluminium case for transport and storage is available. It measures the refractive index (nD).

The main scope of applications is:

- Sugar industry: for example cane sugar
- Pharmacy
- Beverage industry
- Food industry
- Chemical industry
- · Oil industry/Refinery
- Laboratories
- Training

Model KERN	Scales	Measuring range	Accuracy	Division	Price excl. of VAT ex works €
ORT 1RS	Brix Refractive index	0 - 95 % 1,3000 - 1,7000 nD	± 0,1 % ± 0,0002 nD	0,25 % 0,0005 nD	690,-

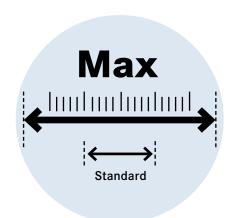



#### Accessory parts: Abbe refractometer - ORT

Model KERN	Description	Price excl. of VAT ex works €
ORA-A1102	Aluminium suitcase Dimension: 310×120×240 mm, weight: 1300 g	80,-
ORA-A2266	Digital thermometer (0 °C/50 °C) (Spare part)	60,-
ORA-A2267	Calibration block for ORT 1RS	40,-
ORA-A1107	Contact liquid – Alpha-Bromonaphthalene (Refractive index: 1,65 nD) Volume: 2,5 ml	25,-
ORA-A3001	Contact liquid - Diiodomethane "Pro" (Refractive index: 1,79 nD) Volume: 2 ml	35,-



Transport and storage case ORA-A1102



Calibration block ORA-A2267

Relationship overview – refractometer calibration (Abbe)					
Model refractometer	Calibration value	Calibration liquid	Article number liquid	Calibration block	Article number calibration block
ORT 1RS	engraved on the block (nD value)	Alpha-Bromnaphthalene CAS 90-11-9	ORA-A1107	yes	ORA-A2267





## Your partner for calibration services, management of test equipment and support

#### **Features**

- Any analogue or digital refractometer will only give correct results if it is checked regularly, i.e. calibrated correctly and adjusted when required. A refractometer or another measuring device is only a reliable measuring and checking tool if it is calibrated and this calibration is documented as part of a quality procedure
- Measuring "correctly" is of elementary significance, as it is not unusual for inaccurate or "wrong" measurements to have expensive economic consequences.
   Calibration or establishing the accuracy of checking equipment must therefore be carried out by laboratories throughout the world
- In the context of standard requirements for monitoring checking equipment, every company with a Quality Management system is obliged to test and document its measuring equipment at regular intervals
- The refractometer calibration certificate documents the intended measuring functionality and confirms the measuring accuracy of your refractometer to you

#### Important

- Refractive index standard traceable to SRM<sup>1</sup> of NIST<sup>2</sup> and PTB<sup>3</sup>
- This service is not possible for the following refractometer models:
  - ORA 6HA
- ORA 1GG
- Calibration of products from other manufacturers is possible on request

<sup>&</sup>lt;sup>3</sup>Physikalisch-Technische Bundesanstalt (German metrology institute)

Mode		Description	Price excl. of VAT ex works €
961-29	90	Calibration certificate for refractometers on initial calibration	99,-
961-29	90R	Calibration certificate for refractometers on recalibration	99,-

<sup>&</sup>lt;sup>1</sup>Standard reference material


<sup>&</sup>lt;sup>2</sup>National Institute of Standards and Technology

# www.kern-sohn.com

## Straight to the target!

Get straight to the right product. Use our new topic area search. Here you will quickly find products which suit your particular area of expertise



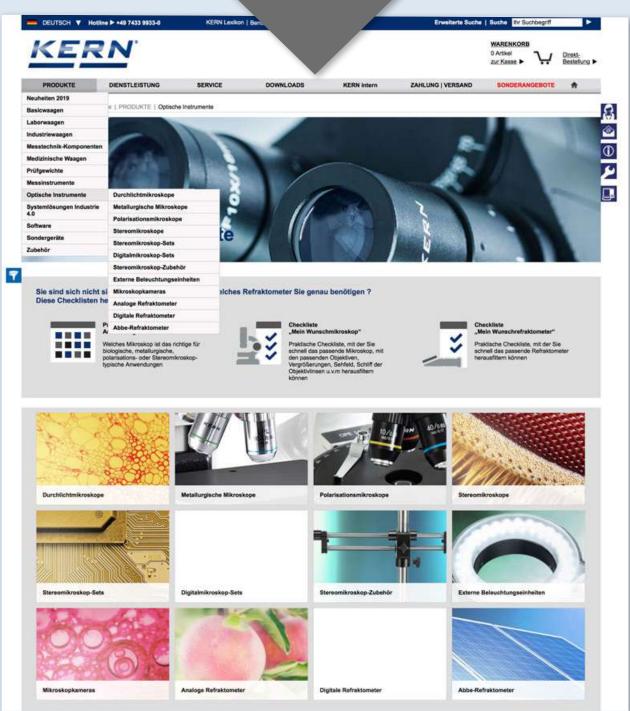

### Checklist for your microscope – your requirements

(1	Which kind of microscope of	do you need?				
			D 447 446			
	Compound microscope	primarily used for transparent/translucent preparation	Page 116 - 118			
Ц	Stereo microscope	surface observation with 3-dimensional optic with slow/medium magnification	Page 119 – 123			
	Phase contrast microscope	Preparations with minimal contrast/very translucent	Page 116 - 118			
	Fluorescence microscope	fluorescent structures, which are specific coloured or auto coloured	Page 116 - 118			
	Polarisation microscope	Preparations with refraction (anisotropic) for example Crystal	Page 116 - 118			
	Metallurgical microscope	surface observation of components, materials and minerals	Page 116 - 118			
	Inverted microscope	used primarily for culture fessel from cell culture/for very thick preparation	Page 116 – 118			
	te your intended use/ scribe your application:					
	te your previous del/manufacturer: (if available)					
_						
	te your min. and max. gnification:					
ma	giiiidatioii.					
(2	What type of eyepiece tube	do you need for your application?				
	Monocular eyepiece tube	view with one eye = 1 eyepiece available				
	Binocular eyepiece tube	view with both eyes = 2 eyepieces available				
	Trinocular eyepiece tube	view with both eyes + additional option to adapt a camera				
	Digital eyepiece tube	view with both eyes + integrated camera				
Att	ention: look also at point 20 - Do you n	eed a camera?				
Ad	ditional comments:					
(3	Which illumination do you r	need for your application?				
	Halogen transmitted illumination	very good illumination/also suitable for dark field and phase contrast				
	LED transmitted illumination	extremely long life time/no heat generation				
	Halogen reflecting illumination	additional illumination, e.g. for Polarisation and metallurgical microscop	oes			
	LED incident illumination	only for stereo microscopes				
	External illumination	external illumination could be ordered additionally, for example ring illu	mination unit, swan neck			
		(cold light source), as Accessories)				
1	avaatina fasta					
	eresting facts	tht microcopy, bosquae they have a hatter bricktone				
		tht microscopy, because they have a better brightness.				
	_	life time and the advantage that there is no heat generation. in our stereo microscopes as standard illumination.				
•	,	,				
Ade	ditional comments:					
-						

bo you need Roemer mummation:			
	no fixed, pre-centred Koehler illumination full Koehler illumination	condenser is centred, can be height-adjusted and focussed, field diaphragm/aperture diaphragm available.  condenser can be fully centred and focussed, field diaphragm/aperture diaphragm available.	
State your intended use/ Describe your application:			
5	How many objectives would	you like to use?	
	4 objectives	Quadplex nosepiece	
	5 objectives	Quintuple nosepiece	
6	What magnification (objective	ves) do you need?	
_			
		when using the 10× magnification eyepiece	
		when using the 10× magnification eyepiece	
		when using the 10× magnification eyepiece	
		when using the 10× magnification eyepiece	
	100× objective (1000× magnification)	when using the 10× magnification eyepiece	
14			
	eresting facts		
► M	agnification formula: objective magnification	ation × eyepiece magnification = Total magnification	
0.			
Sta	te the magnification you require		
۱ ـ ۸	District where continues to the sale of		
Add	litional phase contrast objective		
	NA/lead acceptance de conserva de face		
<u>(</u>	What quality do you need for	r the objective?	
П	Achromatic	DIN standard chicativas	
	Plan achromatic	DIN standard objectives	
	Infinity E-Plan/Semi Plan	DIN standard objectives infinitely corrected objectives for professional methods	
	Achromatic Infinity Plan	infinitely corrected objectives for professional methods	
	Actioniatic mining Flan	minitery corrected objectives for professional methods	
Additional comments:			
, luu	ntional comments.		

8 What eyepiece diameter (visual field) and what eyepiece magnification do you need?				
10× magnification	Dioptre adjustment			
□ Ø 18 mm	☐ Yes, on one side			
□ Ø 18 mm with pointer needle	☐ Yes, on both sides			
□ Ø 18 mm with 0.1 mm scale	□ No			
□ Ø 20 mm				
□ Ø 20 mm with 0.1 mm scale				
2 2 min with 6.1 min Source				
Further magnifications possible:				
(State the magnification you require:)				
(otate the magnification you require.)				
Do you need a camera to say	ve the documents?			
□ yes				
□ no				
Interesting facts				
=	have to very a Compount adopted to adopt a page and			
with a trinocular microscope, you always	have to use a C-mount adapter to adapt a camera!			
Additional comments:				
(e.g. required number of Mpx etc.)				
(**************************************				
Danis and any fourth or form	-1 <sup>1</sup> 0			
Do you need any further fun	ctions?			
☐ Dark field unit				
☐ Polarisation unit				
☐ Fluorescent unit				
☐ Phase-contrast unit				
☐ Colour filter				
☐ Additional objectives				
Additional comments:				
Statement of phase contrast magnification:				
Statement Fluorescence-channel				
(colour UV/V/B/G):				
_				
Further technical characteri	etice:			
i di tilei tecimical characteri	31103.			
State your requirements:				

### Technical requirements of stereo microscope

What type of eyepiece tube do you need for your application?			
_			
	Binocular eyepiece tube	view with both eyes, two eyepieces	
	Trinocular eyepiece tube	view with both eyes and additional option to adapt a camera	
Att	ention: look also at point 20 – Do y	ou need a camera?	
Add	ditional comments:		
13	Please select the require	ed optical system?	
	Greenough	beam paths which are completely separate from each other	
Ц	Parallel/ABBE	beam paths which are completely separate from each other which run parallel	
Δda	litional comments:		
Aut	attional comments.		
14	Which illumination do yo	ou need for your application?	
14			
14 	None	stereo microscope without illumination	
	None Incident illumination	stereo microscope without illumination incident illumination e.g. LED or halogen	
	None Incident illumination Transmitted illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples	
	None Incident illumination Transmitted illumination Coaxial illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus	
	None Incident illumination Transmitted illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples	
	None Incident illumination Transmitted illumination Coaxial illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories	
Add	None Incident illumination Transmitted illumination Coaxial illumination External illumination  ditional comments:  What type of magnificat	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories	
	None Incident illumination Transmitted illumination Coaxial illumination External illumination ditional comments:  What type of magnificat Rotation objective	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories  ion do you need?  changing the magnification by rotating the objective	
Add	None Incident illumination Transmitted illumination Coaxial illumination External illumination  ditional comments:  What type of magnificat	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories	
Add	None Incident illumination Transmitted illumination Coaxial illumination External illumination ditional comments:  What type of magnificat Rotation objective	stereo microscope without illumination incident illumination e.g. LED or halogen additional illumination for translucent samples integrated coaxial illumination for selective depth of focus external illumination could be ordered additionally, for example ring illumination unit, swan neck (cold light source), as Accessories  ion do you need?  changing the magnification by rotating the objective	

What magnification do you need?				
☐ Minimum:	□ Maximum:			
Additional comments:				
Interesting facts  • Magnification formula: Eyepiece magnification × objective magnification (zoom) = Total magnification				
What eyepiece diameter (vis	ual field) and what eyepiece magnification do you need?			
10× magnification  □ Ø 20 mm  □ Ø 22 mm  □ Ø 23 mm	Dioptre adjustment  ☐ Yes, on one side ☐ Yes, on both sides			
Further magnifications possible: (State the magnification you require:)				
What working distance do yo	ou need?			
☐ Minimum:	mm			
Additional comments:				
Interesting facts ► The working distance is the distance between the objective and the sample.				
What size of field of view do you need?				
☐ Minimum:	mm			
Additional comments:				

### Interesting facts

► The field of view is the section which is shown through the magnification. If the magnification (Zoom) is very high, the field of view will be reduced. By magnifying and focusing a specific section, it is not possible to capture the whole sample.

Do you need a camera to save the documents?		
□ yes □ no		
Interesting facts  ► With a trinocular microscope, you always	have to use a C-mount adapter to adapt a camera!	
Additional comments: (e.g. required number of Mpx etc.)		
Do you need any further fun	ctions?	
<ul> <li>□ Dark field unit</li> <li>□ Stand inlays         (preparation-background)</li> <li>□ Universal stand</li> <li>□ Mechanical bench</li> </ul>	e.g. glass, opaque glass, black, white	
Additional comments:		
Further technical characteric	stics:	
State your requirements:		
Please fill in your contact, the	nat we could make you an offer for a suitable microscope	
Customer number:		
Company:		
Surname, first name:		
Street:		
Postcode/Area:		
Country:		
Tel.:		
Fax:		
E-Mail:		

### **Checklist for your refractometer – your requirements**

Which kind of refractometer do you need?				
	Analogue refractometer Digital refractometer Analogue ABBE refractometer	hand-held device for rapid analysis/can be used for mobile applications digital, hand-held device for rapid analysis/can be used for mobile applications stand benchtop device/refraction index & Brix – Measuring for all applications		
State your intended use/ Describe your application:				
2	State your area of applic	ation		
	Sugar/lubricants	determining the sugar content in e.g. fruit, vegetables, juice, drinks which contain sugar, lubricants for drills, milling machines and lathes		
	Honey	determining the sugar content, water content and relative density of liquids		
	Salt	determining the sodium chloride in water (salinity) and determining the salt in water		
	Wine	determining the alcohol content, degree of ripeness through fruit sugar		
	Urine	Determination of specific urine weight (density), the serum content (serum protein in urine) and the refraction index		
	Industry/motor vehicles Refractive index	Determination of glycol concentration (ethylene, propylene), of battery fluid, screen wash determining the refraction index of different substances		
What sort of results are you expecting: (to determine the scale)				
3	Do vou need automatic t	emperature compensation?		
	<u>)                                    </u>			
	yes	ATC allows you take measurements with different environmental, device and sample temperatures between 10 $^{\circ}\text{C}$ and 30 $^{\circ}\text{C}$		
	no	without ATC the environmental, device and sample temperatures must be 20° C in order to get a precise measuring result. If the table is different, the result must be corrected manually.		
	Please fill in your contac	t, that we could make you an offer for a suitable refractometer		
	Customer number:			
	Company:			
	Surname, first name:			
	Street:			
	Postcode/Area:			
	Country:			
	Tel.:			
	Fax:			
	E-Mail:			